

UNCLASSIFIED

AD 295 672

*Reproduced
by the*

**ARMED SERVICES TECHNICAL INFORMATION AGENCY
ARLINGTON HALL STATION
ARLINGTON 12, VIRGINIA**



UNCLASSIFIED

NOTICE: When government or other drawings, specifications or other data are used for any purpose other than in connection with a definitely related government procurement operation, the U. S. Government thereby incurs no responsibility, nor any obligation whatsoever; and the fact that the Government may have formulated, furnished, or in any way supplied the said drawings, specifications, or other data is not to be regarded by implication or otherwise as in any manner licensing the holder or any other person or corporation, or conveying any rights or permission to manufacture, use or sell any patented invention that may in any way be related thereto.

SSD-TDR-62-203

295 672


**AN ACETYLATION METHOD FOR
DETERMINING MIXTURES OF EITHER
HYDRAZINE-1,1-DIMETHYLHYDRAZINE
OR
MONOMETHYLHYDRAZINE-1,1-DIMETHYLHYDRAZINE**

TECHNICAL DOCUMENTARY REPORT NO. SSD-TDR-62-203

DECEMBER 1962

**ROCKET RESEARCH LABORATORIES
CHEMICAL AND MATERIALS BRANCH
EDWARDS, CALIFORNIA**

**HEADQUARTERS
SPACE SYSTEMS DIVISION
AIR FORCE SYSTEMS COMMAND
UNITED STATES AIR FORCE**

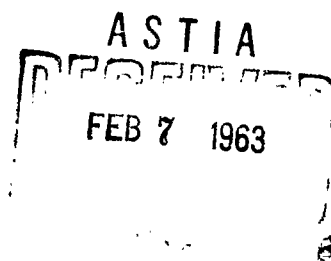

NOTICES

When Government drawings, specifications, or other data are used for any purpose other than in connection with a definitely related Government procurement operation, the United States Government thereby incurs no responsibility nor any obligation whatsoever; and the fact that the Government may have formulated, furnished, or in any way supplied the said drawings, specifications, or other data, is not to be regarded by implication or otherwise as in any manner licensing the holder or any other person or corporation, or conveying any rights or permission to manufacture, use, or sell any patented invention that may in any way be related thereto.

Qualified requesters may obtain copies of this report from the Armed Services Technical Information Agency, (ASTIA), Arlington Hall Station, Arlington 12, Virginia.

SSD-TDR-62-203

AN ACETYLATION METHOD FOR DETERMINING
MIXTURES OF EITHER HYDRAZINE-1, 1-DIMETHYLHYDRAZINE
OR MONOMETHYLHYDRAZINE-1, 1-DIMETHYLHYDRAZINE

TECHNICAL DOCUMENTARY REPORT
NO. SSD-TDR-62-203


DECEMBER 1962

by

HUGH E. MALONE

ROBERT A. BIGGERS

This report has been reviewed and approved.

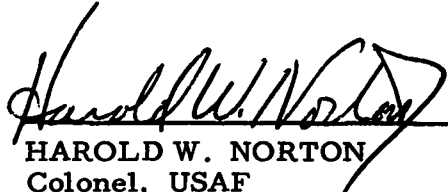

HAROLD W. NORTON
Colonel, USAF
Commander

Rocket Research Laboratories
Chemical and Materials Branch
Edwards, California

ABSTRACT

A method is presented for determining mixtures of either hydrazine-1, 1-dimethylhydrazine (N_2H_4 /UDMH) or monomethylhydrazine-1, 1-dimethylhydrazine (MMH/UDMH). This method is based on the difference in reaction rates for acetylating hydrazines. Both N_2H_4 and MMH react immediately while the acetylation of UDMH proceeds slowly. The reactions are conducted in acetonitrile-acetic acid medium. Two titrations are required: one measures the total basicity of the hydrazine mixtures; the other titration measures UDMH after the N_2H_4 or MMH has reacted with Ac_2O .

TABLE OF CONTENTS

	<u>Page</u>
INTRODUCTION	1
BACKGROUND	1
EXPERIMENTAL	2
Effect of Acetic Anhydride	3
Effect of Water	4
Potentiometric and Colorimetric Endpoints	5
PROCEDURE	6
Preparation of the Sample	6
Determination of Total Basicity of the Hydrazine- UDMH Mixture	6
Determination of 1, 1-Dimethylhydrazine	6
CALCULATIONS	6
REAGENTS	7
RESULTS	7
SUMMARY	9
BIBLIOGRAPHY	14

LIST OF ILLUSTRATIONS

Figure

I.	Effect of Time and Temperature on the Reaction Rate of UDMH with Acetic Anhydride	11
II.	The Effect of Acetic Anhydride on Hydrazines	12
III.	Typical Titration Curve for Both Titrations	13

AN ACETYLATION METHOD FOR DETERMINING MIXTURES OF EITHER HYDRAZINE-1, 1-DIMETHYLHYDRAZINE OR MONOMETHYLHYDRAZINE-1, 1-DIMETHYLHYDRAZINE

INTRODUCTION.

Hydrazine (N_2H_4), monomethylhydrazine (MMH), and 1, 1-dimethylhydrazine (UDMH) are important both individually and in mixtures as rocket fuels. A 50-50 percent mixture of N_2H_4 /UDMH is used as the fuel blend for the Titan II. Though mixtures of MMH/UDMH are discussed herein, they are not presently being used in any weapon system.

The purpose of this report is to present a method for determining mixtures of either N_2H_4 /UDMH and MMH/UDMH based on observed differences in reaction times with acetic anhydride.

BACKGROUND.

Several methods have been developed for analysis of mixtures of hydrazine-1, 1-dimethylhydrazine. These methods use gas chromatographic techniques (1, 5), standard iodate titration in an aqueous medium (9), or the selective reaction of hydrazine with salicylaldehyde and subsequent titration of basicity with perchloric acid in non-aqueous medium (6). This latter method is presently used in Specification MIL-P-27402. It has the end point disadvantage encountered in most weak base-strong acid titrations.

Acetic anhydride used herein for acetylating hydrazines was used primarily as a solvent for titrations of amides (10), weak bases (3), and nitrogen bases (8). It has been used also to improve the end point in non-aqueous titrations (10). Fritz and Schenk (4) used acetic anhydride for acid catalyzed acetylation of hydroxyl compounds in ethyl acetate.


Although the reactions for acetylating hydrazines are known (2) and reported (7), the authors are unaware of any efforts to determine hydrazine and substituted hydrazines by this method.

EXPERIMENTAL.

Reactions of hydrazine, monomethylhydrazine and 1, 1-dimethylhydrazine with acetic anhydride. Reaction products are shown in the following equations:


With N_2H_4 -1, 2 diacetylhydrazine:

Equation (1)

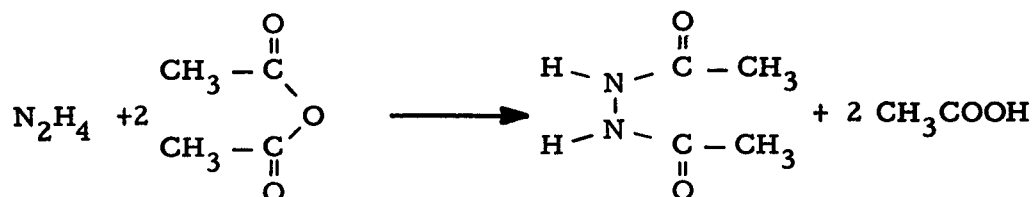


With MMH-Monomethyl 1, 2 diacetylhydrazine:

Equation (2)

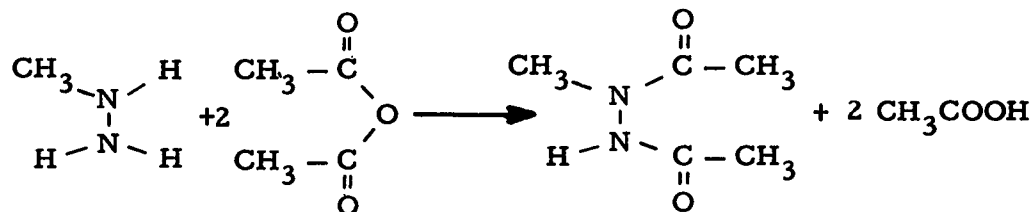

With UDMH-1, 1-Dimethyl 2 acetylhydrazine:

Equation (3)

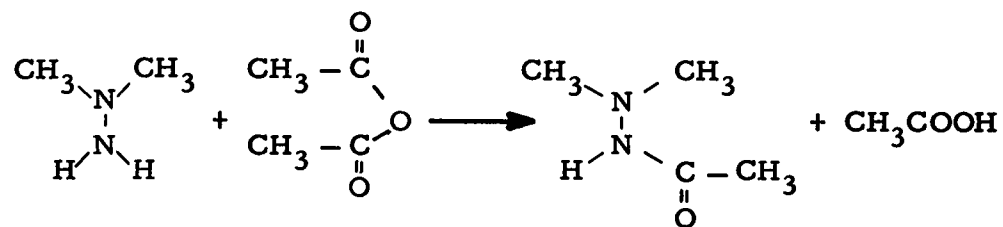

The rates of the reaction differ greatly as both the N_2H_4 and MMH are rendered neutral to perchloric acid titration almost immediately upon UDMH

reacts slowly in acetic acid medium. The product of equation (3) occurs only after a considerable time lapse at room temperature. Figure I shows the effect of time and temperature on the UDMH-acetic anhydride reaction. The reaction, as indicated by a decrease in milliliters of HClO_4 is slow until the temperature is increased.

Effect of Acetic Anhydride.

As the quantity of acetic anhydride is increased in samples of N_2H_4 , MMH, and UDMH, the milliliters of HClO_4 titrant required for neutralization decreases much faster with N_2H_4 and MMH. The UDMH is slow to react even when a considerable excess of acetic anhydride is added. This effect is shown in Figure II. These curves were plotted from titration, with perchloric acid, of individual samples containing a 1 ml aliquot of each hydrazine compound per 50 ml acetic acid. Duplicate samples were treated with 1 ml aliquots of acetic anhydride samples containing 1.0 ml - 6.0 ml acetic anhydride in 50 ml acetic acid. The two curves for both MMH and N_2H_4 substantiate the molar ratio presented in equation (1) and (2).

An analysis, after acetylation, of the individual hydrazines taken from stored drums is shown in Table I.

TABLE I.
CORRECTION VALUES FOR HClO_4 TITRATIONS OF HYDRAZINES

Sample No.	N_2H_4 ml HClO_4 after Ac_2O	UDMH ml HClO_4 after Ac_2O^1	MMH ml HClO_4 after Ac_2O
1	0.20	0.04	0.14
2	0.23	0.03	0.10
3	0.20	0.06	
4	0.21		
5	0.20		
6	0.17		
7	0.05		

¹ml HClO_4 less than total alkalinity value. Sample size 0.025 g.

N_2H_4 drums 1, 2, and 3 were stored for over five years, while drums 4 and 6 were stored for over two years. Drum 5 was stored only six months. Sample marked 7 was redistilled over a 75 plate distillation column. The UDMH and MMH drums were stored for one year.

The values presented above result from the presence of contaminants in the hydrazines. Aniline and ammonia are the usual titratable contaminants in N_2H_4 . Nitroso-dimethylamine, ammonia, dimethylamine, methylene dimethylhydrazine, and methylhydrazine are the titratable contaminants in UDMH (12). Of these contaminants, aniline and methylhydrazine react similarly to N_2H_4 and are nontitratable after acetylation. Methylene dimethylhydrazine, dimethylamine, ammonia, and nitrosodimethylamine will be acetylated at different rates and could titrate as UDMH. Further work is necessary to establish the rates of reaction for these contaminants. However, their presence in the hydrazines are not significant except for dimethylamine (0.7%) and nitrosodimethylamine (0.3%) in UDMH and ammonia (0.3%) in N_2H_4 .

The contaminants usually present in MMH have not been identified as yet. Until they are, titrations of the individual hydrazines are necessary.

Preliminary observations using other acetylation compounds such as acetyl chloride, benzoyl chloride, maleic anhydride, and benzoic anhydride with UDMH indicate that only the two anhydrides can be used as possible substitutes for acetic anhydride.

Effect of Water.

The water content in N_2H_4 -UDMH samples was increased from 0.99% to 7.84% and the samples were analyzed by the procedure presented. Table II shows that increasing the water content does not interfere with the titration.

TABLE II

Percent Water	Sample Weight N_2H_4 /UDMH before H_2O added	Total Weight of Sample	ml HClO_4 Total Alkalinity	ml HClO_4 after Acetylation
0.99	0.8918	0.9007	4.24	1.42
2.31	0.8913	0.9124	4.24	1.42
4.65	0.8903	0.9337	4.24	1.42
7.84	0.8818	0.9570	4.19	1.40

The large excess of acetic anhydride in the N_2H_4 reaction reacts with any water to form acetic acid and eliminates the water interference problem encountered in most non-aqueous titration procedures.

Potentiometric and Colorimetric Endpoints.

Both the potentiometric and quinaldine red indicator endpoints were evaluated. A Beckman zeromatic pH meter with glass-calomel electrodes was used to determine the stoichiometric endpoints for both the total alkalinity titration and the acetic anhydride titration. Figure III shows a typical curve obtained in these titrations. The quinaldine red endpoint changes from red (actually pink) to colorless approximately 0.01 ml after the greatest change in potential occurs, dictated by the millivolt change plot shown in Figure IV. The indicator and potentiometer endpoints can be improved appreciably by using acetonitrile as the solvent in both titrations. For the plots shown an acetonitrileacetic acid mixture (4:1 V/V) was used.

The endpoints can also be greatly improved by using either dioxane or nitromethane as solvents. However, we have encountered considerable difficulty in obtaining dioxane with a low blank value. The nitromethane is such a good monopropellant that we avoid using it as a solvent except when absolutely necessary.

PROCEDURE.

Preparation of Sample.

Pipette 0.4 ml of the N_2H_4 -UDMH mixture into a tared, cooled 50-ml volumetric flask containing 40 ml of acetic acid. Weigh to the nearest 0.1 mg, obtaining the sample weight by difference. Dilute to the 50-ml mark with acetic acid and mix thoroughly.

Determination of the Total Basicity of the Hydrazine-UDMH Mixture -

Pipette a 5 ml aliquot from the prepared sample into a 50-ml beaker containing 20 ml of acetonitrile. Add 5 drops of quinaldine red indicator and titrate with 0.1 N perchloric acid in acetic acid until the red color of the indicator changes to colorless. Indicate this amount for value "A" in calculations. Titrate a blank of acetonitrile and quinaldine red indicator to the colorless endpoint.

Determination of 1, 1-Dimethylhydrazine - Pipette a 5 ml aliquot of the prepared sample into a 50-ml beaker containing 20 ml of acetonitrile-acetic acid (4:1 V/V). Add 2 ml of acetic anhydride and 5 drops of quinaldine red indicator. Titrate with 0.1 N perchloric acid in acetic until the red color of the indicator changes to colorless. Indicate this amount for value "B" in calculations. Titrate a blank of acetonitrile-acetic acid, acetic anhydride and quinaldine red indicator to the colorless endpoint.

CALCULATIONS.

$$\% \text{ Hydrazine} = \frac{(A-a) - (B-b) \times N \times M_1}{w} \times 10$$

$$\% \text{ UDMH} = \frac{(B-b) \times N \times M_2}{w} \times 10$$

$$\% H_2O = 100 - (\% \text{ Hydrazine} + \% \text{ UDMH})$$

where:

- A = ml perchloric acid for total basicity.
- B = ml perchloric acid for the anhydride addition.
- a = ml perchloric acid for blank of acetonitrile and indicator.
- b = ml perchloric acid for blank of acetonitrile-acetic acid, acetic anhydride and indicator.
- N = normality of perchloric acid.
- M_1 = milliequivalent weight of hydrazine (0.03205).
- M_2 = milliequivalent weight of UDMH (0.06010).
- w = weight of sample

REAGENTS

Acetic acid - glacial reagent grade.

Acetic anhydride - ACS reagent grade.

Acetonitrile - ACS reagent grade.

Quinaldine red - 0.2 percent in acetic acid.

0.1 N perchloric acid in acetic acid - dissolve 8.5 ml of 72% perchloric acid in 1 liter of glacial acetic acid and standardize against primary reagent potassium acid phthalate.

RESULTS.

Several samples of N_2H_4 -UDMH mixtures were prepared and analyzed. Table III presents the theoretical, experimental, and variation from calculated percent obtained for samples ranging from 10% N_2H_4 -90% UDMH to 90% N_2H_4 -10% UDMH.

TABLE III

<u>Experimental</u>		<u>Theoretical</u>		<u>Variation in %</u>	
<u>N₂H₄</u>	<u>UDMH</u>	<u>N₂H₄</u>	<u>UDMH</u>	<u>N₂H₄</u>	<u>UDMH</u>
92.18	7.58	92.31	7.33	-0.13	+0.23
96.34	13.03	86.61	13.05	-0.27	-0.02
73.98	25.43	74.31	25.38	-0.33	+0.05
65.59	34.36	65.29	34.67	+0.30	-0.31
55.99	44.35	55.97	44.03	+0.02	+0.32
57.79	42.03	58.18	41.82	-0.39	+0.21
55.12	45.13	55.20	44.80	-0.08	+0.33
55.77	43.80	55.91	44.09	-0.14	-0.29
50.40	48.05	50.30	48.25	+0.10	-0.20
20.65	78.36	20.25	78.80	+0.40	-0.44
12.20	87.63	11.69	88.17	+0.51	-0.54

Standard Deviation N₂H₄ 0.28

UDMH 0.32

The percent error is greater for samples containing larger UDMH content. This can be attributed to the fact that most of the impurities in the mixtures will react similarly to UDMH. In the above samples of N₂H₄/UDMH mixtures, the correction values were not used. The correction values for both the N₂H₄ and UDMH were 0.04 ml.

Several samples of MMH-UDMH mixtures were prepared and analyzed.

Table IV presents the theoretical, experimental, and variation from calculated value.

TABLE IV

<u>Experimental</u>		<u>Theoretical</u>		<u>Variation in %</u>	
<u>MMH</u>	<u>UDMH</u>	<u>MMH</u>	<u>UDMH</u>	<u>MMH</u>	<u>UDMH</u>
41.56	56.54	41.73	57.04	-0.17	-0.50
50.42	49.24	51.06	48.19	-0.64	+1.05
40.07	55.88	40.47	56.06	-0.40	-0.18
47.01	49.84	47.64	51.19	-0.63	-0.35
62.97	34.47	63.12	34.92	-0.15	-0.45
63.57	36.07	63.89	35.34	-0.32	+0.73

Standard Deviation MMH 0.43
UDMH 0.61

When the correction values discussed in Table I are made to MMH/UDMH samples, the following results are obtained:

TABLE V

<u>Experimental</u>		<u>Theoretical</u>		<u>Variations in %</u>	
<u>MMH</u>	<u>UDMH</u>	<u>MMH</u>	<u>UDMH</u>	<u>MMH</u>	<u>UDMH</u>
41.85	57.10	41.73	57.04	+0.12	+0.06
50.78	48.06	51.06	48.19	-0.18	0.13

We encounter difficulty in the analysis of the hydrazine mixtures when the correction values for the individual hydrazines are unknown. The correction values are applied by equating the value for the individual hydrazine proportional to the total concentration hydrazines in the mixture. The UDMH correction value (resulting from MMH type contaminants) is added to the milliliters acid required for either the N_2H_4 or MMH titration. The N_2H_4 or MMH correction value (resulting from UDMH type contaminants) is subtracted from the milliliters acid required for the UDMH titration.

When the correction values become excessive, the results will be unsatisfactory. Values this high are not usually encountered in stored hydrazines. We presume that the high-value MMH samples reported previously were contaminated accidentally and are not representative correction values normally encountered.

SUMMARY.

Detailed procedures for potentiometric or colorimetric titration in acetonitrile/acetic acid solvent with perchloric acid titrant are presented. Mixtures of either UDMH/ N_2H_4 or UDMH/MMH can be analyzed in this system. The procedural basis is the difference in reaction rates for acetylating hydrazines with acetic anhydride.

Analysis of N_2H_4 /UDMH mixtures gave standard deviations of 0.28 for N_2H_4 and 0.32 for UDMH. Standard deviations found for analysis of MMH/UDMH were 0.43 for MMH and 0.61 for UDMH. The standard deviation decreases when correction values are applied.

RECOMMENDATION.

The acetylation procedure for analysis of N_2H_4 /UDMH is recommended to supersede the aldehyde procedure currently required in Specification MIL-P-27402.

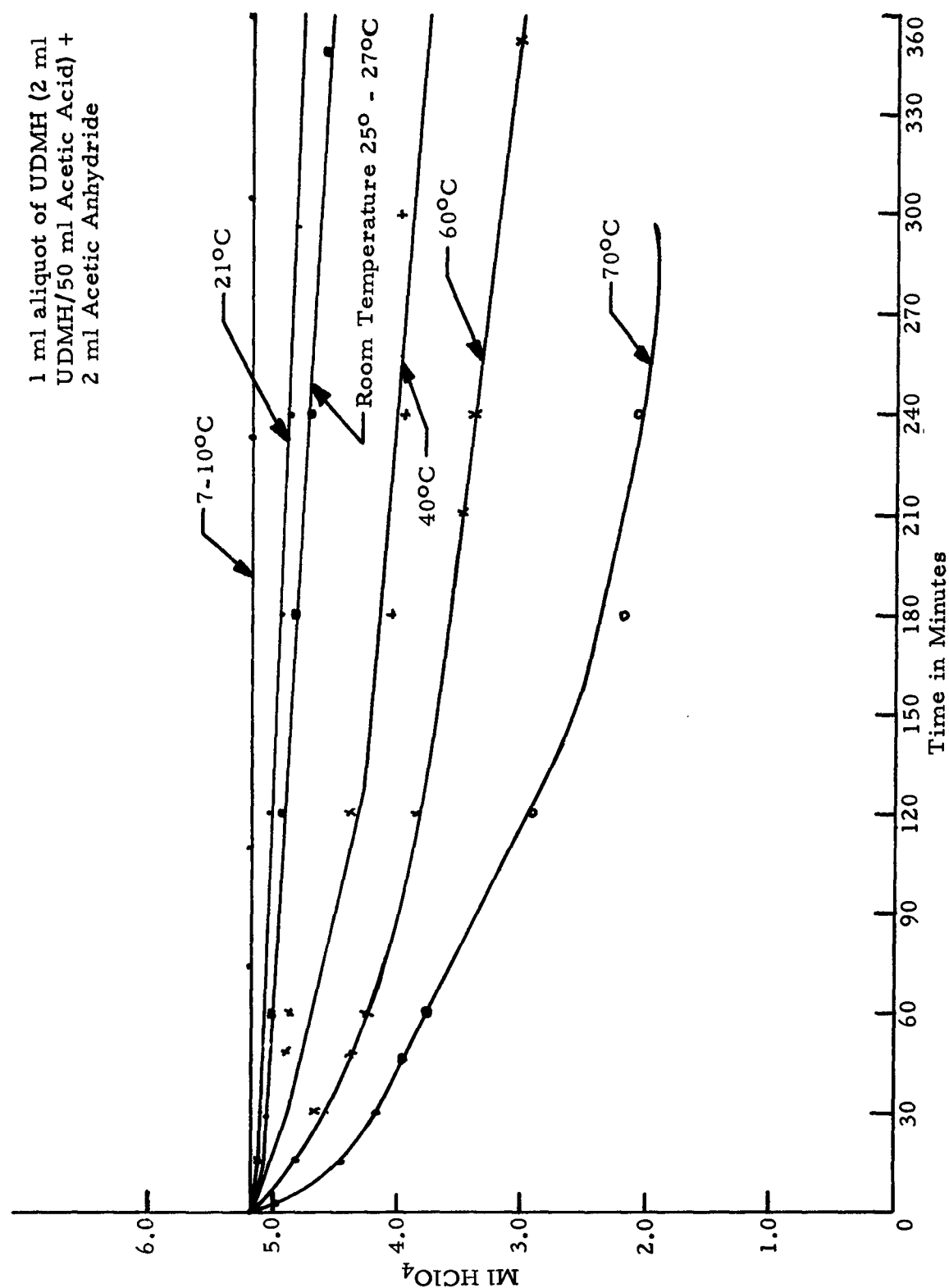


Figure 1. Effect of Time and Temperature on the Reaction Rate of UDMH with Acetic Anhydride.

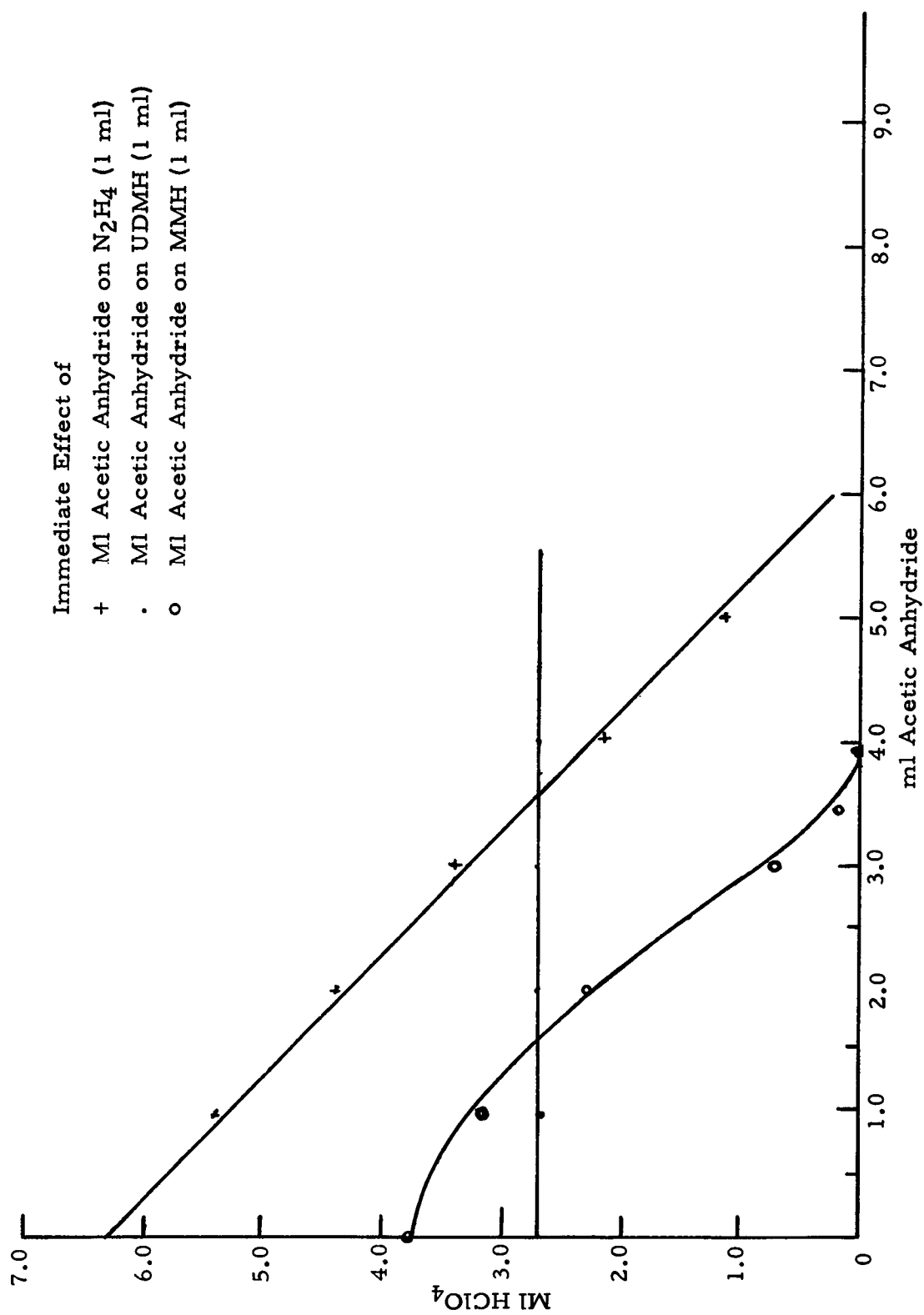


Figure II. The Effect of Acetic Anhydride on Hydrazines

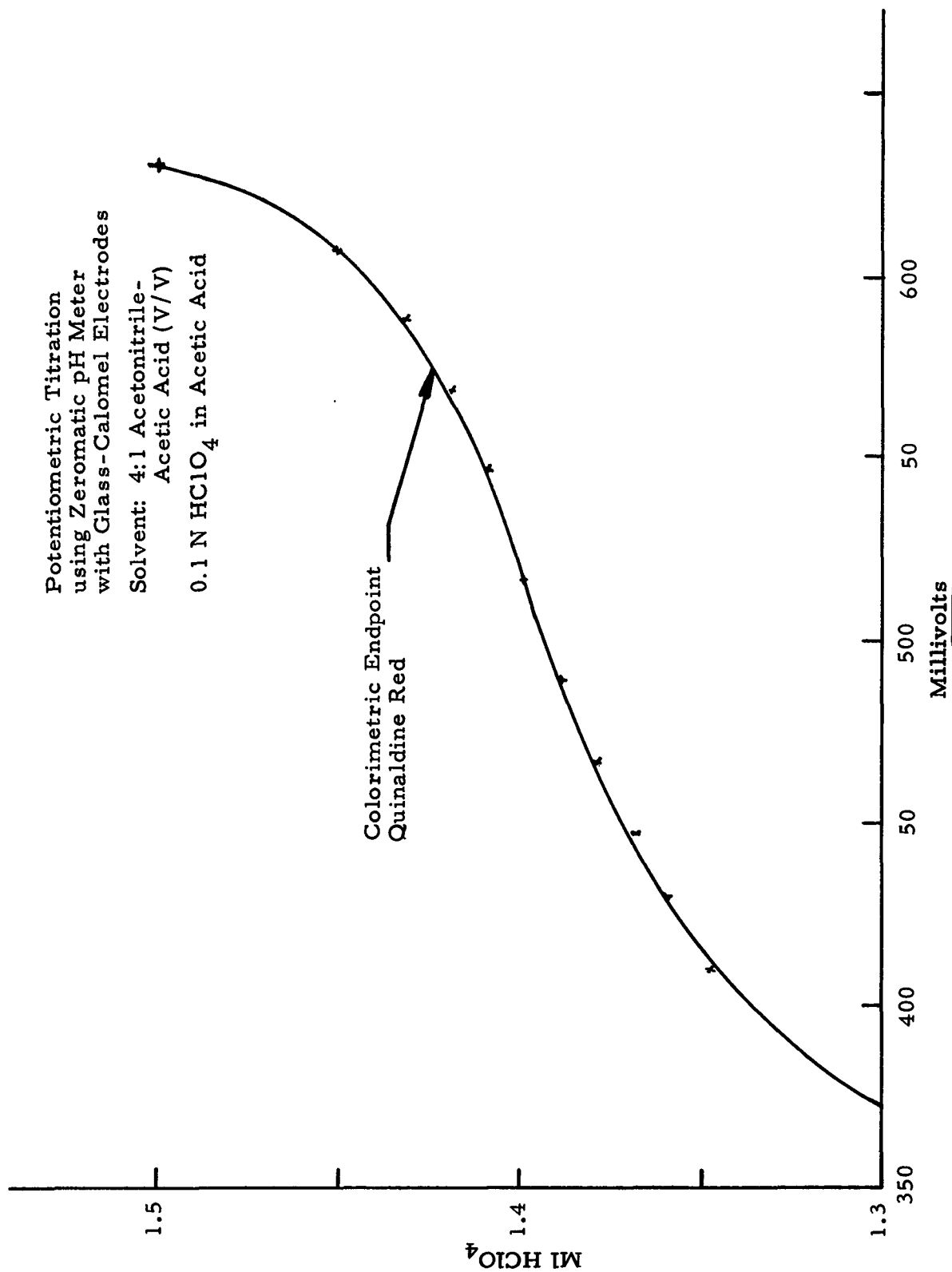


Figure III. Typical Titration Curve for Both Titrations.

BIBLIOGRAPHY

1. Cain, E. F. C. and Stevens, M. R., paper from 2nd International Biannual Symposium on Gas Chromatography (June 1959).
2. Clark, C. C., "Hydrazine," p.43, Copyright (1953) by Mathieson Chemical Corporation.
3. Fritz, J. S. and Fulda, M. O., Anal. Chemistry 25, 1837 (1953).
4. Fritz, J. S. and Schenk, G. H., Anal. Chem. 31, 1808 (1959).
5. Mackey, M. R., unpublished paper, 28 May 1962.
6. Malone, H. E., Anal. Chemistry, 33, 574 (1961).
7. Sadtler, Samuel P. and Son, Inc., Standard Spectra, No. 19955, Aldrich Chemical Co.; 11889, Iowa State College; 11890, Iowa State College.
8. Streuli, C. A., Anal. Chemistry 30, 997 (1958).
9. Weed, A. F., Analytical Method CA-106 by Food Machinery and Chemical Corporation, May (1960).
10. Wimer, D. C., Anal. Chemistry 30, 77 (1958).
11. Wimer, D. C., *ibid*, 30, 2060 (1958).
12. Zarembo, J. D., Analytical Method CA-196 by Food Machinery and Chemical Corporation, February (1962).

DISTRIBUTION

Rocket Research Laboratories
Attn: DGPCS, Mr. Hugh Malone
6593 Test Group (Dev)
Edwards, Calif (50)

Commander
AFSC
Foreign Technology Div
Attn: FTD (TD-E3b)
Wright-Patterson AFB, Ohio (1)

Director, Air University Library
Maxwell AFB, Alabama (2)

Air Force Systems Command
Attn: SCR-2
Andrews AFB
Wash 25, DC (2)

Bureau of Ordnance
Department of the Navy
Wash 25, DC (1)

U. S. Naval Research Laboratory
Attn: Dr. H. W. Carhart
Wash 25, DC (1)

LPIA
Applied Physics Laboratory
Johns Hopkins University
8621 Georgia Ave
Silver Spring, Maryland (10)

Department of the Air Force
Hq USAF, DCS/D
Attn: AFDRT/AS
Wash 25, DC (1)

Bureau of Aeronautics
Department of the Navy
Attn: PP-731
Wash 25, DC (2)

MAAMA (MAOQ/Mr. Alvin Traynor)
Olmsted AFB, Pa. (10)

Office of Naval Research
Attn: Mr. B. Hornstein
Washington 25, DC (3)

Office of Chief of Ordnance
Department of the Army
Attn: Dr. C. M. Hudson (ORDTU)
Wash 25, DC (3)

Commanding General
Army Chemical Center
Attn: Dr. E. H. Krackow
Edgewood, Maryland (2)

Picatinny Arsenal
Attn: Mr. L. H. Rosen
Dr. J. D. Clark
Dover, New Jersey (2)

AFMTC (Mr. Normile)
Patrick AFB, Fla (1)

U. S. Naval Ordnance Test Station
Attn: Mr. E. M. Bens
China Lake, Calif (2)

Commander
Research and Test Squadron
Attn: Mr. Quirmback
Mr. M. D. Weber
Holloman AFB, New Mexico (2)

Naval Powder Factory
Attn: Mr. W. O. Brimijoin
Indianhead, Maryland (2)

New York University
Department of Chemical Engineering
Attn: Dr. P. F. Winternitz
New York 53, NY (1)

Distribution p2 of 3

Westvaco Chlor-Alkali Division
Food Machinery
Attn: Dr. Paul Derr
Dr. Louis Diamond
Princeton, NJ (2)

Commanding General
Engineering Branch
Technical and Engineering Div, O. M. L.
Attn: Mr. D. I. Graham, Jr.
Redstone Arsenal
Huntsville, Alabama (1)

Commanding General
Mechanical Branch
Technical and Engineering Div, O. M. L.
Attn: Mr. P. E. Redding
Redstone Arsenal
Huntsville, Alabama (2)

Aerojet-General Corporation
Attn: Dr. D. L. Armstrong
Mr. E. M. Wilson
Azusa, Calif (2)

Bell Aircraft Corporation
Attn: Dr. S. A. Long
Mr. Henry Heubusch
PO Box 1
Buffalo, NY (3)

Phillips Petroleum Company
Attn: Mr. Homer Fix
Bartlesville, Okla (2)

Commercial Solvents Co.
Research and Development Div
Attn: Dr. R. S. Egly
Dr. J. A. Riddick
Terre Haute, Indiana (2)

Olin-Mathieson Chemical Corporation
Attn: Mr. D. Griffin
New Haven, Conn (2)

Olin-Mathieson Chemical Corporation
Attn: Dr. B. Hill
Baltimore, Maryland (2)

General Chemical Division
Allied Chemical and Dye Corporation
Attn: Mr. T. A. Wallace
40 Rector St
New York, NY (1)

Southern Research Institute
Attn: Dr. W. J. Barrett
917 S Twentieth Ave
Birmingham, Ala (1)

California Institute of Technology
Jet Propulsion Laboratory
2800 Oak Grove Dr
Pasadena, Calif (1)

Commander
Air Proving Grounds, Chemical Section
Climatic Project Air Corps
Attn: Mr. E. G. Taylor
Eglin AFB, Fla (1)

General Electric Company
Rocket Propulsion Units
Attn: Mr. E. St Clair Gantz
Building 300
Cincinnati 15, Ohio (1)

Aerojet-General Corporation
Attn: Dr. W. R. Fish
Sacramento, Calif (1)

Aerojet-General Corporation
Attn: Mr. Harry York
Sacramento, Calif (1)

Rohm and Haas Co.
Redstone Div
Attn: Dr. Keith A. Booman
Huntsville, Ala (1)

Rocketdyne
Attn: Mr. E. F. C. Cain, Code
D/591-356
6633 Canoga Ave
Canoga Park, Calif (1)

Distribution p3 of 3

New York University
College of Engineering, Bldg 3
Attn: Dr. Arrigo A. Carotti
New York 53, NY (1)

Naval Research Lab, Code 6180
Attn: Dr. Robert Hazlett
Wash 25, DC (1)

Thiokol Chemical Corporation
Reaction Motors Div
Attn: Mr. Fred G. Hoffman
Denville, NY (1)

George C. Marshall
Space Flight Center
NASA M-S and M-MA
Attn: Mr. John Nunnelley
Huntsville, Ala (1)

Lewis Research Center
Propellant Chemistry Div
Analytical Chemistry Sec
Attn: Mr. Adolph Spakowski
Cleveland, Ohio (1)

General Chemical Co.
Research Laboratory
Attn: Mr. A. D. Turissini
Morristown, NJ (1)

Army Rocket and Guided Missile Agency
Research Laboratory, OMLD
Attn: Dr. Walter Wharton,
Code ORDXR-ORC
Huntsville, Ala (1)

Allied Chemical Corporation
Nitrogen Div
Attn: Dr. R. E. Ogden
Hopewell, Virginia (1)

Air Products and Chemicals Inc
Attn: Mr. Willard Ent
Allentown, Pa. (1)

Texaco Research Center
Attn: Dr. H. E. Vermillion
Beacon, NY (1)
AFSC (STLO)
Lewis Research Center (NASA)
Cleveland 35, Ohio (1)

United Technology Corporation
Attn: Mr. Paul Marti
PO Box 358
Sunnyvale, Calif (1)

The Martin Company (C-124)
Research Dept
Attn: Dr. James A. Bowman
PO Box 1176
Denver 1, Colo (1)

Stanford Research Institute
Poulter Laboratories
Attn: Dr. Eugene Burns
Menlo Park, Calif (1)

Office of Naval Research
Department of the Navy
Attn: Mr. Thomas Dowd
Wash 25, DC (1)

AFLC (MCRQD, Mr. Russell Dye)
Wright-Patterson AFB, Ohio (1)

AFSC (MCQDT/Mr. Stephen Dzuroff)
Wright-Patterson AFB, Ohio (1)

Armed Service Tech Information
Agency
Arlington Hall Station
Attn: TIPCR
Arlington 12, Virginia (10)

Air Force Flight Test Center
Attn: FTOOT
Edwards AFB, Calif (10)

Rocket Research Laboratories
6593d Test Group (Dev), (DGLP)
Edwards, Calif (10)

DCAS (DCLMT/Mrs. Densmore)
AF Unit Post Office
Los Angeles 45, Calif (1)

BSD (BSREQ, Dr. Paul Leatherman)
Norton AFB, Calif (1)

Robert A. Biggers (DGPCS)
Rocket Research Laboratories
6593 Test Group (Dev)
Edwards, Calif (1)

Rocket Research Laboratories, Edwards, Calif. Rpt No. SSD-TDR-62-203. AN ACETYLATION METHOD FOR DETERMINING MIXTURES OF EITHER HYDRAZINE-1, 1-DIMETHYLHYDRAZINE OR MONOMETHYL-HYDRAZINE-1, 1-DIMETHYLHYDRAZINE. 14 p. incl illus. Dec 62. Unclassified Report A method is presented for determining mixtures of either hydrazine-1, 1-dimethylhydrazine ($N_2H_4/UDMH$) or monomethylhydrazine-1, 1-dimethylhydrazine ($MMH/UDMH$). This method is based on the difference in reaction rates for acetylating hydrazines. Both N_2H_4 and MMH react immediately while the acetylation of $UDMH$ proceeds slowly. The reactions are conducted in	1. Acetylation I. AFSC Project No. 3148, Task No. 30197 II. Malone, Hugh E., Biggers, R. A. III. In ASTIA collection
--	--

Rocket Research Laboratories, Edwards, Calif. Rpt No. SSD-TDR-62-203. AN ACETYLATION METHOD FOR DETERMINING MIXTURES OF EITHER HYDRAZINE-1, 1-DIMETHYLHYDRAZINE OR MONOMETHYL-HYDRAZINE-1, 1-DIMETHYLHYDRAZINE. 14 p. incl illus. Dec 62. Unclassified Report A method is presented for determining mixtures of either hydrazine-1, 1-dimethylhydrazine ($N_2H_4/UDMH$) or monomethylhydrazine-1, 1-dimethylhydrazine ($MMH/UDMH$). This method is based on the difference in reaction rates for acetylating hydrazines. Both N_2H_4 and MMH react immediately while the acetylation of $UDMH$ proceeds slowly. The reactions are conducted in	1. Acetylation I. AFSC Project No. 3148, Task No. 30197 II. Malone, Hugh E., Biggers, R. A. III. In ASTIA collection
--	--

acetonitrile-acetic acid medium. Two titrations are required; one measures the total basicity of the hydrazine mixtures; the other titration measures $UDMH$ after the N_2H_4 or MMH has reacted with Ac_2O.	
---	--

acetonitrile-acetic acid medium. Two titrations are required; one measures the total basicity of the hydrazine mixtures; the other titration measures $UDMH$ after the N_2H_4 or MMH has reacted with Ac_2O.	
---	--