

AD-A042 200

ARMY MISSILE RESEARCH AND DEVELOPMENT COMMAND REDSTO--ETC F/6 7/4
SPECTROSCOPY OF LASER IRRADIATED MIXTURES OF BORON TRICHLORIDE --ETC(U)
APR 77 G A TANTON, H C MEYER, R I GREENBERG
DRDMI-TR-77-7

UNCLASSIFIED

NL

| OF |
AD
A042200



END

DATE
FILMED
8-77

ADA 042200



**U.S. ARMY
MISSILE
RESEARCH
AND
DEVELOPMENT
COMMAND**



Redstone Arsenal, Alabama 35809

DDC FILE COPY

12

TECHNICAL REPORT TR-77-7

SPECTROSCOPY OF LASER IRRADIATED
MIXTURES OF BORON TRICHLORIDE AND
HYDROGEN

Physical Sciences Directorate
Technology Laboratory

8 April 1977

Approved for public release; distribution unlimited.



DISPOSITION INSTRUCTIONS

**DESTROY THIS REPORT WHEN IT IS NO LONGER NEEDED. DO NOT
RETURN IT TO THE ORIGINATOR.**

DISCLAIMER

**THE FINDINGS IN THIS REPORT ARE NOT TO BE CONSTRUED AS AN
OFFICIAL DEPARTMENT OF THE ARMY POSITION UNLESS SO DESIG-
NATED BY OTHER AUTHORIZED DOCUMENTS.**

TRADE NAMES

**USE OF TRADE NAMES OR MANUFACTURERS IN THIS REPORT DOES
NOT CONSTITUTE AN OFFICIAL INDORSEMENT OR APPROVAL OF
THE USE OF SUCH COMMERCIAL HARDWARE OR SOFTWARE.**


UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE (When Data Entered)

| REPORT DOCUMENTATION PAGE | | READ INSTRUCTIONS BEFORE COMPLETING FORM |
|--|-----------------------------|---|
| 1. REPORT NUMBER TR-77-7 | 2. GOVT ACCESSION NO. 14 | 3. RECIPIENT'S CATALOG NUMBER DRDMI-TR-77-7 |
| 4. TITLE (and Subtitle) SPECTROSCOPY OF LASER IRRADIATED MIXTURES OF BORON TRICHLORIDE AND HYDROGEN | | 5. TYPE OF REPORT & PERIOD COVERED Technical Report |
| 7. AUTHOR(s) G. A. Tanton, H. C. Meyer, R. I. Greenberg | | 6. PERFORMING ORG. REPORT NUMBER |
| 9. PERFORMING ORGANIZATION NAME AND ADDRESS Commander US Army Missile Research and Development Command Attn: DRDMI-TR Redstone Arsenal, Alabama 35809 | | 8. CONTRACT OR GRANT NUMBER(s) |
| 11. CONTROLLING OFFICE NAME AND ADDRESS Commander US Army Missile Research and Development Command Attn: DRDMI-TI Redstone Arsenal, Alabama 35809 | | 10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS (DA) 1W362303A214 AMCMS Code 632303.2140411 |
| 14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office) | | 12. REPORT DATE 8 Apr 77 |
| 16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited. | | 13. NUMBER OF PAGES 13 |
| 17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report) | | 15. SECURITY CLASS. (of this report) Unclassified |
| 18. SUPPLEMENTARY NOTES | | 15a. DECLASSIFICATION/DOWNGRADING SCHEDULE |
| 19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Laser induced chemistry Fourier transform IR spectrometer Low level impurities Reaction products | | |
| 20. ABSTRACT (Continue on reverse side if necessary and identify by block number) Preparation of missile propellants is one of the areas expected to be significantly impacted by laser induced chemistry. The synthesis by laser induced chemistry of specific boron compounds used as ingredients in propellants has been the purpose of our recent investigations. | | |


410211

mt

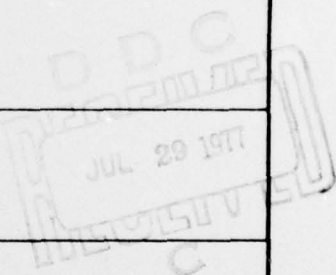

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

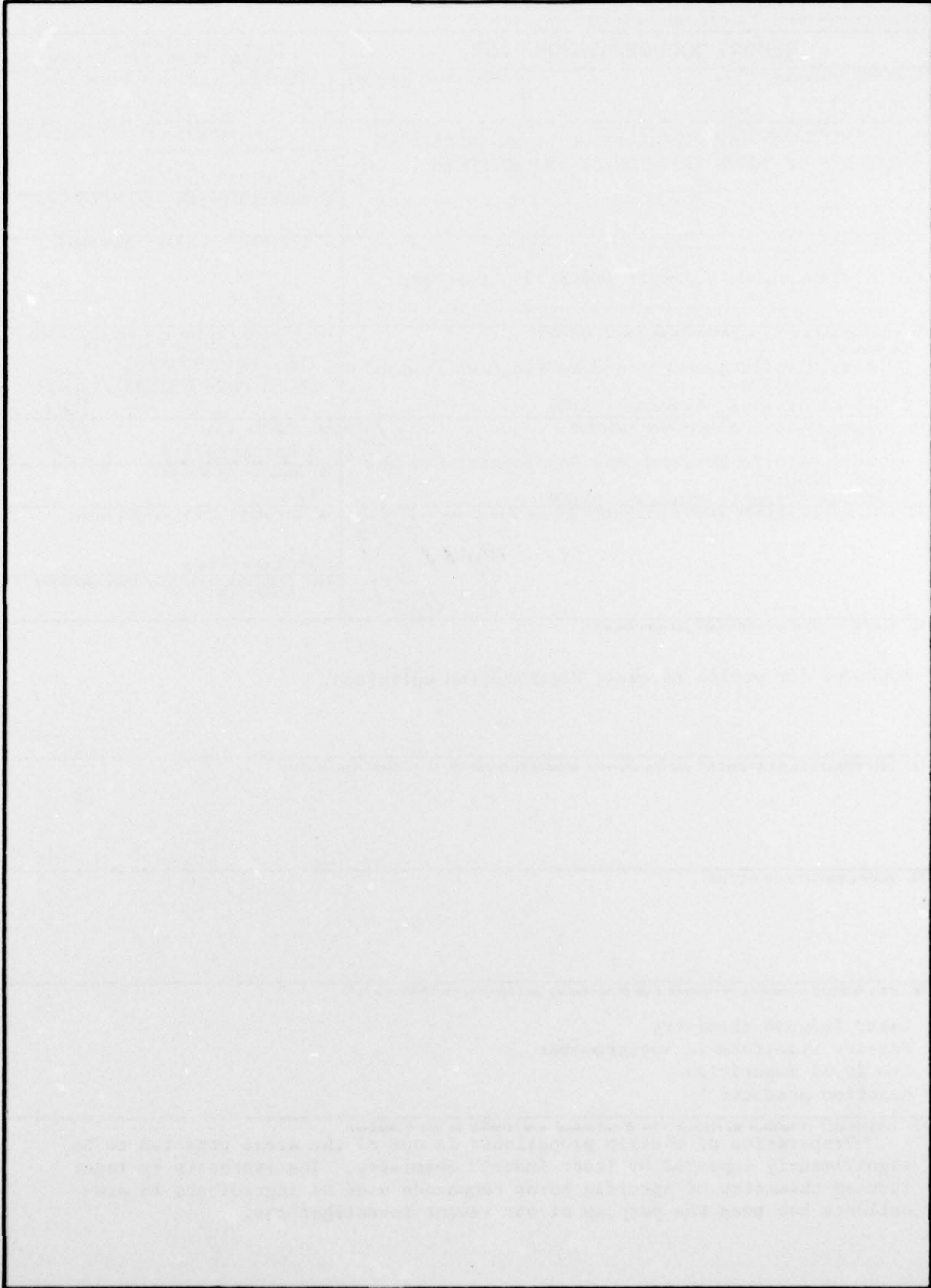

UNCLASSIFIED

SECURITY CLASSIFICATION OF THIS PAGE(When Data Entered)

CONTENTS

| | Page |
|---|------|
| I. INTRODUCTION | 3 |
| II. EXPERIMENTAL | 3 |
| III. TECHNIQUE FOR IDENTIFYING LUMINESCENCE SPECTRA | 3 |
| IV. RESULTS. | 5 |
| V. CONCLUSIONS. | 8 |
| REFERENCES | 10 |

| | |
|---------------------------------|---|
| ACCESSION FOR | |
| RTTS | White Section <input checked="" type="checkbox"/> |
| DIS | Buff Section <input type="checkbox"/> |
| UNANNOUNCED | <input type="checkbox"/> |
| JUSTIFICATION | |
| BY | |
| DISTRIBUTION/AVAILABILITY CODES | |
| Dist. | AVAIL. and/or SPECIAL |
| A | |

I. INTRODUCTION

Preparation of missile propellants is one of the areas expected to be significantly impacted by laser induced chemistry. The synthesis by laser induced chemistry of specific boron compounds used as ingredients in propellants has been the purpose of the recent investigations.

Initially mixtures of boron trichloride (BCl_3) and hydrogen have been irradiated by a pulsed CO_2 laser and the gas was observed [1,2] to emit a pale yellow luminescence. We have obtained the spectrum of this visible luminescence and identified spectral lines. Use of a Fourier transform IR spectrometer (FTS) enabled us to detect low level impurities in both the initial mixture of BCl_3 and H_2 and the reaction products. These IR spectra showed the presence of the impurity phosgene (COCl_2) in the initial gas mixture and CO in the reaction products in addition to the expected HBCl_2 and HCl .

II. EXPERIMENTAL

Typically, a 10-cm brass or stainless steel cell fitted with 5 cm diameter NaCl or KBr windows was used. Three of those cells can be seen in the upper middle part of Figure 1. It was passivated for several hours with BCl_3 then evacuated and refilled with 25 to 100 Torr of BCl_3 and 25 to 400 Torr of H_2 in a desired ratio. This mixture was irradiated with 200 ns wide 0.7 J pulses from a Lumonics Model TEA 203 tunable multigas laser operating with a 2-second repetition rate at the P20 CO_2 (10.591 μm) wavelength. During irradiation, visible luminescence was detected through a window in the cell perpendicular to the beam direction and recorded on Polaroid type 107 film after passing through a prism spectrograph. Absorption spectra were determined with a Digilab Model 20-B Fourier transform spectrometer before and after irradiation. This instrument is shown in Figure 2.

III. TECHNIQUE FOR IDENTIFYING LUMINESCENCE SPECTRA

A quick and accurate method for identifying optical spectra has been developed. In this technique, positions corresponding to given wavelengths are measured from an arbitrarily chosen reference line on a spectrographic plate of known wavelengths taken with our prism spectrograph. These distances and wavelengths are used to calculate the coefficients of a polynomial using a least-squares fit. We use a basic program for polynomial regression of any desired degree up to 9th degree written for the 9830A Hewlett-Packard programmable calculator for this fit. This polynomial is then used to plot a curve of linear distance versus wavelength. This step completes the calibration.



Figure 1. Gas cells for laser irradiation.

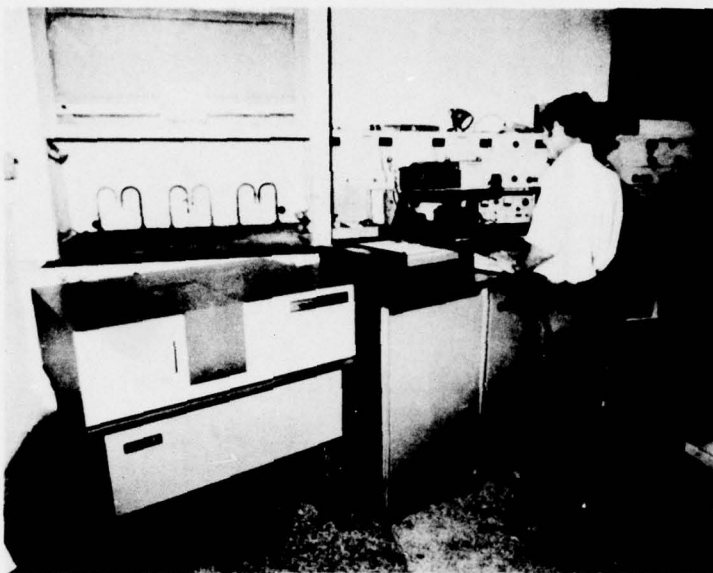


Figure 2. Fourier transform spectrometer.

The calculator will then plot, in our particular spectrographic plate format, the position of any wavelength programmed into it. These calculated spectra are plotted as a column of lines of given equal horizontal length and vertical separation calculated from the polynomial obtained from the calibration points.

Identification of molecular and atomic spectra from a spectrographic plate of an unknown is made by overlapping it with calculator drawn spectra plotted on a transparency. Speed and accuracy are increased over conventional methods since several spectral lines of a series can be simultaneously compared, thereby eliminating the usual tedious method of measuring the position of each line individually. Deviations between data and standard spectral line positions, usually due to the geometry of a particular experimental arrangement or to the mode of excitation, can be evaluated immediately by inspection. In this manner, confusion arising between measured and standard spectra is reduced. Relative intensity in addition to wavelength can be indicated on the calculator plot. This is illustrated in the flame spectrum of boron (Figure 3). The principal limitation of this method is the availability of known spectra.

IV. RESULTS

The absorption spectrum after a 300-pulse irradiation of a 50-Torr mixture of boron trichloride and hydrogen in a 1:1 ratio is shown in Figure 4. An increase in HCl absorption band after irradiation can be seen at 2900 cm^{-1} , accompanied by a decrease of BCl_3 absorption


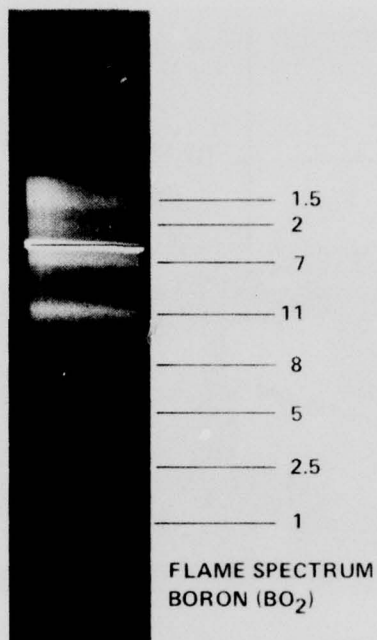


Figure 3. Flame spectrum of boron compared with computed spectrum.

at 1400, 995, and 955 cm^{-1} . Creation of bands at 2620, 1090, and 900 cm^{-1} shows that HBCl_2 was formed as an end product of laser induced reactions according to the reaction

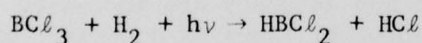

as previously reported by Rockwood and Hudson [1]. Subsequently the optimum starting ratio of BCl_3/H_2 was found to be 1:4 for HBCl_2 formation. The optical emission spectra of BCl_3 irradiated with a pulsed CO_2 laser was found to have a continuous spectral distribution [3,4] from 5300 Å to > 6300 Å. When a mixture of BCl_3 with an excess of H_2 was used emission band spectra which fit the $\text{B}^1\Sigma \rightarrow \text{A}^1\Pi$ and $\text{d}^3\pi \rightarrow \text{a}^3\pi$ series of CO were observed. Bands due to Cl^+ and Cl_2^+ overlapped by the two most prominent bands of the $^3\pi \rightarrow ^3\pi$ series of C_2 were detected as shown in Figure 5. COCl_2 , a main impurity in the c.p. grade BCl_3 used, was identified as the probable source of CO. There were also bands just above 5000 Å which correspond to BO bands previously observed [5].

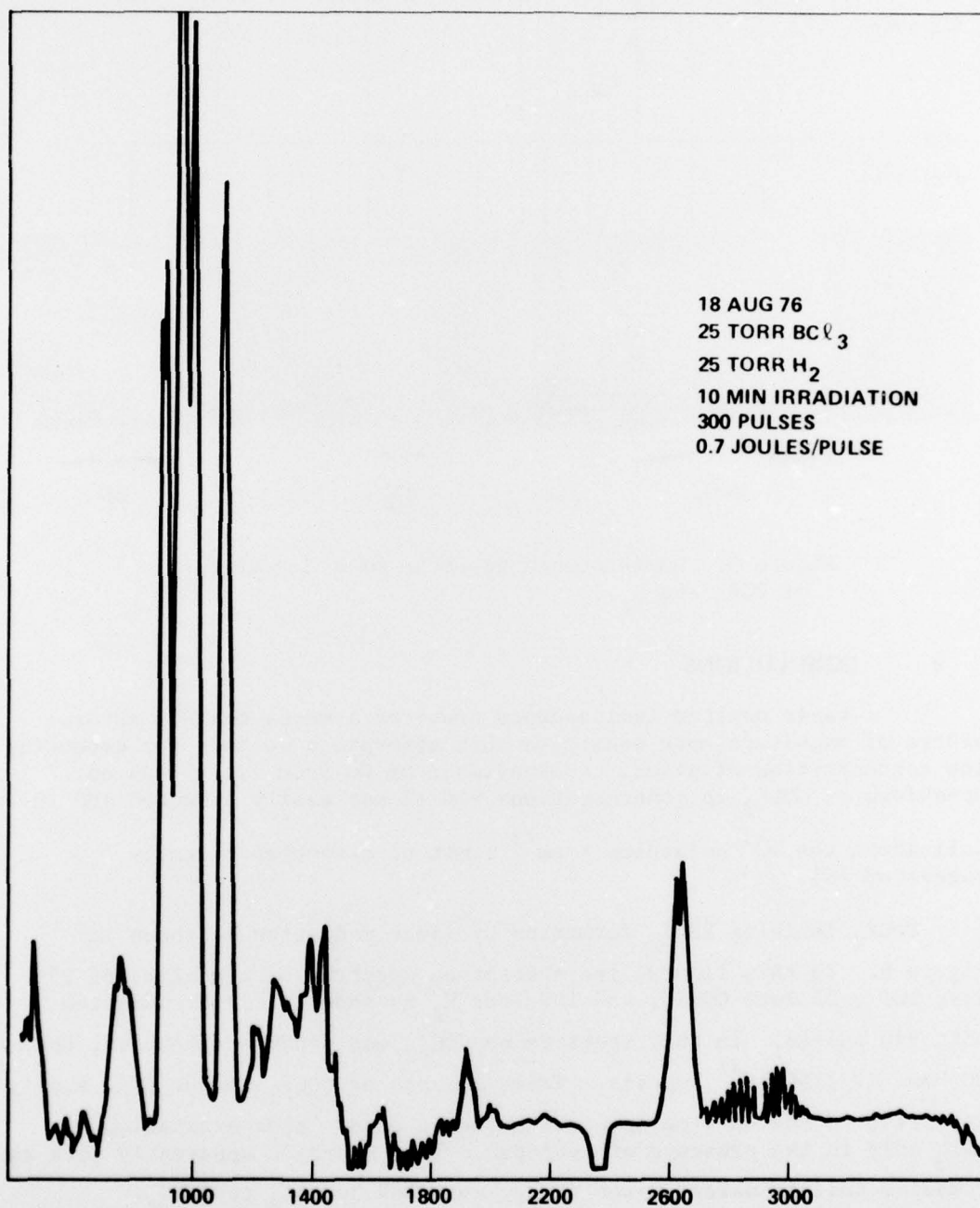


Figure 4. Absorption spectrum of a 50 Torr mixture of BCl_3 and H_2 in a 1:1 ratio, after irradiation with 300 pulses.

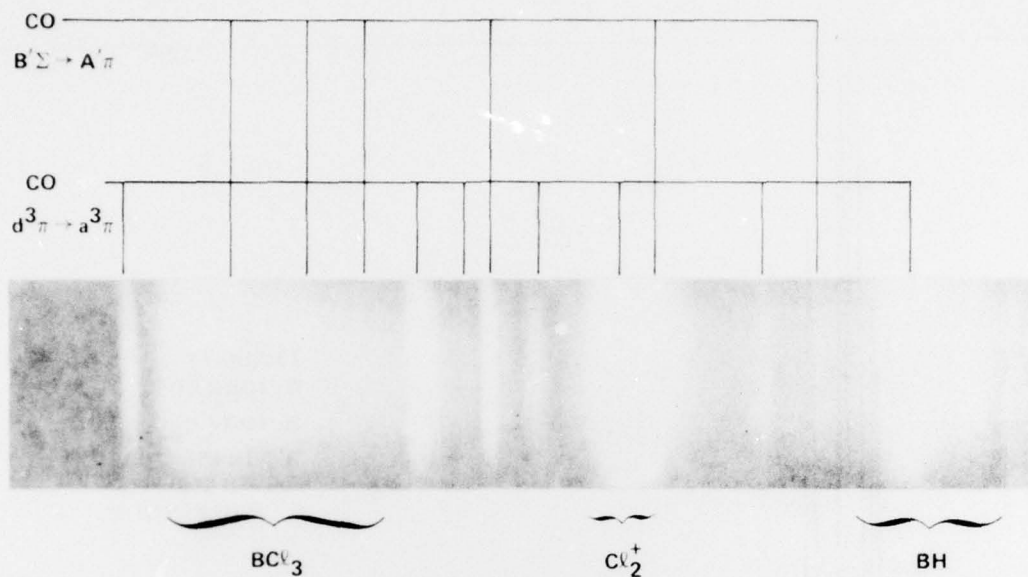


Figure 5. Luminescence spectrum of a 1:4 mixture of BCl_3 and H_2 .

V. CONCLUSIONS

Laser excited luminescence provides a means that is several orders of magnitude more sensitive than absorption methods for detecting low concentration of gases. Luminescence of CO from laser induced breakdown of COCl_2 in concentrations $\sim 0.1\%$ was easily detected and is well above the 10^5 molecules \cdot cm^{-3} limit of detection recently suggested [6].

COCl_2 inhibits HBCl_2 formation by laser radiation as shown in Figure 6. In this figure, the absorption spectrum of a mixture of 25 Torr BCl_3 , 25 Torr COCl_2 , and 100 Torr H_2 is shown after irradiation with 810 pulses. In this spectrum no HBCl_2 was produced; however, the CO band at 2150 cm^{-1} appears. Trace amounts of COCl_2 in our BCl_3 supply gave rise to the CO line spectrum observed under laser excitation of BCl_3 only in the presence of hydrogen. Thus hydrogen apparently acts as a medium that transfers laser energy absorbed by BCl_3 to COCl_2 .

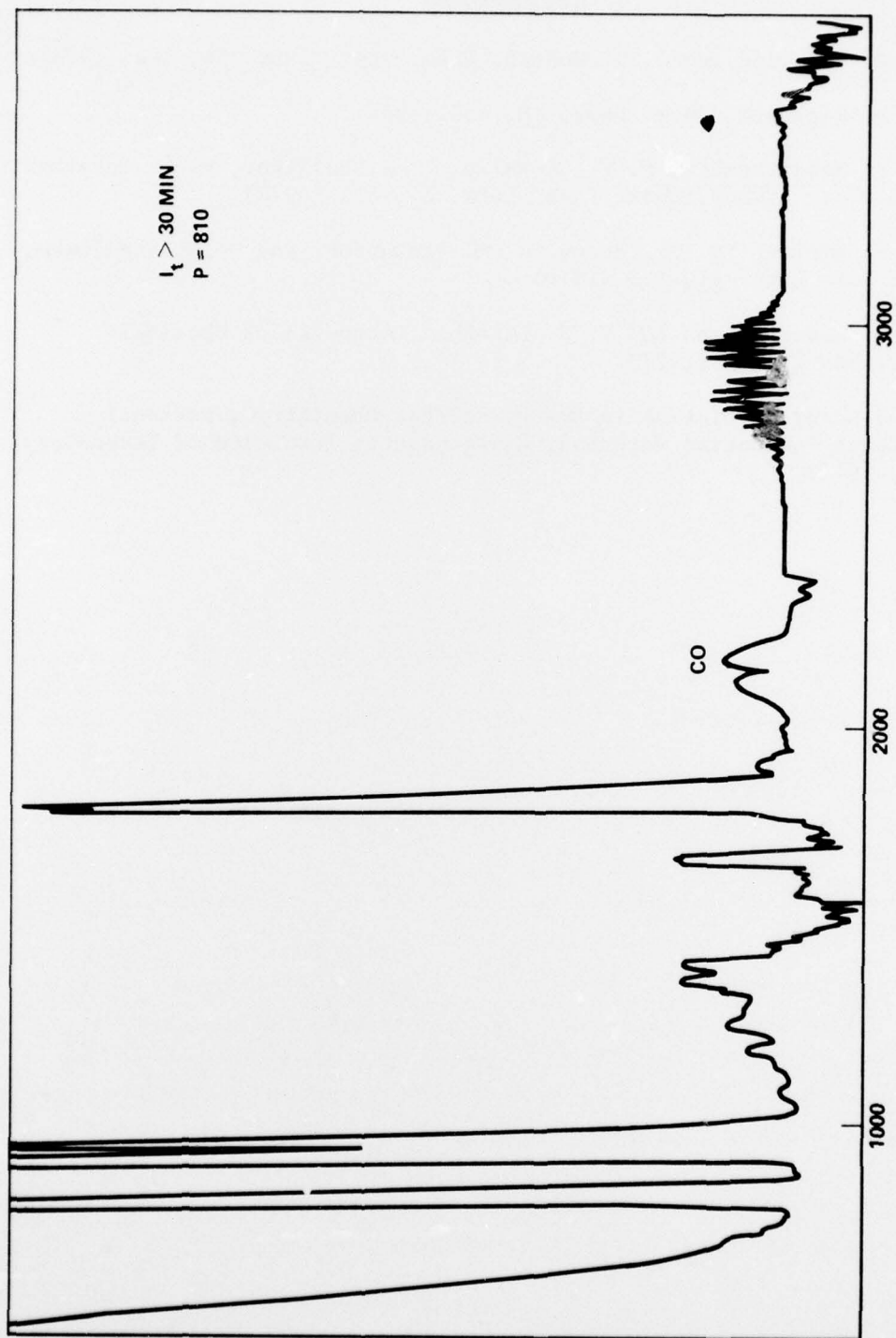


Figure 6. Absorption spectrum of a mixture of 25 Torr BCl_3 , 25 Torr COCl_2 , and 100 Torr H_2 after irradiation with 810 pulses.

REFERENCES

1. S. D. Rockwood and J. W. Hudson, Chem. Phys. Lett. 34, 542 (1975).
2. S. D. Rockwood, Chem. Phys. 10, 453 (1974).
3. R. V. Ambartzumian, N. V. Chekalin, V. S. Doljikov, V. S. Tetokhov, and E. A. Ryabov, Chem. Phys. Lett. 25, 515 (1974).
4. N. V. Karlov, Yu. N. Petrov, A. M. Prokhorov, and O. M. Stel'makh, J.E.T.P. Lett. 11, 135 (1970).
5. R. V. Ambartzumian and V. S. Letokhov, Accounts of Chemical Research 10, 61 (1977).
6. "The Laser Revolution in Energy-Related Chemistry", National Science Foundation Workshop, Massachusetts Institute of Technology, May 1976.

DISTRIBUTION

| | No. of Copies | | No. of Copies |
|--|---------------|---|-----------------------|
| Commander Defense Documentation Center Attn: DDC-TGA Cameron Station Alexandria, Virginia 22314 | 12 | Commander US Army Combined Arms Combat Development Activity Fort Leavenworth, Kansas 66027 | 1 |
| Commander US Army Research Office Attn: Dr. R Lontz P.O. Box 12211 Research Triangle Park, North Carolina 27709 | 2 | Commander US Army Frankford Arsenal Philadelphia, Pennsylvania 19137 | 1 |
| US Army Research and Standardization Group (Europe) Attn: DRXSN-E-RX, Dr. Alfred K. Nedoluha Box 65 FPO, New York 09510 | 2 | Commander US Army Picatinny Arsenal Dover, New Jersey 07801 | 1 |
| US Army Materiel Development and Readiness Command Attn: Dr. Gordon Bushy Dr. James Bender Dr. Edward Sedlak 5001 Eisenhower Avenue Alexandria, Virginia 22333 | 1 1 1 | Commander US Army Tank Automotive Development Command Attn: DRDTA-RWL Warren, Michigan 48090 | 1 |
| Headquarters Hq DA (DAMA-ARZ) Washington, D. C. 20310 | 2 | Commander US Army Mobility Equipment Research and Development Command Fort Belvoir, Virginia 22060 | 1 |
| Director of Defense Research and Engineering Engineering Technology Attn: Mr. L. Weisberg Washington, D. C. 20301 | 2 | Commander US Army Harry Diamond Laboratories Attn: Dr. Stan Kulpa 2800 Powder Mill Road Adelphi, Maryland 20783 | 1 |
| Director Defense Advanced Research Project Agency 1400 Wilson Boulevard Arlington, Virginia 22209 | 1 | Commander US Army Armament Command Rock Island, Illinois 61202 | 1 |
| Commander US Army Aviation Systems Command 12th and Spruce Streets St. Louis, Missouri 63166 | 1 | Commander US Army Foreign Science and Technology Center Federal Office Building 220 7th Street, NE Charlottesville, Virginia 22901 | 1 |
| Director US Army Air Mobility Research and Development Laboratory Ames Research Center Moffett Field, California 94035 | 1 | Commander US Army Training and Doctrine Command Fort Monroe, Virginia 22351 | 1 |
| Commander US Army Electronics Command Attn: DRSEL-TL-T, Dr. Jacobs DRSEL-CT, Dr. R. Buser Fort Monmouth, New Jersey 07703 | 1 1 | Director Ballistic Missile Defense Advanced Technology Center Attn: ATC-D ATC-O ATC-R ATC-T P.O. Box 1500 Huntsville, Alabama 35808 | 1 1 1 1 1 |
| Director US Army Night Vision Laboratory Attn: John Johnson Fort Belvoir, Virginia 22060 | 1 | Commander US Naval Air Systems Command Washington, D. C. 20360 | 1 |
| Director Atmospheric Sciences Laboratory US Army Electronics Command White Sands Missile Range, New Mexico 88002 | 1 | Chief of Naval Research Department of the Navy Washington, D. C. 20360 | 1 |
| Director US Army Ballistic Research Laboratories Attn: Ken Richer Aberdeen Proving Ground, Maryland 31005 | 1 | Commander US Naval Air Development Center Attn: Radar Division Warminster, Pennsylvania 18974 | 1 |
| | | Commander US Naval Electronics Lab Center San Diego, California 92152 | 1 |
| | | Commander US Naval Surface Weapons Center Dahlgren, Virginia 22448 | 1 |

| | No. of Copies | | No. of Copies |
|--|------------------|--|---------------|
| Commander US Naval Weapons Center Attn: Mr. Robert Moore China Lake, California 93555 | 1 | Director Calspan Corporation Attn: R. Kell P.O. Box 235 Buffalo, New York 14221 | 1 |
| Director Naval Research Laboratory Attn: Code 5300, Radar Division Dr. Skolnik Code 5370, Radar Geophysics Branch Code 6460, Electromagnetic Prop Branch Washington, D. C. 20390 | 1 1 1 1 | California Institute of Technology Attn: Dr. N. George 1201 E. California Boulevard Pasadena, California 91109 | 1 |
| Commander Rome Air Development Center US Air Force Attn: R. McMillan, OCSA James Wasielewski, TRRC Griffiss Air Force Base, New York 13440 | 1 1 | Optical Science Consultants Attn: Dr. D. L. Fried P.O. Box 388 Yorba Linda, California 92686 | 1 |
| Commander US Air Force, AFOSR/NP Attn: LT COL Gordon Wepfer Bolling Air Force Base Washington, D. C. 20332 | 1 | Dr. D. T. Hodges, Jr. The Aerospace Corporation 2350 East El Segundo Boulevard Los Angeles, California 90009 | 1 |
| Commander US Air Force Avionics Laboratory Attn: D. Rees CPT James D. Pryce, AFAL/WE Dr. B. L. Sowers, AFAL/RWI Wright-Patterson Air Force Base, Ohio 45433 | 1 1 1 | Commander Center for Naval Analyses Attn: Document Control 1401 Wilson Boulevard Arlington, Virginia 22209 | 1 |
| Commander (AFEL) Hanscom Air Force Base, Maryland 01731 | 1 | Hughes Aircraft Company Missile System Group Attn: Dr. J. Glassman Canoga Park, California 91304 | 1 |
| Commander AFATL/LMF Eglin Air Force Base, Florida 32544 | 1 | Goodyear Aerospace Corporation Arizona Division Attn: Mr. Fred Wilcox Mr. P. W. Murray Litchfield Park, Arizona 85340 | 1 1 |
| Hughes Research Laboratory Attn: Mr. Smith John F. Heney J. M. Baird 3001 Malibu Canyon Road Malibu, California 90265 | 1 1 1 | Honeywell, Inc. Systems and Research Division Attn: Dr. Paul Kruse Minneapolis, Minnesota 55413 | 1 |
| Environmental Research Institute of Michigan Radar and Optics Division Attn: Dr. A. Kozma Dr. C. C. Aleksoff P.O. Box 618 Ann Arbor, Michigan 41807 | 1 1 | Raytheon Company Attn: A. V. Jelalian 528 Boston Post Road Sudbury, Massachusetts 01776 | 1 |
| Dr. J. G. Castle Physics Department University of Alabama 4701 University Drive, N.W. Huntsville, Alabama 35807 | 1 | Georgia Institute of Technology Engineering Experiment Station Attn: James Gallaper 347 Ferst Drive Atlanta, Georgia 30332 | 1 |
| Dr. Clyde Riley Department of Chemistry University of Alabama 4701 University Drive, N.W. Huntsville, Alabama 35807 | 1 | The Rand Corporation Attn: Dr. S. J. Dudzinsky, Jr. 1700 Main Street Santa Monica, California 90406 | 2 |
| Science and Technology Division Institute of Defense Analysis Attn: Dr. Vincent J. Corcoran 400 Army-Navy Drive Arlington, Virginia 22202 | 1 | Superior Technical Services 4308 Governors Drive Attn: Ray King Huntsville, Alabama 35805 | 1 |
| | | DRDMI-LP, Mr. Voigt | 1 |
| | | DRDMI-X, Dr. McDaniel | 1 |
| | | -T, Dr. Kobler | 1 |
| | | -H, Dr. Rose | 1 |
| | | -HS, Dr. Roberts | 1 |
| | | -TD | 1 |
| | | -TE | 1 |
| | | -TE, Mr. Pittmann | 1 |
| | | -TC | 1 |
| | | -TK | 1 |

| | No. of Copies |
|-----------------------|---------------|
| DRDMI-TR, Dr. Hartman | 1 |
| Dr. Bennett | 1 |
| Mrs. Davis | 5 |
| -TRD | 50 |
| -TBD | 3 |
| -TI (Record Set) | 1 |
| (Reference Copy) | 1 |