

AD-A190 160

DETERMINATION OF AQUO-COMPLEXES OF CHROMIUM AND IRON IN 1/1
CHROMIUM PLATING A. (U) ARMY ARMAMENT RESEARCH
DEVELOPMENT AND ENGINEERING CENTER WTT. S SOPOK


UNCLASSIFIED

NOV 87 ARCCD-TR-87034

F/8 7/4

ML



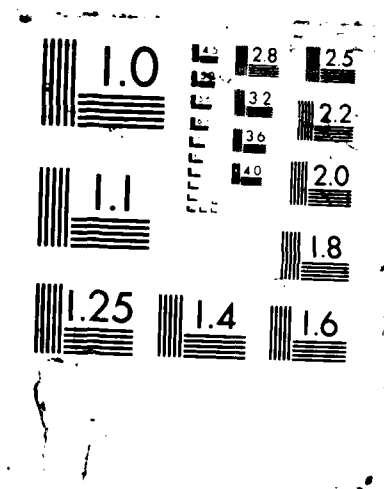

DTIC FILE COPY

AD

TECHNICAL REPORT ARCCB-TR-87034

**DETERMINATION OF AQUO-COMPLEXES OF
CHROMIUM AND IRON IN CHROMIUM PLATING AND
POLISHING SOLUTIONS BY ION CHROMATOGRAPHY**

SAMUEL SOPOK

DTIC
ELECTE
JAN 05 1988
S D

NOVEMBER 1987

AD-A190 168



**US ARMY ARMAMENT RESEARCH, DEVELOPMENT
AND ENGINEERING CENTER**

**CLOSE COMBAT ARMAMENTS CENTER
BENÉT WEAPONS LABORATORY
WATERVLIET, N.Y. 12189-4050**

APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED

87 12 28 032

DISCLAIMER

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorized documents.

The use of trade name(s) and/or manufacturer(s) does not constitute an official indorsement or approval.

DESTRUCTION NOTICE

For classified documents, follow the procedures in DoD 5200.22-M, Industrial Security Manual, Section II-19 or DoD 5200.1-R, Information Security Program Regulation, Chapter IX.

For unclassified, limited documents, destroy by any method that will prevent disclosure of contents or reconstruction of the document.

For unclassified, unlimited documents, destroy when the report is no longer needed. Do not return it to the originator.

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER ARCCB-TR-87034	2. GOVT ACCESSION NO. <i>AD-A190168</i>	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) DETERMINATION OF AQUO-COMPLEXES OF CHROMIUM AND IRON IN CHROMIUM PLATING AND POLISHING SOLUTIONS BY ION CHROMATOGRAPHY		5. TYPE OF REPORT & PERIOD COVERED Final
7. AUTHOR(s) Samuel Sopok		6. PERFORMING ORG. REPORT NUMBER
9. PERFORMING ORGANIZATION NAME AND ADDRESS US Army ARDEC Benet Laboratories, SMCAR-CCB-TL Watervliet, NY 12189-4050		8. CONTRACT OR GRANT NUMBER(s)
11. CONTROLLING OFFICE NAME AND ADDRESS US Army ARDEC Close Combat Armaments Center Picatinny Arsenal, NJ 07806-5000		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS AMCMS No. 6126231BL PRON No. M77F000901M71A
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		12. REPORT DATE November 1987
		13. NUMBER OF PAGES 7
		15. SECURITY CLASS. (of this report) UNCLASSIFIED
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES Submitted to <u>LC-GC Magazine of Liquid and Gas Chromatography.</u>		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Aquo-complexes Chromium Iron Chromium Plating Solutions Polishing Solutions Ion Chromatography		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) Aquo-complexes of chromium and iron in chromium plating and polishing solutions of individual concentration levels above 7.5 g/l result in a poor quality plate condition. An ion chromatographic procedure utilizing UV-visible and atomic absorption detection for determining 0-20 g/l concentration levels of these aquo-complexes is described. Since no derivatizations are necessary, matrix effects are much less common compared to previous chromatographic methods. The relative standard deviation is about two percent and sensitivity is about 1 ppm for these determinations.		

TABLE OF CONTENTS

	<u>Page</u>
INTRODUCTION	1
EXPERIMENTAL PROCEDURES	1
RESULTS AND DISCUSSION	2
REFERENCES	6

TABLES

I. UV-VISIBLE SPECTRA OF DETECTABLE COMPLEXES	3
II. DETERMINATION OF Cr(III) IONS IN CHROMIUM PLATING SOLUTIONS	4
III. DETERMINATION OF Fe(III) IONS IN CHROMIUM PLATING SOLUTIONS	4
IV. DETERMINATION OF Fe(II) IONS IN CHROMIUM POLISHING SOLUTIONS	5

Accession For	
NTIS GRA&I	☒
DTIC TAB	☐
Unannounced	☐
Justification	
By	
Distribution/	
Availability Codes	
Avail and/or	
Dist	Special
A-1	



INTRODUCTION

The separation and subsequent detection of aquo-complexes of chromium and iron in chromium plating and polishing solutions by ion exchange chromatography are necessary to maximize the plating quality of low alloy steels.

These aquo-complexes should not individually exceed 7.5 g/l, but may range from 0-20 g/l. The ions result from metal finishing processes such as polishing, reverse-plating, and normal plating providing Fe(II), Fe(III), and Cr(III) ions, respectively. Additional Cr(III) ions result from accidental oxidation of ethylene glycol due to coolant leaks.

Past methods (refs 1-7) have derivatized these aquo-complexes to enhance separation and detection. Such methods have included high cost-element specific atomic absorption detection (ref 1) and more practical low cost-multi-element UV-visible detection (refs 2-7).

The chemistries of these derivatizations are difficult to control when past methods are applied to chromium plating and polishing solutions. Acid strength and variability affect the derivatization chemistries causing matrix effects. The purpose of this study is to develop an ion chromatographical procedure that eliminates the derivatization chemistries, while providing an analysis range of 0-20 g/l for these aquo-complexes in chromium plating and polishing solutions.

EXPERIMENTAL PROCEDURE

The Dionex model 2020i ion chromatograph is the analysis system used. The system flowstream consists of the element reservoir (0.283 M HNO_3), element delivery system (2.3 ml/min), injector valve (500 μl loop), separator column

References are listed at the end of this report.

(Dionex HPIC-CG2), UV-visible detector (Dionex, 0.5 AUFS), and the atomic absorption detector (Perkin-Elmer, aspirator flow rate is 2.8 ml/min). Air-tight coupling of the two detectors is easily accomplished with the standard tubing provided. For aquo-complexes of iron and chromium, UV-visible detection is at 330 and 570 nm and atomic absorption detection is at 248.3 and 357.9 nm, respectively.

UV-visible wavelengths are chosen from initial spectra (300-600 nm) acquired for these aquo-complexes on a UV-visible spectrophotometer (Bausch and Lomb).

Matrix types include chromium plating (2.50 g/l H_2SO_4 and 250 g/l CrO_3) and polishing (750 g/l H_3PO_4 and 750 g/l H_2SO_4) solutions at typical operating levels. Chlorides of iron and chromium are used to prepare aquo-complex standards. Both standard and sample matrices require a 250 dilution to attain the 0-80 ppm detection range for these aquo-complexes.

The integrator (Dionex) attenuation is 16 and 256 for UV-visible and atomic detection, respectively.

RESULTS AND DISCUSSION

For UV-visible detection, it is necessary to choose wavelengths for Fe(III) and Cr(III) ions that have sufficient absorbance for detection but do not interfere with other matrix components present. For the chromium plating and polishing solution matrices, Table I illustrates that 300 nm and 570 nm are ideal for Fe(III) and Cr(III) ion detection, respectively. There is no interference if the ions elute simultaneously.

Separation of these aquo-complexes is accomplished by cation exchange chromatography using a nitric acid eluent which completely oxidizes sample Fe(II) ions to Fe(III) ions, while leaving Cr(III) ions intact.

TABLE I. UV-VISIBLE SPECTRA OF DETECTABLE COMPLEXES

Wavelength (nm)	Cr(IV) (1000 ppm)	Cr(III) (1000 ppm)	Fe(III) (1000 ppm)
325	1.90	0.03	1.50
330	1.93	0.04	1.57
350	2.00	0.08	1.70
375	2.10	0.21	1.60
400	2.00	0.33	0.73
425	1.90	0.27	0.18
450	1.60	0.14	0.04
475	0.49	0.07	0.01
500	0.08	0.09	0.00
525	0.01	0.16	0.00
550	0.00	0.24	0.00
570	0.00	0.28	0.00
575	0.00	0.29	0.00
600	0.00	0.26	0.00

Tables II through IV show linear calibration data in the concentration range of 0-80 ppm for aquo-complexes of iron and chromium after a 250 dilution of chromium plating and polishing solutions. The percent relative standard deviation and detection limit of the aquo-complex calibration data are two percent and 1 ppm, respectively. For retention times of the aquo-complexes of chromium and iron, UV-visible detection has 3.3 and 3.9 minutes and atomic absorption detection has 3.7 and 4.3 minutes, respectively. In these metal finishing solutions, unretained peaks did not affect determinations of the aquo-complexes using either detector.

TABLE II. DETERMINATION OF Cr(III) IONS IN CHROMIUM PLATING SOLUTIONS

Cr(III) (ppm)	UV-Visible Height (AU)	Atomic Absorption Height (AU)
20	0.053	0.051
40	0.105	0.100
60	0.157	0.150
80	0.215	0.199

TABLE III. DETERMINATION OF Fe(III) IONS IN CHROMIUM PLATING SOLUTIONS

Fe(III) (ppm)	UV-Visible Height (AU)	Atomic Absorption Height (AU)
20	0.063	0.045
40	0.122	0.085
60	0.184	0.129
80	0.243	0.171

TABLE IV. DETERMINATION OF Fe(II) IONS IN CHROMIUM POLISHING SOLUTIONS

Fe(II) (ppm)	UV-Visible Height (AU)	Atomic Absorption Height (AU)
20	0.066	0.038
40	0.126	0.077
60	0.192	0.118
80	0.266	0.147

Ion chromatography, without derivatization, is an effective method for the determination of aquo-complexes of chromium and iron in chromium plating and polishing solutions for low alloy steels. This method can be applied to other chemical matrices with proper precautions.

REFERENCES

1. J. C. VanLoon, B. Radziuk, N. Kahn, J. Lichwa, F. J. Fernandez, and J. D. Kerber, "Metal Speciation Using Atomic Absorption Spectroscopy," At. Absorpt. Newsl., Vol. 16, No. 4, 1977, p. 79-83.
2. R. Cassidy and S. Elchuk, "Dynamically Coated Columns for the Separation of Metal Ions and Anions by Ion Chromatography," Anal. Chem., Vol. 54, No. 9, 1982, pp. 1558-1563.
3. J. O'Laughlin, "Separation of Cationic Metal Chelates of 1,10-Phenanthroline by Liquid Chromatography," Anal. Chem., Vol. 54, No. 2, 1982, pp. 178-181.
4. P. Hobbs, P. Jones, and L. Ebdon, "Determination of Trace Metals by High-Performance Liquid Chromatographic Separation and Spectrophotometric Detection," Anal. Proc., Vol. 20, No. 12, 1983, pp. 613-616.
5. J. Riviello and C. Pohl, "Advances in the Separation and Detection of Transition Metals by Ion Chromatography," American Chemical Society 35th Pittsburgh Conference, 1984.
6. J. Riviello and C. Pohl, "Ion Chromatography of Transition Metals," American Chemical Society 25th Rocky Mountain Conference, 1983.
7. S. Herberling and J. Riviello, "Comparison of High Performance Chromatographic Methods for Determination of Metals," American Chemical Society 36th Pittsburgh Conference, 1985.

TECHNICAL REPORT INTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>
CHIEF, DEVELOPMENT ENGINEERING BRANCH	
ATTN: SMCAR-CCB-D	1
-DA	1
-DC	1
-DM	1
-DP	1
-DR	1
-DS (SYSTEMS)	1
CHIEF, ENGINEERING SUPPORT BRANCH	
ATTN: SMCAR-CCB-S	1
-SE	1
CHIEF, RESEARCH BRANCH	
ATTN: SMCAR-CCB-R	2
-R (ELLEN FOGARTY)	1
-RA	1
-RM	1
-RP	1
-RT	1
TECHNICAL LIBRARY	5
ATTN: SMCAR-CCB-TL	
TECHNICAL PUBLICATIONS & EDITING UNIT	2
ATTN: SMCAR-CCB-TL	
DIRECTOR, OPERATIONS DIRECTORATE	1
ATTN: SMCWV-OD	
DIRECTOR, PROCUREMENT DIRECTORATE	1
ATTN: SMCWV-PP	
DIRECTOR, PRODUCT ASSURANCE DIRECTORATE	1
ATTN: SMCWV-QA	

NOTE: PLEASE NOTIFY DIRECTOR, BENET LABORATORIES, ATTN: SMCAR-CCB-TL, OF ANY ADDRESS CHANGES.

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
ASST SEC OF THE ARMY RESEARCH AND DEVELOPMENT ATTN: DEPT FOR SCI AND TECH THE PENTAGON WASHINGTON, D.C. 20310-0103	1	COMMANDER ROCK ISLAND ARSENAL ATTN: SMCRI-ENM ROCK ISLAND, IL 61299-5000	1
ADMINISTRATOR DEFENSE TECHNICAL INFO CENTER ATTN: DTIC-FDAC CAMERON STATION ALEXANDRIA, VA 22304-6145	12	DIRECTOR US ARMY INDUSTRIAL BASE ENGR ACTV ATTN: AMXIB-P ROCK ISLAND, IL 61299-7260	1
COMMANDER US ARMY ARDEC ATTN: SMCAR-AEE	1	COMMANDER US ARMY TANK-AUTMV R&D COMMAND ATTN: AMSTA-DDL (TECH LIB) WARREN, MI 48397-5000	1
SMCAR-AES, BLDG. 321	1	COMMANDER	
SMCAR-AET-O, BLDG. 351N	1	US MILITARY ACADEMY	1
SMCAR-CC	1	ATTN: DEPARTMENT OF MECHANICS	
SMCAR-CCP-A	1	WEST POINT, NY 10996-1792	
SMCAR-FSA	1		
SMCAR-FSM-E	1	US ARMY MISSILE COMMAND	
SMCAR-FSS-D, BLDG. 94	1	REDSTONE SCIENTIFIC INFO CTR	2
SMCAR-MSI (STINFO)	2	ATTN: DOCUMENTS SECT, BLDG. 4484	
PICATINNY ARSENAL, NJ 07806-5000		REDSTONE ARSENAL, AL 35898-5241	
DIRECTOR US ARMY BALLISTIC RESEARCH LABORATORY ATTN: SLCBR-DD-T, BLDG. 305	1	COMMANDER US ARMY FGN SCIENCE AND TECH CTR ATTN: DRXST-SD	1
ABERDEEN PROVING GROUND, MD 21005-5066		220 7TH STREET, N.E. CHARLOTTESVILLE, VA 22901	
DIRECTOR US ARMY MATERIEL SYSTEMS ANALYSIS ACTV ATTN: AMXSY-MP	1	COMMANDER US ARMY LABCOM	
ABERDEEN PROVING GROUND, MD 21005-5071		MATERIALS TECHNOLOGY LAB	
COMMANDER		ATTN: SLCMT-IML (TECH LIB)	2
HQ, AMCCOM		WATERTOWN, MA 02172-0001	
ATTN: AMSMC-IMP-L	1		
ROCK ISLAND, IL 61299-6000			

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, US ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050, OF ANY ADDRESS CHANGES.

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST (CONT'D)

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
COMMANDER US ARMY LABCOM, ISA ATTN: SLCIS-IM-TL 2800 POWDER MILL ROAD ADELPHI, MD 20783-1145	1	COMMANDER AIR FORCE ARMAMENT LABORATORY ATTN: AFATL/MN EGLIN AFB, FL 32543-5434	1
COMMANDER US ARMY RESEARCH OFFICE ATTN: CHIEF, IPO P.O. BOX 12211 RESEARCH TRIANGLE PARK, NC 27709-2211	1	COMMANDER AIR FORCE ARMAMENT LABORATORY ATTN: AFATL/MNF EGLIN AFB, FL 32542-5000	1
DIRECTOR US NAVAL RESEARCH LAB ATTN: MATERIALS SCI & TECH DIVISION CODE 26-27 (DOC LIB) WASHINGTON, D.C. 20375	1 1	METALS AND CERAMICS INFO CTR BATTELLE COLUMBUS DIVISION 505 KING AVENUE COLUMBUS, OH 43201-2693	1

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, US ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050, OF ANY ADDRESS CHANGES.

END
DATE
FILMED

4-88
DTIC