

DTIC FILE COPY

AD-A221 135

2

AD

MEMORANDUM REPORT ARCCB-MR-90010

**DETERMINATION OF SODIUM CYANIDE IN COPPER
AND CADMIUM CYANIDE PLATING SOLUTIONS
BY PRECIPITATION-FORMATION TITRATION**

SAMUEL SOPOK

DTIC
ELECTE
MAY 02 1990
S B D

MARCH 1990



**US ARMY ARMAMENT RESEARCH,
DEVELOPMENT AND ENGINEERING CENTER
CLOSE COMBAT ARMAMENTS CENTER
BENÉT LABORATORIES
WATERVLIET, N.Y. 12189-4050**



APPROVED FOR PUBLIC RELEASE; DISTRIBUTION UNLIMITED

90 05 00 045

DISCLAIMER

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorized documents.

The use of trade name(s) and/or manufacturer(s) does not constitute an official indorsement or approval.

DESTRUCTION NOTICE

For classified documents, follow the procedures in DoD 5200.22-M, Industrial Security Manual, Section II-19 or DoD 5200.1-R, Information Security Program Regulation, Chapter IX.

For unclassified, limited documents, destroy by any method that will prevent disclosure of contents or reconstruction of the document.

For unclassified, unlimited documents, destroy when the report is no longer needed. Do not return it to the originator.

REPORT DOCUMENTATION PAGE		READ INSTRUCTIONS BEFORE COMPLETING FORM
1. REPORT NUMBER ARCCB-MR-90010	2. GOVT ACCESSION NO.	3. RECIPIENT'S CATALOG NUMBER
4. TITLE (and Subtitle) DETERMINATION OF SODIUM CYANIDE IN COPPER AND CADMIUM CYANIDE PLATING SOLUTIONS BY PRECIPITATION-FORMATION TITRATION		5. TYPE OF REPORT & PERIOD COVERED Final
		6. PERFORMING ORG. REPORT NUMBER
7. AUTHOR(s) Samuel Sopok		8. CONTRACT OR GRANT NUMBER(s)
9. PERFORMING ORGANIZATION NAME AND ADDRESS U.S. Army ARDEC Benet Laboratories, SMCAR-CCB-TL Watervliet, NY 12189-4050		10. PROGRAM ELEMENT, PROJECT, TASK AREA & WORK UNIT NUMBERS AMCMS No. 6126.23.1BLO.OAR PRON No. 1A02ZRAENMSC
11. CONTROLLING OFFICE NAME AND ADDRESS U.S. Army ARDEC Close Combat Armaments Center Picatinny Arsenal, NJ 07806-5000		12. REPORT DATE March 1990
		13. NUMBER OF PAGES 9
14. MONITORING AGENCY NAME & ADDRESS (if different from Controlling Office)		15. SECURITY CLASS. (of this report) UNCLASSIFIED
		15a. DECLASSIFICATION/DOWNGRADING SCHEDULE
16. DISTRIBUTION STATEMENT (of this Report) Approved for public release; distribution unlimited.		
17. DISTRIBUTION STATEMENT (of the abstract entered in Block 20, if different from Report)		
18. SUPPLEMENTARY NOTES		
19. KEY WORDS (Continue on reverse side if necessary and identify by block number) Chemical Analysis Plating Solutions Sodium Cyanide Cadmium Cyanide Plating Solutions Copper Cyanide Precipitation-Formation Titration		
20. ABSTRACT (Continue on reverse side if necessary and identify by block number) The chemical literature lacks an acceptable analytical method to adequately determine and monitor sodium cyanide in copper and cadmium cyanide plating solutions in the plating processes. In this report, an improved method is presented providing acceptable analysis and monitoring of this salt in these plating processes. The optimum operating ranges of sodium cyanide are 58 g/l maximum in copper cyanide plating solutions and 149 g/l maximum in cadmium cyanide plating solutions. The resulting precisions are in the range of (CONT'D ON REVERSE)		

20. ABSTRACT (CONT'D)

0.5 to 4.0 g/l, providing adequate monitoring of these metal finishing solutions supported by six years of testing.

UNCLASSIFIED

TABLE OF CONTENTS

	<u>Page</u>
ACKNOWLEDGEMENTS	ii
INTRODUCTION	1
EXPERIMENTAL PROCEDURE	1
RESULTS AND DISCUSSION	2
REFERENCES	4

TABLES

I. EXPERIMENTAL TITRATION DATA FOR THE SODIUM CYANIDE STANDARD SOLUTION	5
II. EXPERIMENTAL TITRATION DATA FOR SODIUM CYANIDE IN COPPER CYANIDE SAMPLE SOLUTIONS	5
III. EXPERIMENTAL TITRATION DATA FOR SODIUM CYANIDE IN CADMIUM CYANIDE SAMPLE SOLUTIONS	5
IV. PRECISION OF A 2-ml CLASS-A PIPET	6
V. PRECISION OF 50-g/l SODIUM CYANIDE STANDARD SOLUTION	6
VI. PRECISION OF 17-g/l SILVER NITRATE TITRANT SOLUTION	7
VII. PRECISION OF A 50-ml CLASS-A BURET	7

Accession For	
NTIS GRA&I	☒
DTIC TAB	☐
Unannounced	☐
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A-1	

ACKNOWLEDGEMENTS

Special thanks are given to Ellen Fogarty and Rose Neifeld of Benet Laboratories for their respective word processing and technical editing work on this manuscript.

INTRODUCTION

The chemical literature lacks an acceptable analytical method to determine and monitor sodium cyanide in copper and cadmium cyanide plating solutions in the plating processes. Lack of optimization of these plating solutions causes serious problems for the metal finishing industry such as poor quality products and wasted resources.

The only common chemical analysis method to determine sodium cyanide in these solutions is a color comparison test kit such as the one provided by Merck Corporation (ref 1), however, this method provides poor precisions.

The specific method described in this report provides acceptable analysis and monitoring of this salt in the metal finishing processes. The method consists of a precipitation-formation titration (refs 2,3). Relative precisions of this method are in the range of 2 to 3 percent.

EXPERIMENTAL PROCEDURE

Strict analytical chemistry methods and procedures are followed throughout this experimental procedure section. An excellent source of reference for these methods and procedures is by Fritz and Schenk (ref 3).

Two analytical reagent grade standard solutions are required. The first is a 50 ± 0.5 -g/l sodium cyanide solution (1 g/l sodium hydroxide for stabilization), and the second is a 17 ± 0.1 -g/l silver nitrate titrant solution. These two solutions are prepared and standardized as described in References 2 and 3.

Two other analytical reagents are required. The first is potassium iodide indicator powder and the second is concentrated ammonium hydroxide solution.

The procedure to prepare standard and sample solutions is the same for the copper and cadmium cyanide solutions. This titration analysis requires that 2 milliliters (ml) of sample or standard solution is pipetted into a 400-ml



beaker. Then the beaker is filled to about the 100-ml mark with deionized water and a stirring bar is added. Also added to this beaker is 0.1 gram of potassium iodide indicator and 5 ml of concentrated ammonium hydroxide solution.

The solution is titrated using the silver nitrate titrant to a faint yellow turbid permanent endpoint, recording the amount of titrant dispensed in milliliters at the endpoint.

All standard and sample solutions are analyzed in triplicate. Sodium cyanide concentrations in the samples are calculated by normal chemical stoichiometry.

RESULTS AND DISCUSSION

Experimental standard titration data are presented in Table I for these solutions. Experimental sample titration data are presented in Table II for copper cyanide sample solutions one and two and in Table III for cadmium cyanide sample solutions three and four. The titration can be summarized using the following equations. The first two involve the cyanide-silver ion relationship and the third involves the iodide-silver ion endpoint indicator relationship:

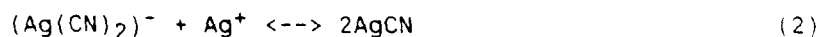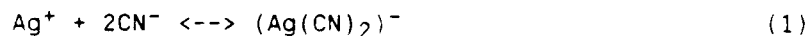

Silver cyanide will not precipitate in these solutions due to the presence of ammonium ions allowing only the silver iodide to precipitate at the endpoint.

The calculation for determining the concentration of sodium cyanide in the original sample solutions is

$$\text{g/l NaCN} = (50)(A)/(B) \quad (4)$$

where

A = ml of titrant for sample

B = ml of titrant for standard

where the constant value in Eq. (4) is the combined result of many constants (refs 2,3).

Using Eq. (4), the respective sodium cyanide values for copper cyanide sample solutions one and two in Table II are 35.1 and 49.9 g/l. Likewise, using the same equation, the respective sodium cyanide values for cadmium cyanide sample solutions three and four in Table III are 84.9 and 140.1 g/l.

It is useful to evaluate the variations in precision for the materials and methods used. Tables IV through VII present the data for the 2-ml class-A pipets, 50-g/l sodium cyanide standard solution, 17-g/l silver nitrate titrant solution, and the 50-ml class-A burets, respectively.

The data obtained by this method are sufficient to adequately monitor the sodium cyanide in the metal finishing processes, thus providing efficient use of resources. The optimum operating ranges of sodium cyanide are 58 g/l maximum in the copper cyanide plating solutions and 149 g/l maximum in the cadmium cyanide plating solutions. The resulting precisions are in the range of 0.5 to 4.0 g/l, providing adequate monitoring of these metal finishing solutions supported by six years of testing.

REFERENCES

1. "EM Quant Chlorine and Cyanide Tests," EM Science Division of Merck Corporation, Cherry Hill, NJ, 1980.
2. D. Peters, J. Hayes, and G. Hieftje, Chemical Separations and Measurements: Theory and Practice of Analytical Chemistry, W. B. Saunders Company, Philadelphia, PA, 1974.
3. J. Fritz and G. Schenk, Quantitative Analytical Chemistry, Fifth Edition, Allyn and Bacon, Inc., Boston, MA, 1987.

TABLE I. EXPERIMENTAL TITRATION DATA FOR THE
SODIUM CYANIDE STANDARD SOLUTION

Replicates	Titrant Used (ml)
1	10.25
2	10.20
3	10.20
$\bar{x}$ (avg)	10.22

TABLE II. EXPERIMENTAL TITRATION DATA FOR SODIUM
CYANIDE IN COPPER CYANIDE SAMPLE SOLUTIONS

Replicates	Sample One Titrant Used (ml)	Sample Two Titrant Used (ml)
1	7.15	10.20
2	7.15	10.20
3	7.20	10.20
$\bar{x}$ (avg)	7.17	10.20

TABLE III. EXPERIMENTAL TITRATION DATA FOR SODIUM
CYANIDE IN CADMIUM CYANIDE SAMPLE SOLUTIONS

Replicates	Sample Three Titrant Used (ml)	Sample Four Titrant Used (ml)
1	17.35	28.60
2	17.35	28.65
3	17.35	28.65
$\bar{x}$ (avg)	17.35	28.63

TABLE IV. PRECISION OF A 2-ml CLASS-A PIPET

Replicate	2-ml Pipet volume (ml)*
1	2.00
2	2.01
3	1.99
4	1.99
5	1.99
6	2.01
$\bar{x}$ (avg)	2.00
S_n	0.01

*Volumes are calculated from the weight-volume relationship of a pipetted deionized water solution corrected for temperature.

TABLE V. PRECISION OF A 50-g/l SODIUM CYANIDE STANDARD SOLUTION

Replicate	Sodium Cyanide (g/l)*
1	49.97
2	49.99
3	49.99
4	50.01
5	50.00
6	50.02
$\bar{x}$ (avg)	50.00
S_n	0.02

*Sodium cyanide concentrations are calculated by methods given in References 2 and 3.

TABLE VI. PRECISION OF A 17-g/l SILVER NITRATE TITRANT SOLUTION

Replicate	Silver Nitrate (g/l)*
1	17.1
2	17.0
3	17.2
4	16.8
5	16.9
6	16.8
X(avg)	17.0
Sn	0.2

*Silver nitrate concentrations are calculated by methods given in References 2 and 3.

TABLE VII. PRECISION OF A 50-ml CLASS-A BURET

Replicate	Volume (ml)*
1	24.94
2	24.98
3	25.02
4	25.05
5	24.98
6	25.05
X(avg)	25.00
Sn	0.04

*Volumes are calculated from the weight-volume relationship of a contained deionized water solution corrected for temperature at a point one-half full--25 ml.

TECHNICAL REPORT INTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>
CHIEF, DEVELOPMENT ENGINEERING DIVISION	
ATTN: SMCAR-CCB-D	1
-DA	1
-DC	1
-DM	1
-DP	1
-DR	1
-DS (SYSTEMS)	1
CHIEF, ENGINEERING SUPPORT DIVISION	
ATTN: SMCAR-CCB-S	1
-SE	1
CHIEF, RESEARCH DIVISION	
ATTN: SMCAR-CCB-R	2
-RA	1
-RM	1
-RP	1
-RT	1
TECHNICAL LIBRARY	5
ATTN: SMCAR-CCB-TL	
TECHNICAL PUBLICATIONS & EDITING SECTION	3
ATTN: SMCAR-CCB-TL	
DIRECTOR, OPERATIONS DIRECTORATE	1
ATTN: SMCWV-OD	
DIRECTOR, PROCUREMENT DIRECTORATE	1
ATTN: SMCWV-PP	
DIRECTOR, PRODUCT ASSURANCE DIRECTORATE	1
ATTN: SMCWV-QA	

NOTE: PLEASE NOTIFY DIRECTOR, BENET LABORATORIES, ATTN: SMCAR-CCB-TL, OF ANY ADDRESS CHANGES.

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
ASST SEC OF THE ARMY RESEARCH AND DEVELOPMENT ATTN: DEPT FOR SCI AND TECH THE PENTAGON WASHINGTON, D.C. 20310-0103	1	COMMANDER ROCK ISLAND ARSENAL ATTN: SMCRI-ENM ROCK ISLAND, IL 61299-5000	1
ADMINISTRATOR DEFENSE TECHNICAL INFO CENTER ATTN: DTIC-FDAC CAMERON STATION ALEXANDRIA, VA 22304-6145	12	DIRECTOR US ARMY INDUSTRIAL BASE ENGR ACTV ATTN: AMXIB-P ROCK ISLAND, IL 61299-7260	1
COMMANDER US ARMY ARDEC ATTN: SMCAR-AEE	1	COMMANDER US ARMY TANK-AUTMV R&D COMMAND ATTN: AMSTA-ODL (TECH LIB) WARREN, MI 48397-5000	1
SMCAR-AES, BLDG. 321	1	COMMANDER US MILITARY ACADEMY ATTN: DEPARTMENT OF MECHANICS WEST POINT, NY 10996-1792	1
SMCAR-AET-O, BLDG. 351N	1		
SMCAR-CC	1		
SMCAR-CCP-A	1		
SMCAR-FSA	1		
SMCAR-FSM-E	1	US ARMY MISSILE COMMAND REDSTONE SCIENTIFIC INFO CTR ATTN: DOCUMENTS SECT, BLDG. 4484 REDSTONE ARSENAL, AL 35898-5241	2
SMCAR-FSS-D, BLDG. 94	1		
SMCAR-IMI-I (STINFO) BLDG. 59	2		
PICATINNY ARSENAL, NJ 07806-5000			
DIRECTOR US ARMY BALLISTIC RESEARCH LABORATORY ATTN: SLCBR-DD-T, BLDG. 305 ABERDEEN PROVING GROUND, MD 21005-5066	1	COMMANDER US ARMY FGN SCIENCE AND TECH CTR ATTN: DRXST-SD 220 7TH STREET, N.E. CHARLOTTESVILLE, VA 22901	1
DIRECTOR US ARMY MATERIEL SYSTEMS ANALYSIS ACTV ATTN: AMXSY-MP ABERDEEN PROVING GROUND, MD 21005-5071	1	COMMANDER US ARMY LABCOM MATERIALS TECHNOLOGY LAB ATTN: SLCMT-IML (TECH LIB) WATERTOWN, MA 02172-0001	2
COMMANDER HQ, AMCCOM ATTN: AMSMC-IMP-L ROCK ISLAND, IL 61299-6000	1		

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, US ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050, OF ANY ADDRESS CHANGES.

TECHNICAL REPORT EXTERNAL DISTRIBUTION LIST (CONT'D)

	<u>NO. OF COPIES</u>		<u>NO. OF COPIES</u>
COMMANDER US ARMY LABCOM, ISA ATTN: SLCIS-IM-TL 2800 POWDER MILL ROAD ADELPHI, MD 20783-1145	1	COMMANDER AIR FORCE ARMAMENT LABORATORY ATTN: AFATL/MN EGLIN AFB, FL 32542-5434	1
COMMANDER US ARMY RESEARCH OFFICE ATTN: CHIEF, IPO P.O. BOX 12211 RESEARCH TRIANGLE PARK, NC 27709-2211	1	COMMANDER AIR FORCE ARMAMENT LABORATORY ATTN: AFATL/MNF EGLIN AFB, FL 32542-5434	1
DIRECTOR US NAVAL RESEARCH LAB ATTN: MATERIALS SCI & TECH DIVISION CODE 26-27 (DOC LIB) WASHINGTON, D.C. 20375	1 1	METALS AND CERAMICS INFO CTR BATTELLE COLUMBUS DIVISION 505 KING AVENUE COLUMBUS, OH 43201-2693	1

NOTE: PLEASE NOTIFY COMMANDER, ARMAMENT RESEARCH, DEVELOPMENT, AND ENGINEERING CENTER, US ARMY AMCCOM, ATTN: BENET LABORATORIES, SMCAR-CCB-TL, WATERVLIET, NY 12189-4050, OF ANY ADDRESS CHANGES.