

USAARL Report No. 92-11

DTIC

ELECTRONIC

APR 28 1992

S

C

AD-A249 633



**Test and Evaluation Report
of the SSCOR PacVac
Portable Suction Pump Model 10002**

By

**James E. Bruckart (Project Officer)
and
Joseph R. Licina (Project Officer)**

**Bill Olding (UES, Inc.)
Martin Quattlebaum (UES, Inc.)**

Biodynamics Research Division

February 1992

92-10637



Approved for public release; distribution unlimited.

92 4 24 110

**United States Army Aeromedical Research Laboratory
Fort Rucker, Alabama 36362-5292**

Notice

Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.


Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

Disclaimer

The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:

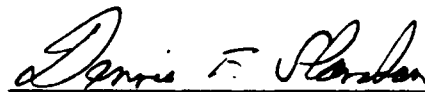

DENNIS F. SHANAHAN
LTC, MC, MFS
Director, Biodynamics
Research Division

Released for publication:



ROGER W. WILEY, O.D., Ph.D.
Chairman, Scientific
Review Committee



DAVID H. KARNEY
Colonel, MC, SFS
Commanding

REPORT DOCUMENTATION PAGE

Form Approved
OMB No. 0704-0188

1a. REPORT SECURITY CLASSIFICATION UNCLASSIFIED			1b. RESTRICTIVE MARKINGS		
2a. SECURITY CLASSIFICATION AUTHORITY			3. DISTRIBUTION/AVAILABILITY OF REPORT Approved for public release, distribution unlimited		
2b. DECLASSIFICATION/DOWNGRADING SCHEDULE					
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 92-11			5. MONITORING ORGANIZATION REPORT NUMBER(S)		
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		6b. OFFICE SYMBOL (If applicable) SGRD-UAD-IE	7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command		
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 577 Fort Rucker, AL 36362-5292			7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012		
8a. NAME OF FUNDING/SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (If applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER Partial effort, under contract No. DAMD 17-86-C-6215		
8c. ADDRESS (City, State, and ZIP Code)			10. SOURCE OF FUNDING NUMBERS		
PROGRAM ELEMENT NO. 0603807A		PROJECT NO. 3M463807D836	TASK NO. 6 LC	WORK UNIT ACCESSION NO. 201	
11. TITLE (Include Security Classification) Test and Evaluation Report of the SSCOR PacVac Portable Suction Pump Model 10002					
12. PERSONAL AUTHOR(S) James E. Bruckart, Joseph R. Licina, Bill Olding, and Martin Quattlebaum					
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____		14. DATE OF REPORT (Year, Month, Day) 1992 February	
15. PAGE COUNT 45					
16. SUPPLEMENTARY NOTATION					
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)		
FIELD	GROUP	SUB-GROUP	Electromagnetic compatibility, test and evaluation, aeromedical equipment		
19. ABSTRACT (Continue on reverse if necessary and identify by block number) The SSCOR PacVac Portable Suction Pump Model 10002 was tested for electromagnetic interference/compatibility in the UH-60A helicopter under the U.S. Army Program for Testing and Evaluation of Equipment for Aeromedical Operations. The tests were conducted using current military and industrial standards and procedures for electromagnetic interference/compatibility and human factors. Since in-flight testing of the PacVac portable suction unit could not be completed, its compatibility with the U.S. Army medical evacuation UH-60A Blackhawk could not be determined.					
20. DISTRIBUTION/AVAILABILITY OF ABSTRACT ☒ UNCLASSIFIED/UNLIMITED ☐ SAME AS RPT. ☐ DTIC USERS			21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED		
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center			22b. TELEPHONE (Include Area Code) (205) 255-6907		22c. OFFICE SYMBOL SGRD-UAX-SI

Acknowledgment

Laboratory testing was accomplished at the U.S. Army Aero-medical Research Laboratory using government furnished equipment by Universal Energy Systems, Inc. (UES), under contract No. DAMD 17-86-C-6215.

Accession For	
DD Form 891	☒
DD Form 748	☐
Form needed	☐
Classification	
Description/	
Availability Codes	
Special and/or	
Dist	Special
A-1	



Table of contents

Section	Page
1. <u>EXECUTIVE DIGEST</u>	
1.1 Test objectives	1-1
1.2 Testing authority	1-2
1.3 Scope	1-2
1.4 Material description	1-2
1.5 Summary	1-3
1.6 Conclusion	1-3
2. <u>SUBTESTS</u>	
2.1 Initial inspection	2-1
2.2 Battery life evaluation	2-1
2.3 Human factors evaluation (laboratory)	2-2
2.4 Altitude (low pressure) test	2-3
2.5 Vibration test	2-4
2.6 High temperature test	2-5
2.7 Low temperature test	2-6
2.8 Humidity test	2-8
2.9 Electromagnetic characteristics test	2-9
2.10 In-flight human factors evaluation	2-10
2.11 In-flight EMI/EMC characteristics test	2-10

Table of contents (Continued)

3. SUPPORTING DOCUMENTATION

3.1	Detailed test information	3-1
3.2	Test data	3-2
3.3	Criteria, significant problems, and suggested improvements	3-25
3.4	References	3-27
3.5	Abbreviations	3-28
3.6	List of manufacturers.	3-30
3.7	Distribution list	3-31

Section 1. Executive digest

The Army program for Test and Evaluation of Aeromedical Equipment uses existing military standards (MIL-STD) and collective professional expertise to test and evaluate selected medical equipment proposed for use aboard Army aircraft. Equipment meeting these standards ensures the safety of the crew, patients, and aircraft by eliminating risks due to: (1) Interference by the medical equipment with aircraft systems/subsystems operation, (2) by the aircraft system's interference with the operation of the medical equipment, (3) the medical equipment's susceptibility to environmental exposure, or (4) physical and/or functional incompatibility while in use on board selected rotary-wing aircraft. This program tests both developmental and nondevelopmental (off the shelf) medical equipment destined for use aboard Army medical evacuation aircraft.

1.1 TEST OBJECTIVES

1.1.1 To determine if the medical equipment is complete and operational per the manufacturer's operating instructions.

1.1.2 To ensure the electrical safety of the medical equipment.

1.1.3 To ensure the equipment will function as designed throughout the rated battery operation time.

1.1.4 To ensure the safety of the operator, the patient, and the aircrew.

1.1.5 To assess design considerations which, potentially, could contribute to an operator error.

1.1.6 To determine if the medical equipment can function as designed in a low-pressure environment.

1.1.7 To determine the ability of the medical equipment to withstand the vibrational stresses expected in a rotary-wing flight environment without degradation or malfunction.

1.1.8 To determine the ability of the medical equipment to be stored and operated in a high-temperature environment.

1.1.9 To determine the ability of the medical equipment to be stored and operated in a low-temperature environment.

1.1.10 To determine the ability of the medical equipment to operate satisfactorily for short periods of time during exposure to highly humid conditions.

1.1.11 To assess the levels of electromagnetic emissions produced by the medical equipment within selected frequency ranges.

1.1.12 To assess the minimum electromagnetic susceptibility levels of the medical equipment within selected frequency ranges.

1.1.13 To assess the physical and/or functional compatibility of the medical equipment while in use on board the aircraft.

1.1.14 To assess the electromagnetic interference (EMI) and electromagnetic compatibility (EMC) characteristics of the medical equipment with the host aircraft and its installed systems.

1.2 TESTING AUTHORITY

Research and Technology Work Unit Summary, dated 5 October 1989. Project number 3M463807D836, titled, Army Program for Testing and Evaluation of Equipment for Aeromedical Operations.

1.3 SCOPE

1.3.1 This test was conducted at the United States Army Aeromedical Research Laboratory (USAARL), Cairns Army Airfield (CAAF), and designated test flight areas in and around Fort Rucker, Alabama.

1.3.2 The USAARL UH-60A aircraft, serial number 88-26069, with subsystems delineated in paragraph 3.2.2, was available for flight testing of the SSCOR PacVac portable suction pump*. The PacVac suction pump was not available for flight testing.

1.3.3 Laboratory testing was accomplished at USAARL using government furnished equipment (GFE) by Universal Energy Systems, Inc. (UES), under contract No. DAMD 17-86-C-6215.

1.3.4 Prior to flight testing the following tests were accomplished: Acceptance inspection, equipment training, electromagnetic compatibility, human factors and safety, environmental compatibility, and in-flight compatibility.

1.3.5 An airworthiness release (AWR) dated 16 August 1990 was received from the U.S. Army Aviation Systems Command (AVSCOM).

1.4 MATERIAL DESCRIPTION

The SSCOR PacVac portable suction pump is a battery-operated suction pump for pharyngeal suctioning during resuscitation emergencies. Batteries are recharged by an ac charger included with the unit. It is designed to be ready for immediate suctioning and to provide 30 minutes of suction in the event of a

* See list of manufacturers

resuscitation emergency. Suction power is switchable between two power settings. The pump also is equipped with a disposable collection canister containing a hydrophobic filter element to prevent aerosolization of particles down to 0.45 micron diameter.

1.5 SUMMARY

1.5.1 Laboratory testing

1.5.1.1 Battery Life Evaluation: The battery in the PacVac portable suction unit provided power for an average of 45 minutes in three cycles of charging and operation. This is the maximum time specified in the operation manual for the suction unit.

1.5.1.2 Human Factors Evaluation: The PacVac suction unit was found satisfactory in all major categories of the evaluation. Standards referenced include MIL-STD-1472D, AAMI Human Factors Engineering Guidelines, and UL-544.

1.5.1.3 Environmental Tests: The PacVac portable suction can be expected to perform in a variety of environmental conditions. Its performance was found to be satisfactory in all stages of the environmental testing.

1.5.1.4 Radiated Emissions Tests (RE02): The PacVac portable suction unit may be unsatisfactory for use in certain EMI sensitive environments. Broadband emissions were detected in the test frequency ranges that exceeded the test limits. Emission limits are set forth in MIL-STD-461A, Notice 4.

1.5.1.5 Radiated Susceptibility Test (RS03): The PacVac portable suction unit was not found to be susceptible to radiated emissions.

1.5.2 In-flight testing

In-flight testing of the PacVac portable suction unit could not be accomplished. The manufacturer of the PacVac suction unit, Model 10002, discontinued production of this unit and none was available from the manufacturer or local military medical units for in-flight testing.

1.6 CONCLUSIONS

Since in-flight testing of the PacVac portable suction unit could not be completed, its compatibility with the U.S. Army medical evacuation UH-60A Blackhawk with the subsystems listed in paragraph 3.2.2 could not be determined.

Section 2. Subtests

2.1 INITIAL INSPECTION

2.1.1 Objective

To determine if the PacVac portable suction unit is complete and operational for testing per the manufacturer's operating instructions.

2.1.2 Criteria

2.1.2.1 The physical inventory was conducted solely for investigation and documentation.

2.1.2.2 The PacVac suction will suction 550 mL of water through a vacuum line in 7.5 seconds or less. This is equivalent to 550 mmHg suction or 25 liters per minute airflow.

2.1.3 Test procedure

2.1.3.1 A complete physical inventory of the PacVac suction unit was completed per the manufacturer's equipment list.

2.1.3.2 An operational validation test of the PacVac suction unit was conducted per the manufacturer's operating instructions by USAARL's medical maintenance personnel.

2.1.4 Test findings

2.1.4.1 The PacVac suction unit was inventoried and found to be complete.

2.1.4.2 The PacVac suction unit operated as prescribed in the manufacturer's operating manual, P/N not available. Criterion met.

2.2 BATTERY LIFE EVALUATION (Laboratory)

2.2.1 Objective

To ensure the equipment will function as designed throughout the rated battery operation time.

2.2.2 Criterion

Verify manufacturer's specified full power battery life expectancy of 30 to 45 minutes continuous operation.

2.2.3 Test procedure

2.2.3.1 Charging and operation cycles were conducted in ambient room conditions of 23°C, 40-60 percent relative humidity (RH).

2.2.3.2 The PacVac suction unit was operated continuously at full speed with no load until no suction was detected.

2.2.4 Test findings

The suction unit operated continuously for an average of 45 minutes before a low battery condition was noted. This performance is consistent with the manufacturer's specification of 30 to 45 minutes operation on a fully charged battery. Criterion met.

2.3 HUMAN FACTORS EVALUATION (Laboratory)

2.3.1 Objectives

2.3.1.1 To assure the safety of the operator, the potential patient, and the aircrew.

2.3.1.2 To assess the design considerations which could potentially contribute to an operator error.

2.3.2 Criterion

The PacVac portable suction unit must be rated satisfactory in all major categories of the evaluation. These include visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.3.3 Test procedure

2.3.3.1 The evaluation was conducted in a laboratory under fluorescent lighting and ambient room conditions.

2.3.3.2 The PacVac suction unit was operated according to prescribed instructions through its full range of functions.

2.3.4 Test finding

The PacVac suction unit was found to be satisfactory in all major areas of the evaluation. Criterion met.

2.4 ALTITUDE (LOW PRESSURE) TEST [IAW METHOD 500.2, MIL-STD-810D]

2.4.1 Objective

To determine if the PacVac suction unit can function as designed in a low-pressure environment.

2.4.2 Criterion

The PacVac suction unit will function as designed while exposed to an altitude equivalency of 15,000 feet above sea level.

2.4.3 Test procedure

2.4.3.1 A pretest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.4.3.2 The Altitude Test was performed in a Tenney Engineering model 64S altitude chamber*. This test is based on MIL-STD-810D, Method 500.2. The PacVac suction unit was placed in operation near the center of the floor of the chamber. The PacVac suction unit was turned on by means of a remote arm through the chamber wall and monitored for ability to suction 550 mL water in less than 7.5 seconds. Chamber pressure was decreased to 420 mmHg (15,000 ft equivalent altitude) over a 15-minute period, held constant for 60 minutes, then raised, at 1500 fpm, to ambient conditions (760 mmHg) over a 10-minute period. There are no provisions for the control of temperature or humidity inside this chamber.

2.4.3.3 A posttest performance check was conducted to ensure proper operation of the PacVac suction unit after the exposure to low pressure.

2.4.4 Test findings

2.4.4.1 The pretest performance check met criterion 2.1.2.2.

2.4.4.2 No failures in the PacVac suction unit's performance were noted before, during, or after the altitude test. Criterion met.

2.4.4.3 The posttest performance check met criterion 2.1.2.2.

2.5 VIBRATION TEST [IAW METHOD 514.3, MIL-STD-810D]

2.5.1 Objective

To determine the ability of the PacVac suction unit to withstand the vibrational stresses expected in a rotary-wing environment without degradation or malfunction.

2.5.2 Criterion

While exposed to vibrational stresses, the PacVac suction unit will remain operational.

2.5.3 Test procedure

2.5.3.1 A pretest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.5.3.2 The vibration test was performed using an Unholtz-Dickey model TA115-40/CSTA vibration test system*. It is a single-axis system with an electromagnetic driver unit. The test consisted of sinusoidal vibrations, superimposed on random vibrations over a frequency range of 500 Hz, as shown below. These vibrations are derived from measurements taken on the floor under the co-pilot's seat in a UH-1 helicopter traveling at 120 knots. The reference spectrum breakpoints are from MIL-STD-810D, Method 514.3; reference spectrum levels are based on field measurements with a conservatism factor of 1.5. Independent tests were conducted in the X, Y, and Z axes.

Z-axis

duration: 60 minutes

broadband intensity: 0.4506 G_{rms}

random vibration: initial slope : 99.00 dB/Hz

5 Hz level: 0.00006210 $G_{sqr/Hz}$

100 Hz level: 0.0006210 $G_{sqr/Hz}$

300 Hz level: 0.0006210 $G_{sqr/Hz}$

500 Hz level: 0.00006210 $G_{sqr/Hz}$

final slope: -99.00 dB/oct

sinusoidal vibration: .5450 G_{pk} at 11.25 Hz

.1690 G_{pk} at 22.50 Hz

.1200 G_{pk} at 33.75 Hz

.0310 G_{pk} at 45.00 Hz

.0530 G_{pk} at 56.25 Hz

X and Y axes

duration: 60 minutes each
broadband intensity: 0.3099 G_{rms}
random vibration: initial slope: 99.00 dB/oct
5 Hz level: 0.00002920 $G_{sqr/Hz}$
100 Hz level: 0.0002920 $G_{sqr/Hz}$
300 Hz level: 0.0002920 $G_{sqr/Hz}$
500 Hz level: 0.00002920 $G_{sqr/Hz}$
final slope: -99.00 dB/oct
sinusoidal vibration: .3200 G_{pk} at 11.25 Hz
.0670 G_{pk} at 22.50 Hz
.0950 G_{pk} at 33.75 Hz
.0350 G_{pk} at 45.00 Hz
.0770 G_{pk} at 56.25 Hz

The PacVac suction unit was strapped to the vibration table fixture, and its performance was evaluated before, during, and after exposure to vibration.

2.5.3.3 A posttest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.5.4 Test findings

2.5.4.1 The pretest performance check met criterion 2.1.2.2.

2.5.4.2 No failures in the PacVac suction unit's performance occurred before, during, or after exposure to vibration. Criterion met.

2.5.4.3 The posttest performance check met criterion 2.1.2.2.

2.6 HIGH TEMPERATURE TEST [IAW METHOD 501.2, MIL-STD-810D]

2.6.1 Objective

To determine the ability of the PacVac suction unit to be stored and operated in a high-temperature environment.

2.6.2 Criteria

2.6.2.1 During the high-temperature operation check, the PacVac suction unit must display consistent performance.

2.6.2.2 After the high-temperature storage cycle, the PacVac suction unit must be able to display consistent performance.

2.6.3 Test procedure

2.6.3.1 A pretest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.6.3.2 The high-temperature test was conducted in a Tenney Engineering model ZWUL-10107D Walk-in Controlled Environment Chamber*. This test is based on MIL-STD-810D, Method 501.2. For the high-temperature operation test, the PacVac suction unit was placed in operation on a wire test stand near the center of the environmental chamber. The chamber temperature was raised to 49°C and the humidity was stabilized at a maximum of 20 percent RH within 15 minutes. The environmental control system is capable of regulating temperature within $\pm 2^\circ\text{C}$ and humidity within ± 5 percent RH. Temperature and humidity were held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to minimize the change in chamber conditions during performance checks. After the operational test, the PacVac suction unit was allowed to return to ambient conditions over a 30-minute period.

2.6.3.3 A posttest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.6.3.4 The PacVac suction unit was stored (not operated) at temperatures of 63°C for 1 hour, 71°C for 4 hours, then again at 63°C for 1 hour. The chamber and PacVac suction unit then were returned to ambient conditions over a 30-minute period.

2.6.3.5 A poststorage performance check was conducted to ensure proper performance of the PacVac suction unit.

2.6.4 Test findings

2.6.4.1 The pretest performance check met criterion 2.1.2.2.

2.6.4.2 No operational failures occurred during the high-temperature test. Criterion met.

2.6.4.3 The posttest performance check met criterion 2.1.2.2.

2.6.4.4 The PacVac suction unit functioned properly after the high-temperature storage test. Criterion met.

2.7 LOW TEMPERATURE TEST [IAW METHOD 502.2, MIL-STD-810D]

2.7.1 Objective

To determine the ability of the PacVac suction unit to be stored and operated in a low-temperature environment.

2.7.2 Criteria

2.7.2.1 During the low-temperature operation check, the PacVac suction unit must display consistent operation.

2.7.2.2 After the low-temperature storage cycle, the PacVac suction unit must be able to display consistent operation.

2.7.3 Test procedure

2.7.3.1 A pretest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.7.3.2 The PacVac suction unit was placed on the floor of the environmental chamber and the temperature was lowered to 0°C within 25 minutes. The environmental control system is capable of regulating temperature within 2°C. Humidity cannot be controlled in the chamber at freezing temperatures. The temperature was held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to minimize the change in chamber conditions, and a performance check was conducted. The chamber temperature then was raised to ambient temperature within a 30-minute period.

2.7.3.3 A posttest performance check was conducted to ensure proper operation of the PacVac suction unit.

2.7.3.4 The PacVac suction unit was "stored" in a nonoperational mode with the power cord coiled and placed on top of the PacVac suction unit. The PacVac suction unit was placed on the floor of the environmental test chamber and the temperature was lowered to -46°C for 6 hours. The chamber then was raised to ambient temperature over a 30-minute period.

2.7.3.5 A poststorage performance check was conducted to ensure proper operation of the PacVac suction unit.

2.7.4 Test findings

2.7.4.1 The pretest performance check met criterion 2.1.2.2.

2.7.4.2 No operational failures occurred during the low-temperature test. Criterion met.

2.7.4.3 The posttest performance check met criterion 2.1.2.2.

2.7.4.4 The PacVac suction unit functioned properly after the low-temperature storage test. Criterion met.

2.8 HUMIDITY TEST [IAW METHOD 507.2, MIL-STD-810D]

2.8.1 Objective

To determine the ability of the PacVac suction unit to operate satisfactorily for short periods of time during exposure to highly humid conditions.

2.8.2 Criterion

While exposed to a high humidity environment, the PacVac suction unit must display consistent operation.

2.8.3 Test procedure

2.8.3.1 A pretest performance check was conducted to ensure the proper operation of the PacVac suction unit.

2.8.3.2 The humidity test was conducted in a Tenney Engineering model ZWUL-10107D Walk-in Controlled Environment Chamber*. This test is based on MIL-STD-810D, Method 507.2. For the humidity test, the PacVac suction unit was placed in operation on a wire test stand near the center of the environmental chamber. The chamber temperature was raised to a temperature of 29.5°C and a relative humidity of 95 percent within 25 minutes. Temperature and relative humidity were maintained for 4 hours. The environmental control system is capable of regulating temperature within $\pm 2^{\circ}\text{C}$ and humidity within ± 5 percent RH. At 45-minute intervals, the suction performance was checked. The chamber door was opened briefly to minimize the change in chamber conditions. The chamber and the PacVac suction unit were returned to ambient conditions before the posttest performance validation check was conducted.

2.8.3.3 A posttest performance check was conducted to ensure the proper operation of the PacVac suction unit.

2.8.4 Test findings

2.8.4.1 The pretest performance check met criterion 2.1.2.2.

2.8.4.2 No operational failures occurred during the humidity test. Criterion met.

2.8.4.3 The posttest performance check met criterion 2.1.2.2.

2.9 ELECTROMAGNETIC CHARACTERISTICS TEST [IAW MIL-STD-461A, Notice 4, and MIL-STD-462, Notice 3]

2.9.1 Objectives

2.9.1.1 To assess the maximum levels of radiated electromagnetic emissions produced by the PacVac suction unit in the 14 kHz to 1.0 GHz frequency range.

2.9.1.2 To assess the tolerances of radiated electromagnetic susceptibility of the PacVac suction unit within the 10 kHz to 10 GHz broadband electric field and the 14 kHz to 12.4 GHz narrowband.

2.9.2 Criteria

2.9.2.1 The PacVac suction unit shall not produce emissions in excess of the limits set forth in paragraph 6.13, MIL-STD-461A, Notice 4.

2.9.2.2 The PacVac suction unit shall not malfunction when it is subjected to radiated emissions as specified in paragraph 6.20, MIL-STD-461A, Notice 4.

2.9.3 Test procedure

2.9.3.1 The radiated emissions test was performed according to MIL-STD-462, Notice 3, Method RE02. The PacVac suction unit was positioned on a wooden test stand 1 meter tall, 0.18 meters wide, and 0.21 meters long, inside the EMI chamber. The unit was directly in line with, and at a horizontal distance of 1 meter from the receiving antennas. The antennas were mounted for both vertical and horizontal polarities and connected to the appropriate EMI receivers. Electrometrics EMC-25 and EMC-50 receivers were used for this test. Their frequency ranges in testing are 14 kHz to 1 GHz and 1 to 12.4 GHz. Broadband and narrowband detection methods were used from 14 kHz to 1 GHz. Narrowband detection methods were used from 1 to 12.4 GHz. The suction unit operated continuously when tested.

2.9.3.2 The radiated susceptibility test was performed according to MIL-STD-462, Notice 3, Method RS03. The PacVac suction unit was positioned on a wooden test stand 1 meter tall, 0.18 meters wide, and 0.21 meters long, inside the EMI chamber. The unit was directly in line with, and at a horizontal distance of 1 meter from the transmitting antennas. The antennas were mounted for both vertical and horizontal polarities and connected to radio frequency (RF) transmitters. The suction unit was exposed to fields of 10 V/m from 200 MHz to 2 GHz, and 5 V/m from 2 to 10 GHz. All RF carrier waves were 50 percent amplitude modulated with a 1000 Hz tone.

2.9.4 Test findings

2.9.4.1 During the radiated emissions test, broadband emissions from 14 kHz to 718 MHz were detected which exceeded specification limits of MIL-STD-461A, Notice 4, by 0.2 to 24.8 dB. Criterion partially met.

2.9.4.2 The PacVac suction unit was not susceptible to radiated emissions. Criterion met.

2.10 IN-FLIGHT HUMAN FACTORS EVALUATION

2.10.1 Objective

To assess the physical and/or functional compatibility of the PacVac suction unit while in use on board the aircraft.

2.10.2 Criterion

The flight surgeon shall be able to operate the PacVac suction unit without physical or functional restrictions aboard the aircraft. Major areas of concern include: Proper operation, visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

2.10.3 Test procedure

Not applicable since unit not available for in-flight testing.

2.10.4 Test findings

Not available since suction unit not available for in-flight testing. Criterion not evaluated.

2.11 IN-FLIGHT EMI/EMC CHARACTERISTICS TEST

2.11.1 Objective

To assess the EMI/EMC characteristics of the PacVac suction unit with the host aircraft and its installed systems.

2.11.2 Criteria

2.11.2.1 The PacVac suction unit shall not radiate EMI to disrupt or interfere with other equipment or systems aboard the aircraft.

2.11.2.2 The aircraft shall not radiate EMI to disrupt or interfere with the PacVac suction unit's operation.

2.11.3 Test procedure

A qualitative EMI/EMC assessment will be performed with both the PacVac suction unit and the aircraft operating as source and victim. The PacVac suction unit and applicable aircraft instruments and systems will be monitored for unusual operation, readings, surges, or power anomalies.

2.11.4 Test findings

Not available since suction unit not available for in-flight testing. Criterion not evaluated.

Section 3. Supporting documentation

3.1 DETAILED TEST INFORMATION

3.1.1 General information

3.1.1.1 PacVac suction unit testing is not considered a major action significantly affecting the quality of the human environment and, therefore, qualifies for categorical exclusion A-28, AR 200-1, Appendix A.

3.1.1.2 A safety pilot will be designated for each flight. Flight operations will be conducted IAW the aircraft operator's manual, appropriate aircrew training manuals, and test item technical data.

3.1.2 Material description

3.1.2.1 The PacVac portable suction unit is a rechargeable, battery-operated pump for pharyngeal suctioning during resuscitation emergencies. The pump operates on 12 volts dc, supplied by an rechargeable internal sealed lead acid battery. The battery can be recharged in 4 to 6 hours using the 120 Vac recharger supplied with the unit. It is designed to be available for immediate operation and to operate on battery power for 30 to 45 minutes. Suction power is switchable between two power settings. This provides 550 mmHg suction with full power or 80 mmHg with the regulator valve open. Secretions are collected in a collection canister. Standard accessories include an 800 mL collection canister, vacuum connecting line, battery charger, and cordura carrying case.

3.1.2.2 Dimensions: 10 x 8 x 15 inches (25.4 x 20.3 x 38.1 cm).

3.1.2.3 Weight: 11 lbs (5 kg).

3.1.2.4 Power requirements: 12 volts dc supplied by internal sealed lead-acid battery. Battery recharged by 120 Vac recharger. The unit will operate the pump for 30 to 45 minutes on a fully charged battery and requires 4 to 6 hours to recharge the battery.

3.2 TEST DATA

3.2.1 Photographic description



3.2.2 Aircraft equipment list

Item No.	Nomenclature
1	Receiver radio -- R-1496A/ARN-89 (automatic direction finder)
2	Displacement gyro -- CN-1314/A
3	Gyro directional -- CN-998/ASN-43
4	Signal data converter -- CV-3338/ASN-128
5	Receiver -- R-2139/ARN-123 (VOR/LOC/MB/GS)
6	Command instrument system processor -- 70600-01038-101
7	SAS amplifier -- 70901-02908-104 (flight control stability augmentation system)
8	Rate gyro -- TRU-2A/A
9	Amplifier, impedance -- AM-4859A/ARN-89
10	Cargo hook -- FE-7590-145
11	Receiver, radar -- RT-1193/ASN-128 (doppler navigation receiver)
12	Barometric altimeter -- AAU-31/A-1
13	Barometric altimeter -- AAU-32A
14	Receiver/transmitter -- RT-1300/ARC-186 (VHF-AM and/or FM radio)
15	UHF-AM radio set -- RT-1518/ARC-164
16	Interphone control -- C6533/ARC (aircraft intercom control)
17	Receiver/transmitter -- RT-1115D/APN-209 (radar altimeter)
18	Indicator altimeter -- ID-1917C/APN-209 (radar altimeter)
19	Control radio set -- C-7392A/ARN-89 (automatic direction finder)
20	Comparator signal data -- CM-482/ARC-186 (comparator for ARC-186)
21	Receiver/transmitter -- RT-1296A/APX-100 (transponder with IFF)
22	Computer display unit -- CP-1252/ASN-128 (doppler navigation system)
23	Compass set controller -- C-8021E/ASN75
24	Magnetic compass - standby -- MS-17983-4

3.2.3 Battery life evaluation

Battery Life Evaluation Report Form

Nomenclature: PacVac portable suction

Manufacturer: SSCOR

Model number: 10002

Serial number: 105

Military item number: None

Options installed: None

Manufacturer battery life specification: 30 to 45 minutes continuous operation.

Specified battery recharge time: 4 to 6 hours to fully charge a depleted battery.

Specified mode of operation under battery power: Suction with no load.

Overall performance: Pass

Measurements: The unit averaged 45 minutes continuous operation.

Comments: None

3.2.4 Human factors evaluation

Human Factors Evaluation Report Form

Nomenclature: PacVac portable suction

Manufacturer: SSCOR

Model number: 10002

Serial number: 105

Military item number: None

Options installed: None

Date of test: 7 Sep 88

Item configuration during test: Item prepared for operation,
sitting on a counter top.

Checklist for HFE

RESULTS

VISUAL DISPLAYS:

n/a

display type, format, content
location of displays
indicator lights
scalar displays
color coding
legends and labels
cathode ray tubes
counters
flags, go-no-go, center-null indicators

Comments: None

CONTROLS:

Satisfactory

location
characteristics of controls
labeling
control - display relationships

Comments: None

TIME REQUIRED TO PREPARE FOR OPERATION (list in comment)

Comments: Less than 2 minutes.

MAINTAINABILITY:

Satisfactory

- component location
- component characteristics
- rests and stands
- covers, cases, access doors
- handles
- lubrication
- component mounting
- cord storage provisions
- external accessibility
- internal accessibility
- list special tools required
- list realistic inspection requirements
- list realistic inspection intervals

Comments: None

CONDUCTORS:

Satisfactory

- binding and securing
- length
- protection
- routing
- conductor coding
- fabrication
- connectors

Comments: None

FASTENERS:

Satisfactory

- access through inspection panel covers
- enclosure fasteners
- device mounting bolts and fasteners

Comments: None

TEST POINTS:

n/a

- general
- location and mounting
- test point labeling and coding

Comments: None

TEST EQUIPMENT:

Satisfactory

general
equipment self-test
indicators (list in comments)
controls
positive indication of proper operation

Comments: None

FUSES AND CIRCUIT BREAKERS:

Satisfactory

external accessibility
easy replacement or reset by operator

Comments: None

LABELS AND CODING:

Satisfactory

placed above controls and displays
near or on the items they identify
not obscured by other equipment components
describe the function of the items they identify
readable from normal operating distance
conspicuous placards adjacent to hazardous items

Comments: None

SAFETY:

Satisfactory

manual
materials
fire and explosive protection
operator protection from mechanical hazards
patient protection from mechanical hazards
electrical safety (operator and patient)

Comments: None

3.2.5 Altitude test

Altitude Test Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 26 Aug 88

Item configuration during test: Item turned on in the standby mode, operating on dc (battery) power, sitting on chamber floor.

Performance test criteria: Pump suction 550 mL water in less than 7.5 seconds.

Ambient conditions outside chamber:

Temperature	100°F
Humidity	50% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None (battery)
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Battery charger

IN-TEST DATA

Time of test start: 1420

POSTTEST DATA

Posttest performance check (complete check of item and accessories):

Time of test end: 1545

Item functional (based on performance test criteria): Yes

Deviation from pretest : None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.6 Vibration test

Vibration Test Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 26 Aug 88

Item configuration during test: Item strapped down on vibration table fixture; ac and dc operation.

Performance test criteria: Pump suction 550 mL water in less than 7.5 seconds.

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None (battery)
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Battery charger

IN-TEST DATA

Data and performance checks during test:

Times of test start: 0827

Time at first check:

X: 0835

Y: 0950

Z: 1245

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Time at second check:

X: 0925

Y: 1045

Z: 1330

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

POSTTEST DATA

Posttest performance check (complete check of item and accessories):

Item functional (based on performance test criteria): Yes

Item intact: Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): High frequency vibrations, probably resulting from loose batteries inside the unit, caused momentary shutdown of vibration table circuits.

Comments on other data: None

3.2.7 High temperature test

High Temperature Test (Equipment Operating) Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 1 Sep 88

Item configuration during test: Unit was sitting on chamber floor, ready for operation.

Performance test criteria: Consistent and accurate displays and measurements.

Ambient conditions outside chamber:

Temperature	23°C
Humidity	55% RH
Barometric Pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	Reservoir of water
list unconnected terminals	None
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.16
distance from floor (meters)	0.0

IN-TEST DATA

Time of test start: 0900

Performance checks during test:

First check:

Time: 1000
Temperature: 48.8°C
Humidity: 16% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria)
Yes, all OK
Deviation from pretest: None

Second check:

Time: 1030
Temperature: 49°C
Humidity: 16% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria)
Yes, all OK
Deviation from pretest: None

Third check:

Time: 1100
Temperature: 49°C
Humidity: 15% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria)
Yes, all OK
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1130
Item functional (based on performance test criteria)
Yes, all OK
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.8 High temperature storage test

High Temperature Test (Equipment in Storage) Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 2 Sep 88

Item configuration during test: Sitting on chamber floor, in storage, not operating.

Performance test criteria: Suctions 550 mL water in less than 7.5 seconds.

Ambient conditions outside chamber:

Temperature	25°C
Humidity	91% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Battery charger
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

Time of test start: 0800

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1500

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.9 Low temperature test

Low Temperature Test (Equipment Operating) Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 1 Sep 88

Item configuration during test: Sitting on chamber floor, ready for operation.

Performance test criteria: Suctions 550 mL water in less than 7.5 seconds.

Ambient conditions outside chamber:

Temperature	23°C
Humidity	50% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	None
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.0
distance from floor (meters)	0.0

Time of test start: 1310

Performance checks during test:

First check:

Time: 1340
Temperature: 0°C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Second check:

Time: 1415
Temperature: 0°C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Third check:

Time: 1445
Temperature: 0°C
Humidity: n/a
Barometric pressure: 1 atm
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1515
Item functional (based on performance test criteria): Yes
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions):

Condensation on unit was allowed to dry before final performance check.

Comments on other data: None

3.2.10 Low temperature storage test

Low Temperature Test (Equipment in Storage) Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 6 Sep 88

Item configuration during test: Sitting on chamber floor, in storage, not operating.

Performance test criteria: Pumps 550 mL water in less than 7.5 seconds.

Ambient conditions outside chamber:

Temperature	25°C
Humidity	50% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Battery charger
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

Time of test start: 0745

POSTTEST DATA

Posttest performance check:
(complete check of item and accessories)

Time of test end: 1445

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: The unit was allowed to return to ambient conditions overnight before final performance check.

3.2.11 Humidity test

Humidity Test Report Form

Nomenclature: PacVac portable suction
Manufacturer: SSCOR
Model number: 10002
Serial number: 105
Military item number: None

Options installed: None

Date of test: 7 Sep 88

Item configuration during test: The unit was sitting on the chamber floor, ready for operation.

Performance test criteria: Pumps 550 mL water in less than 7.5 seconds.

Ambient conditions outside chamber:

Temperature	23°C
Humidity	87% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	None
list connections to simulators	None
list connections to dummy loads	None
list unconnected terminals	Battery charger
distance from north wall (meters)	0.75
distance from south wall (meters)	0.75
distance from east wall (meters)	2.0
distance from west wall (meters)	2.0
distance from ceiling (meters)	2.6
distance from floor (meters)	0.0

IN-TEST DATA

Time of test start: 0800

Performance checks during test:

First check:

Time: 0915
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): All OK
Deviation from pretest: None

Second check:

Time: 1000
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): All OK
Deviation from pretest: None

Third check:

Time: 1045
Temperature: 29.5°C
Humidity: 95% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): All OK
Deviation from pretest: None

Fourth check:

Time: 1130
Temperature: 29.5°C
Humidity: 94% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): All OK
Deviation from pretest: None

Fifth check:

Time: 1215
Temperature: 29.5°C
Humidity: 94% RH
Barometric pressure: 1 atm
Item functional (based on performance test criteria): All OK
Deviation from pretest: None

POSTTEST DATA

Posttest performance check:

(complete check of item and accessories)

Time of test end: 1430

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.12 Electromagnetic characteristics test

Electromagnetic Characteristics Testing
Evaluation of Performance

T & E Item Number: 07

Date: 1 Sep 88

Nomenclature: PacVac portable suction

Manufacturer: SSCOR

Model number: 10002

Serial number: 105

Military item number: n/a

Conducted Emissions Tests

CE01 Testing configuration(s): n/a
 Performance (pass/fail): n/a

Comments: None

CE02 Testing configuration(s): n/a
 Performance (pass/fail): n/a

Comments: None

CE04 Testing configuration(s): n/a
 Performance (pass/fail): n/a

Comments: None

Conducted Susceptibility Tests

CS02 Testing configuration(s): n/a
 Performance (pass/fail): n/a

Comments: None

CS06 Testing configuration(s): n/a
 Performance (pass/fail): n/a

Comments: None

Radiated Emissions Tests

RE02 Testing configuration(s): Operating on wooden
 test stand in the EMC chamber
 Performance (pass/fail): Fail

 Comments: BB failures 0.2 to 24.8 dB over
 specification, 14 kHz to 718.2 MHz.

Radiated Susceptibility Tests

RS03 Testing configuration(s): Operating on the wooden
 test stand in the EMC chamber.

 Performance (pass/fail): Pass

 Comments: Item connected to a pressure meter to
 allow continuous assessment of pump suction.

3.3 CRITERIA, SIGNIFICANT PROBLEMS, AND SUGGESTED IMPROVEMENTS

3.3.1 Criteria

Item		<u>Applicable</u>	
<u>No.</u>	<u>Criteria (Source)</u>	<u>Remarks</u>	<u>subparagraph</u>
1	The physical inventory is conducted solely for investigation and documentation.	n/a	2.1.2.1
2	The PacVac suction unit will display consistent and accurate measurements.	met	2.1.2.2
3	Verify manufacturer's specified full power internal battery life expectancy of 30 to 45 minutes.	met	2.2.2
4	The PacVac suction unit will be rated satisfactory in all major categories of the evaluation. These include: Visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.	par- tially met	2.4.2
5	The PacVac suction unit will display consistent and accurate measurements while exposed to an altitude equivalency of 15,000 feet above sea level.	met	2.5.2
6	The PacVac suction unit will remain operational and display consistent and accurate measurements while exposed to vibrational stresses.	met	2.6.2
7	The PacVac suction unit will display consistent and accurate measurements during the high temperature operation check.	met	2.7.2.1
8	The PacVac suction unit will display consistent and accurate measurements after the high temperature storage.	met	2.7.2.2

9	The PacVac suction unit will display consistent and accurate measurements during the low temperature operation check.	met	2.8.2.1
10	The PacVac suction unit will display consistent and accurate measurements after the low temperature storage.	met	2.8.2.2
11	The PacVac suction unit will display consistent and accurate measurements while exposed to a high humidity.	met	2.9.2
12	The PacVac suction unit will not produce emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.13.	partially met	2.10.2.1
13	The PacVac suction unit will not malfunction when it is subjected to radiated fields as specified in MIL-STD-461A, Notice 4, paragraph 6.20.	met	2.10.2.2
14	The flight surgeon will be able to operate the PacVac suction unit without physical or functional restrictions aboard the aircraft.	not tested	2.11.2.1
15	The PacVac suction unit will not radiate EMI to disrupt or interfere with the other equipment or systems aboard the aircraft.	not tested	2.12.2.2
16	The aircraft will not radiate EMI to disrupt or interfere with the PacVac suction unit.	not tested	2.12.2.3

3.3.2 Significant problems which require corrective action

None

3.3.3 Suggested improvements

None

3.4 REFERENCES

- 3.4.1 Department of Defense. 1971. EMI characteristics, requirements for equipment. Washington, D.C. MIL-STD-461A, Notice 4. February.
- 3.4.2 Department of Defense. 1971. EMI characteristics, measurement of. Washington, D.C. MIL-STD-462, Notice 3. February.
- 3.4.3 Department of Defense. 1983. Environmental test methods and engineering guidelines. Washington, D.C. MIL-STD-810D. July.
- 3.4.4 Department of the Army. 1982. Environmental protection and enhancement. Washington, D.C. Army Regulation 200-1. June.
- 3.4.5 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, D.C. TB 38-750-2. April.
- 3.4.6 Department of Defense. 1985. Standard general requirements for electronic equipment. Washington, D.C. MIL-STD-454K. February.
- 3.4.7 Underwriters Laboratory's, Inc. 1978. Standard for safety, medical and dental equipment. Chicago, Illinois. UL-544.
- 3.4.8 Department of Defense. 1989. Human engineering design criteria for military systems, equipment, and facilities. Washington, D.C. MIL-STD-1472D. March.
- 3.4.9 Association for the Advancement of Medical Instruments. 1988. Human factors engineering guidelines and preferred practices for the design of medical devices. Arlington, Virginia. AAMI-HE-1988. February.
- 3.4.10 Department of the Army. 1978. Operator's manual, UH-60 and EH-60 helicopter, with changes 1-5. Washington, D.C. TM 55-1520-237-10. January.
- 3.4.11 Mitchell, G. W., and Adams, J. E. 1988. Technical test and evaluation of aeromedical equipment. Fort Rucker, AL: U.S. Army Aeromedical Research Laboratory. USAARL Letter Report LR-88-16-1-2.

3.5 ABBREVIATIONS

ac	alternating current
AEST	aeromedical equipment suitability test
atm	atmosphere
AVSCOM	U.S. Army Aviation Systems Command
AWR	airworthiness release
BB	broadband
BPM	beats per minute
CAAF	Cairns Army Airfield
CRT	cathode ray tube
dB	decibel
dc	direct current
EMC	electromagnetic compatibility
EMI	electromagnetic interference
fpm	feet per minute
GFE	government furnished equipment
GHz	gigahertz
Gpk	gravity, peak
G(rms)	gravity (root mean square)
Hz	hertz
IAW	in accordance with
IGE	in-ground effect
ITOP	in-flight test operating procedure
kHz	kilohertz
KIAS	knots indicated airspeed
LCD	liquid crystal display
LED	light emitting diode
MEDEVAC	medical evacuation
MHz	megahertz
MIL-STD	military standard
mL	milliliter
mm	millimeter
mmHg	millimeters of Mercury
m.s.l.	mean sea level
NB	narrowband
NBC	nuclear, biological and chemical
NiCad	nickel cadmium
NVG	night vision goggle

PacVac	SSCOR PacVac portable suction pump
RAM	random access memory
RF	radio frequency
RH	relative humidity
ROM	read only memory
TB	technical bulletin
TFT	technical feasibility testing
T & E	test and evaluation
UES	Universal Energy Systems, Inc.
USAARL	U.S. Army Aeromedical Research Laboratory
V/m	volts per meter

3.6 LIST OF MANUFACTURERS

- 3.6.1 SSCOR
 1210 El Vago Street
 La Canada, CA 91011
- 3.6.2 Sikorsky Aircraft
 6900 Main Street
 Stratford, CT 06601
- 3.6.3 Tenney Engineering, Inc.
 1090 Springfield Road
 Post Office Box 3142
 Union, NJ 07083
- 3.6.4 Unholtz-Dickey Corporation
 6 Brookside Drive
 Wallingford, CT 06492

3.7 DISTRIBUTION LIST

Commander, U.S. Army Natick Research,
Development and Evaluation Center
ATTN: STRNC-MIL (Documents
Librarian)
Natick, MA 01760-5040

Commander
U.S. Army Aviation Systems Command
ATTN: AMSAV-ECU
4300 Goodfellow Bouvelard
St. Louis, MO 63120-1790

Commander/Director
U.S. Army Combat Surveillance
and Target Acquisition Lab
ATTN: DELCS-D
Fort Monmouth, NJ 07703-5304

Commander
10th Medical Laboratory
ATTN: Audiologist
APO New York 09180

Naval Air Development Center
Technical Information Division
Technical Support Detachment
Warminster, PA 18974

Commanding Officer, Naval Medical
Research and Development Command
National Naval Medical Center
Bethesda, MD 20814-5044

Deputy Director, Defense Research
and Engineering
ATTN: Military Assistant
for Medical and Life Sciences
Washington, DC 20301-3080

Commander, U.S. Army Research
Institute of Environmental Medicine
Natick, MA 01760

U.S. Army Avionics Research
and Development Activity
ATTN: SAVAA-P-TP
Fort Monmouth, NJ 07703-5401

U.S. Army Communications-Electronics
Command
ATTN: AMSEL-RD-ESA-D
Fort Monmouth, NJ 07703

Library
Naval Submarine Medical Research Lab
Box 900, Naval Sub Base
Groton, CT 06349-5900

Commander
Man-Machine Integration System
Code 602
Naval Air Development Center
Warminster, PA 18974

Commander
Naval Air Development Center
ATTN: Code 602-B (Mr. Brindle)
Warminster, PA 18974

Commanding Officer
Harry G. Armstrong Aerospace
Medical Research Laboratory
Wright-Patterson
Air Force Base, OH 45433

Director
Army Audiology and Speech Center
Walter Reed Army Medical Center
Washington, DC 20307-5001

Commander, U.S. Army Institute
of Dental Research
ATTN: Jean A. Setterstrom, Ph. D.
Walter Reed Army Medical Center
Washington, DC 20307-5300

Naval Air Systems Command
Technical Air Library 950D
Room 278, Jefferson Plaza II
Department of the Navy
Washington, DC 20361

Naval Research Laboratory Library
Shock and Vibration
Information Center, Code 5804
Washington, DC 20375

Director, U.S. Army Human
Engineering Laboratory
ATTN: Technical Library
Aberdeen Proving Ground, MD 21005

Commander, U.S. Army Test
and Evaluation Command
ATTN: AMSTE-AD-H
Aberdeen Proving Ground, MD 21005

Director
U.S. Army Ballistic
Research Laboratory
ATTN: DRXBR-OD-ST Tech Reports
Aberdeen Proving Ground, MD 21005

Commander
U.S. Army Medical Research
Institute of Chemical Defense
ATTN: SGRD-UV-AO
Aberdeen Proving Ground,
MD 21010-5425

Commander, U.S. Army Medical
Research and Development Command
ATTN: SGRD-RMS (Ms. Madigan)
Fort Detrick, Frederick, MD 21702-5012

Director
Walter Reed Army Institute of Research
Washington, DC 20307-5100

HQ DA (DASG-PSP-O)
5109 Leesburg Pike
Falls Church, VA 22041-3258

Naval Research Laboratory
Library Code 1433
Washington, DC 20375

Harry Diamond Laboratories
ATTN: Technical Information Branch
2800 Powder Mill Road
Adelphi, MD 20783-1197

U.S. Army Materiel Systems
Analysis Agency
ATTN: AMXSY-PA (Reports Processing)
Aberdeen Proving Ground
MD 21005-5071

U.S. Army Ordnance Center
and School Library
Simpson Hall, Building 3071
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental
Hygiene Agency
Building E2100
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research
and Development Center
Aberdeen Proving Ground, MD
21010-5423

Commander
U.S. Army Medical Research
Institute of Infectious Disease
SGRD-UIZ-C
Fort Detrick, Frederick, MD 21702

Director, Biological
Sciences Division
Office of Naval Research
600 North Quincy Street
Arlington, VA 22217

Commander
U.S. Army Materiel Command
ATTN: AMCDE-XS
5001 Eisenhower Avenue
Alexandria, VA 22333

Commandant
U.S. Army Aviation
Logistics School ATTN: ATSQ-TDN
Fort Eustis, VA 23604

Headquarters (ATMD)
U.S. Army Training
and Doctrine Command
Fort Monroe, VA 23651

Structures Laboratory Library
USARTL-AVSCOM
NASA Langley Research Center
Mail Stop 266
Hampton, VA 23665

Naval Aerospace Medical
Institute Library
Building 1953, Code 03L
Pensacola, FL 32508-5600

Command Surgeon
HQ USCENTCOM (CCSG)
U.S. Central Command
MacDill Air Force Base FL 33608

Air University Library
(AUL/LSE)
Maxwell Air Force Base, AL 36112

U.S. Air Force Institute
of Technology (AFIT/LDEE)
Building 640, Area B
Wright-Patterson
Air Force Base, OH 45433

Henry L. Taylor
Director, Institute of Aviation
University of Illinois-Willard Airport
Savoy, IL 61874

COL Craig L. Urbauer, Chief
Office of Army Surgeon General
National Guard Bureau
Washington, DC 50310-2500

Commander
U.S. Army Aviation Systems Command
ATTN: SGRD-UAX-AL (MAJ Gillette)
4300 Goodfellow Blvd., Building 105
St. Louis, MO 63120

U.S. Army Aviation Systems Command
Library and Information Center Branch
ATTN: AMSAV-DIL
4300 Goodfellow Boulevard
St. Louis, MO 63120

Federal Aviation Administration
Civil Aeromedical Institute
Library AAM-400A
P.O. Box 25082
Oklahoma City, OK 73125

Commander
U.S. Army Academy
of Health Sciences
ATTN: Library
Fort Sam Houston, TX 78234

Commander
U.S. Army Institute of Surgical Research
ATTN: SGRD-USM (Jan Duke)
Fort Sam Houston, TX 78234-6200

AAMRL/HEX
Wright-Patterson
Air Force Base, OH 45433

John A. Dellinger,
Southwest Research Institute
P. O. Box 28510
San Antonio, TX 78284

Product Manager
Aviation Life Support Equipment
ATTN: AMCPM-ALSE
4300 Goodfellow Boulevard
St. Louis, MO 63120-1798

Commander
U.S. Army Aviation
Systems Command
ATTN: AMSAV-ED
4300 Goodfellow Boulevard
St. Louis, MO 63120

Commanding Officer
Naval Biodynamics Laboratory
P.O. Box 24907
New Orleans, LA 70189-0407

Assistant Commandant
U.S. Army Field Artillery School
ATTN: Morris Swott Technical Library
Fort Sill, OK 73503-0312

Commander
U.S. Army Health Services Command
ATTN: HSOP-SO
Fort Sam Houston, TX 78234-6000

Director of Professional Services
HQ USAF/SGDT
Bolling Air Force Base, DC 20332-6188

U.S. Army Dugway Proving Ground
Technical Library, Building 5330
Dugway, UT 84022

U.S. Army Yuma Proving Ground
Technical Library
Yuma, AZ 85364

AFFTC Technical Library
6510 TW/TSTL
Edwards Air Force Base,
CA 93523-5000

Commander
Code 3431
Naval Weapons Center
China Lake, CA 93555

Aeromechanics Laboratory
U.S. Army Research and Technical Labs
Ames Research Center, M/S 215-1
Moffett Field, CA 94035

Sixth U.S. Army
ATTN: SMA
Presidio of San Francisco, CA 94129

Commander
U.S. Army Aeromedical Center
Fort Rucker, AL 36362

U.S. Air Force School
of Aerospace Medicine
Strughold Aeromedical Library Technical
Reports Section (TSKD)
Brooks Air Force Base, TX 78235-5301

Dr. Diane Damos
Department of Human Factors
ISSM, USC
Los Angeles, CA 90089-0021

U.S. Army White Sands
Missile Range
ATTN: STEWS-IM-ST
White Sands Missile Range, NM 88002

U.S. Army Aviation Engineering
Flight Activity
ATTN: SAVTE-M (Tech Lib) Stop 217
Edwards Air Force Base, CA 93523-5000

Ms. Sandra G. Hart
Ames Research Center
MS 262-3
Moffett Field, CA 94035

Commander, Letterman Army Institute
of Research
ATTN: Medical Research Library
Presidio of San Francisco, CA 94129

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200

Commander
U.S. Army Medical Materiel
Development Activity
Fort Detrick, Frederick, MD 21702-5009

Commander
U.S. Army Aviation Center
Directorate of Combat Developments
Building 507
Fort Rucker, AL 36362

U. S. Army Research Institute
Aviation R&D Activity
ATTN: PERI-IR
Fort Rucker, AL 36362

Commander
U.S. Army Safety Center
Fort Rucker, AL 36362

U.S. Army Aircraft Development
Test Activity
ATTN: STEBG-MP-P
Cairns Army Air Field
Fort Rucker, AL 36362

Commander U.S. Army Medical Research
and Development Command
ATTN: SGRD-PLC (COL Sedge)
Fort Detrick, Frederick, MD 21702

MAJ John Wilson
TRADOC Aviation LO
Embassy of the United States
APO New York 09777

Netherlands Army Liaison Office
Building 602
Fort Rucker, AL 36362

British Army Liaison Office
Building 602
Fort Rucker, AL 36362

Italian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Directorate of Training Development
Building 502
Fort Rucker, AL 36362

Chief
USAHEL/USAAVNC Field Office
P. O. Box 716
Fort Rucker, AL 36362-5349

Commander U.S. Army Aviation Center
and Fort Rucker
ATTN: ATZQ-CG
Fort Rucker, AL 36362

Commander/President
TEXCOM Aviation Board
Cairns Army Air Field
Fort Rucker, AL 36362

MAJ Terry Newman
Canadian Army Liaison Office
Building 602
Fort Rucker, AL 36362

German Army Liaison Office
Building 602
Fort Rucker, AL 36362

LTC Patrice Cottebrune
French Army Liaison Office
USAAVNC (Building 602)
Fort Rucker, AL 36362-5021

Brazilian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Australian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Dr. Garrison Rapmund
6 Burning Tree Court
Bethesda, MD 20817

Commandant Royal Air Force
Institute of Aviation Medicine
Farnborough Hants UK GU14 6SZ
Dr. A. Kornfield, President
Biosearch Company
3016 Revere Road
Drexel Hill, PA 29026

Commander
U.S. Army Biomedical Research
and Development Laboratory
ATTN: SGRD-UBZ-I
Fort Detrick, Frederick, MD 21702

Defense Technical Information Center
Cameron Station
Alexandra, VA 22313

Commander, U.S. Army Foreign Science
and Technology Center
AIFRTA (Davis)
220 7th Street, NE
Charlottesville, VA 22901-5396

Director,
Applied Technology Laboratory
USARTL-AVSCOM
ATTN: Library, Building 401
Fort Eustis, VA 23604

U.S. Army Training
and Doctrine Command
ATTN: Surgeon
Fort Monroe, VA 23651-5000

Aviation Medicine Clinic
TMC #22, SAAF
Fort Bragg, NC 28305

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

Commander, U.S. Army Missile
Command
Redstone Scientific Information Center
ATTN: AMSMI-RD-CS-R/ILL
Documents Redstone Arsenal, AL 35898

U.S. Army Research and Technology
Laboratories (AVSCOM)
Propulsion Laboratory MS 302-2
NASA Lewis Research Center
Cleveland, OH 44135

Dr. H. Dix Christensen
Bio-Medical Science Building, Room 753
Post Office Box 26901
Oklahoma City, OK 73190

Dr. Christine Schlichting
Behavioral Sciences Department
Box 900, NAVUBASE NLON
Groton, CT 06349-5900

Commandant
Academy of Health Sciences
ATTN: HSHA-COM (LTC Huether)
Fort Sam Houston, TX 78234

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

COL Eugene S. Channing, O.D.
Brooke Army Medical Center
ATTN: HSHE-EAH-O
Fort Sam Houston, TX 78234-6200