

1

USAARL Report No. 93-4

AD-A259 924



An Automated Method for Determining Mass Properties

By

Michael B. Deavers

and

R. Joseph McEntire

DTIC
ELECTE
JAN 26 1993
S C D

93-01372



3308

Biodynamics Research Division

December 1992

Approved for public release; distribution unlimited.

United States Army Aeromedical Research Laboratory
Fort Rucker, Alabama 36362-5292

Notice

Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.

Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

Disclaimer



The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:



JOHN V. BARSON
LTC, MC, SFS
Director, Biodynamics
Research Division

Released for publication:

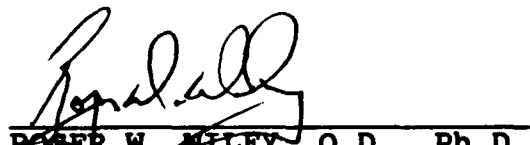
ROGER W. MILEY, O.D., Ph.D.
Chairman, Scientific
Review Committee
DAVID H. KARNEY
Colonel, MC, SFS
Commanding

REPORT DOCUMENTATION PAGE

Form Approved
 OMB No. 0704-0188

1a. REPORT SECURITY CLASSIFICATION Unclassified		1b. RESTRICTIVE MARKINGS	
2a. SECURITY CLASSIFICATION AUTHORITY		3. DISTRIBUTION / AVAILABILITY OF REPORT Approved for public release, distribution unlimited	
2b. DECLASSIFICATION / DOWNGRADING SCHEDULE			
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 93-4		5. MONITORING ORGANIZATION REPORT NUMBER(S)	
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory	6b. OFFICE SYMBOL (if applicable) SGRD-UAD-IE	7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command	
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 77 Fort Rucker, AL 36362-5292		7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012	
8a. NAME OF FUNDING / SPONSORING ORGANIZATION	8b. OFFICE SYMBOL (if applicable)	9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER	
8c. ADDRESS (City, State, and ZIP Code)		10. SOURCE OF FUNDING NUMBERS	
		PROGRAM ELEMENT NO. 0602787A	PROJECT NO. 3M162787A878
		TASK NO. AG	WORK UNIT ACCESSION NO. 138
11. TITLE (Include Security Classification) An automated method for determining mass properties			
12. PERSONAL AUTHOR(S) Michael B. Deavers and B. Joseph McEntire			
13a. TYPE OF REPORT Final	13b. TIME COVERED FROM _____ TO _____	14. DATE OF REPORT (Year, Month, Day) 1992 December	15. PAGE COUNT 22
16. SUPPLEMENTARY NOTATION			
17. COSATI CODES		18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number)	
FIELD	GROUP	center-of-mass, mass moments of inertia, impacts, man-mounted equipment	
20	11		
23	02		
19. ABSTRACT (Continue on reverse if necessary and identify by block number) The U.S. Army Aeromedical Research Laboratory performs mass properties testing with the KSR330-60 Mass Properties Instrument (MPI). For man-mounted equipment, it is important to define and measure the mass properties in order to accurately develop math models, understand human performance impacts, and to conduct comparative evaluations. Determining mass properties requires a thorough theoretical understanding of the center-of-mass and mass moments of inertia. This report provides the operating theory and procedures to measure mass properties with the MPI.			
20. DISTRIBUTION / AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS		21. ABSTRACT SECURITY CLASSIFICATION Unclassified	
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center		22b. TELEPHONE (Include Area Code) (205) 255-6907	22c. OFFICE SYMBOL SGRD-UAX-SI

Acknowledgments

The authors wish to thank Mr. Troy Bowman for recommendations towards the operation of the MPI.

DTIC QUALITY INSPECTED 8

Accession For	
NTIS GRA&I	<input checked="" type="checkbox"/>
DTIC TAB	<input type="checkbox"/>
Unannounced	<input type="checkbox"/>
Justification	
By	
Distribution/	
Availability Codes	
Dist	Avail and/or Special
A-1	

This page intentionally left blank.

Table of Contents

List of figures.	2
Introduction	3
Equipment.	3
Theory	3
Center-of-mass.	3
Moments of inertia.	4
Mass moments of inertia.	4
Mass products of inertia	6
Principal moments and axes of inertia.	9
Procedure of operation	11
Calibration	12
Measurements.	14
Step 1: Tare measurement.	14
Step 2: Part measurement.	15
Step 3: Calculations.	15
Conclusions.	15
References	17
Appendix A: Accuracy of the MPI	18
Appendix B: Determination of torsional pendulum physical constant	19
Appendix C: KSR main menu	20
Appendix D: Mounting of test fixture.	21
Appendix E: Manufacturers' list	22

List of figures

Figure		Page
1	KSR330-60 Mass Properties Instrument setup	4
2	Moments produced during CM testing	5
3	Inverted torsional pendulum.	6
4	Differential mass with respect to the CM	7
5	Arbitrary axis with respect to the coordinate axis.	8
6	Direction cosines for principal axis 1	10
7	Relation between coordinate and principal axes	11
8	Calibration beam and weights	13
9	MPI setup for CM calibration	13
D-1	Mounting pattern of test platform.	21

Introduction

Mass property evaluations are an integral part of the engineering analysis performed by the Biodynamics Research Division at the U.S. Army Aeromedical Research Laboratory (USAARL). For man-mounted equipment, it is important to define and measure the mass properties in order to accurately develop math models, understand human performance impacts, and to conduct comparative evaluations. To measure these properties, the Space Electronics KSR330-60 Mass Properties Instrument (MPI) was selected. The purpose of this paper is to describe the MPI's operating theory and the procedure to measure mass properties with the MPI.

Equipment

The KSR330-60 MPI was designed and produced by Space Electronics, Inc. The MPI is capable of testing a 330-lb load. The maximum full-scale moment allowed due to the offset of the part's CM is 60 in-lb. The accuracies of the MPI are listed in Appendix A. The MPI requires a continuous source of pressurized (70 to 95 PSI) nitrogen or dry clean instrument air. A JUN-AIR 2000 oil-less air compressor with an absorption dryer was selected to supply instrument air to the MPI. The MPI operation is fully automated with a personal computer and a dot matrix printer. The MPI system setup is shown in Figure 1.

Theory

The MPI is designed to obtain accurate measurements of the center-of-mass (CM) and the mass moment of inertia (MOI). The CM is defined as the point in a body at which its entire mass may be assumed to be concentrated. The MOI of a body is defined as a measure of its resistance to rotational acceleration.

Center-of-mass

For each orthogonal axis of a test part, its CM is calculated by measuring the force required to balance the test part on the test platform's pivot axis. This pivot axis serves as the fulcrum for the test platform, which is suspended by a rotary gas bearing. The test part's mass and the unknown distance between the test part's CM location and the test platform's pivot axis results in a moment applied about the pivot axis (Figure 2). The MPI measures the force required to balance

* See manufacturer's list.



Figure 1. KSR330-60 Mass Properties Instrument setup.

this moment caused by the test part. This balancing force (F_1) is sensed by a force transducer located at a known distance (d_1) from the fulcrum. From the moment equation,

$$M_1 = F_1 d_1,$$

the moment (M_1) produced by the transducer relative to the pivot axis can be determined. In order to balance the test platform, this moment must equal the moment resulting from the CM of the test part. Since the force applied by the test part (F_2) is known from its mass, the unknown distance from the fulcrum to the CM of the test part (d_2) is easily calculated from the following equation:

$$F_1 d_1 = F_2 d_2$$

$$\text{or } d_2 = (F_1 d_1) / F_2.$$

Moments of inertia

Mass moments of inertia

The MOI calculations are based on the period of rotational oscillation of the test platform, which is configured as an

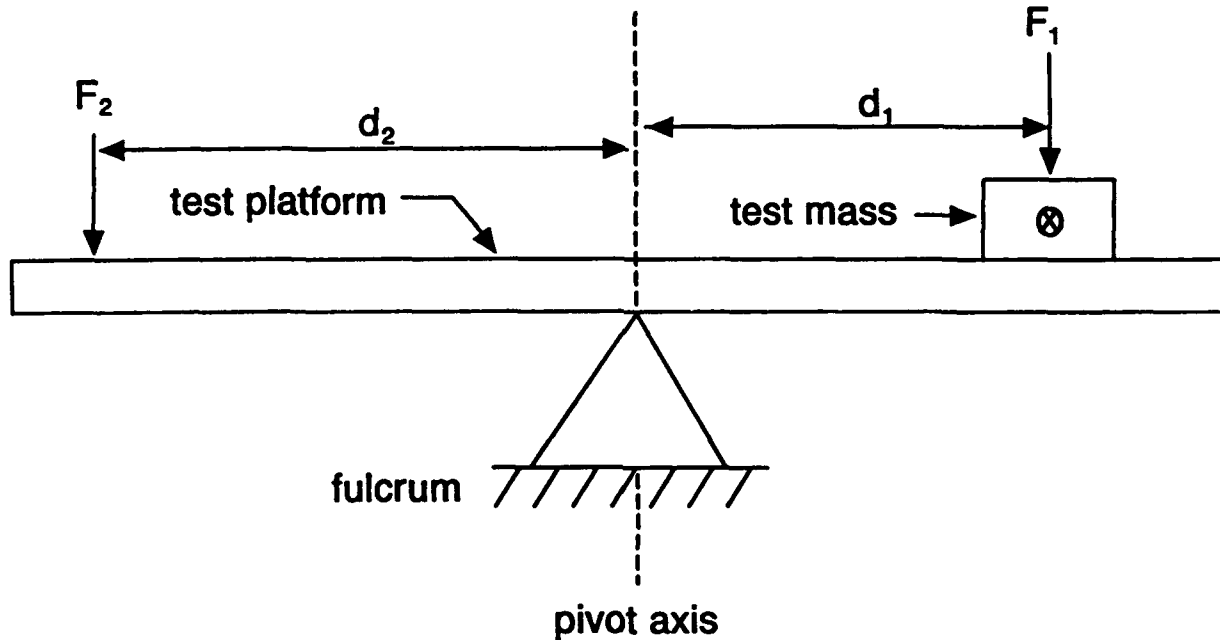


Figure 2. Moments produced during CM testing.

inverted torsional pendulum (Figure 3). To begin the oscillation, the platform is rotated slightly and released. The test platform's MOI is related directly to the oscillation period (T) by the following equation:

$$I = CT^2,$$

where C is a physical constant dependent on the torsional pendulum configuration. This constant is found experimentally as shown in Appendix B.

To determine the MOI of the test part, a tare period (T_i) must first be measured for the platform and any test fixture necessary to support the test part. The test part then is added to the system and another period (T_f) is acquired. Using the calibration constant, C, the test part MOI then is calculated:

$$I = C(T_f^2 - T_i^2).$$

This measurement and calculation provide the part's MOI about a single axis. This procedure must be repeated for the part's two remaining axes. These moments, I_{xx} , I_{yy} , and I_{zz} , are used to determine the part's principal moments of inertia and principal axes of inertia. The MOIs are mathematically represented by:

$$I_{xx} = \int (y^2 + z^2) dm,$$

$$I_{yy} = \int (z^2 + x^2) dm,$$

$$I_{zz} = \int (x^2 + y^2) dm.$$

where x , y , and z are the distances from the basic coordinate axes of the CM to the differential mass (dm) (Figure 4).

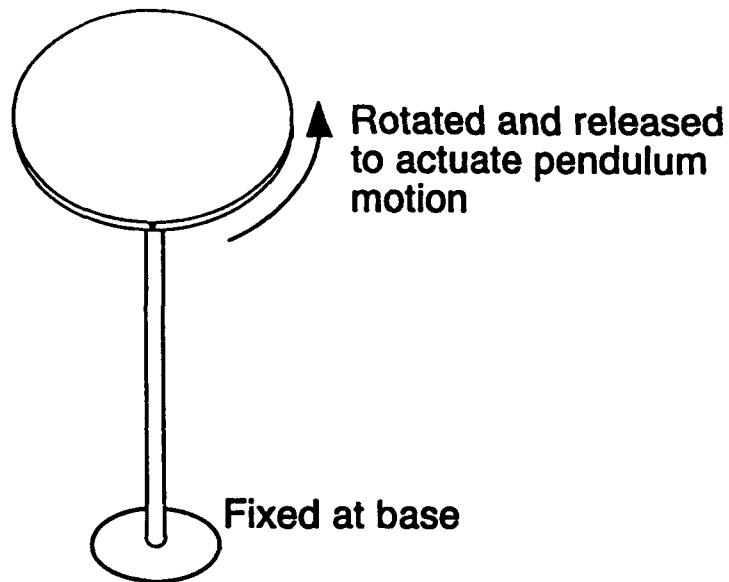


Figure 3. Inverted torsional pendulum.

Mass products of inertia

Along with the MOIs, the products of inertia are required in order to calculate the principal moments of inertia and the principal axes of inertia. The products of inertia, P_{xy} , P_{xz} , and P_{yz} , are mathematically defined by the following integrals:

$$P_{xy} = P_{yx} = \int xy dm$$

$$P_{xz} = P_{zx} = \int xz dm$$

$$P_{yz} = P_{zy} = \int yz dm$$

where x , y , and z are the distances from the basic coordinate axes of the CM to the differential mass (dm) (Figure 4).

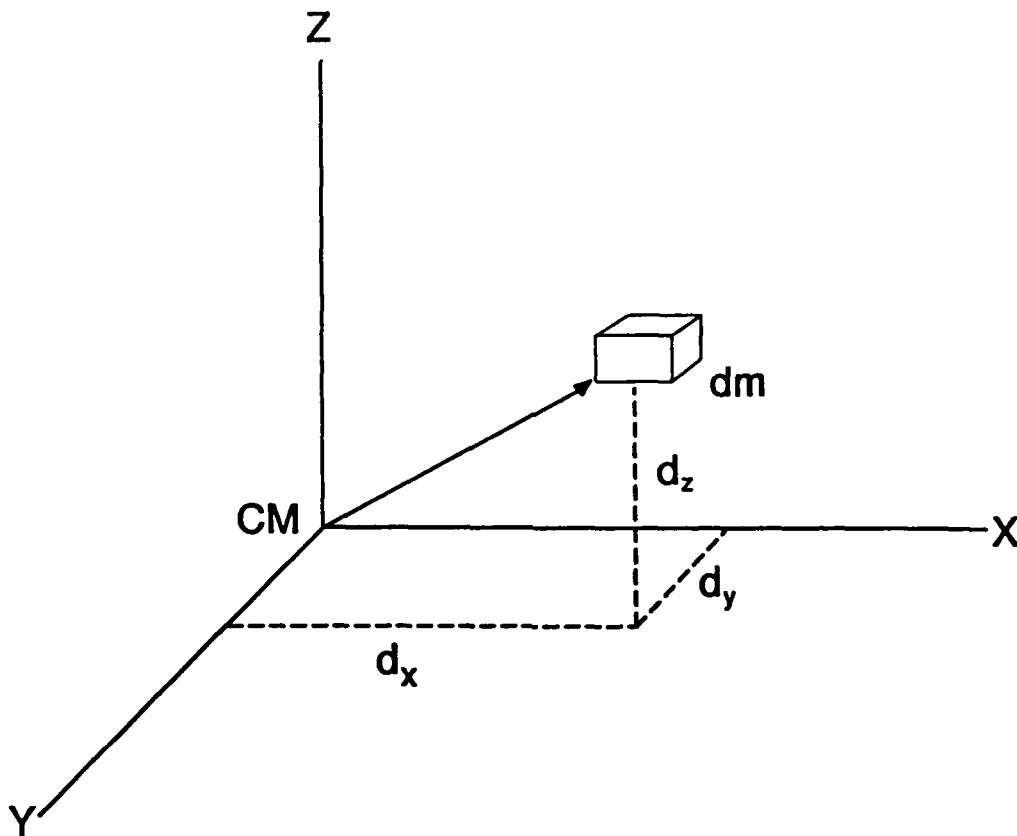


Figure 4. Differential mass with respect to the CM.

To experimentally determine P_{xy} , P_{xz} , and P_{yz} with the MPI, perform the following steps for the x-y plane, the x-z plane, and the y-z plane.

- a. Record the MOIs (I_{xx} , I_{yy} , and I_{zz}) that were determined during the basic MOI testing.
- b. Rotate the test part to an arbitrary axis (remaining within the x-y, y-z, or x-z planes) θ° from the coordinate axis (Figure 5). In order to simplify the calculations, a 45° rotation is suggested.
- c. Run an MOI test with the test part at this orientation. The resulting MOI values (I_{xy} , I_{xz} , or I_{yz}) provide the MOI of the test part with respect to a new coordinate system with axes within the x-y, x-z, and y-z planes of the basic coordinate system.

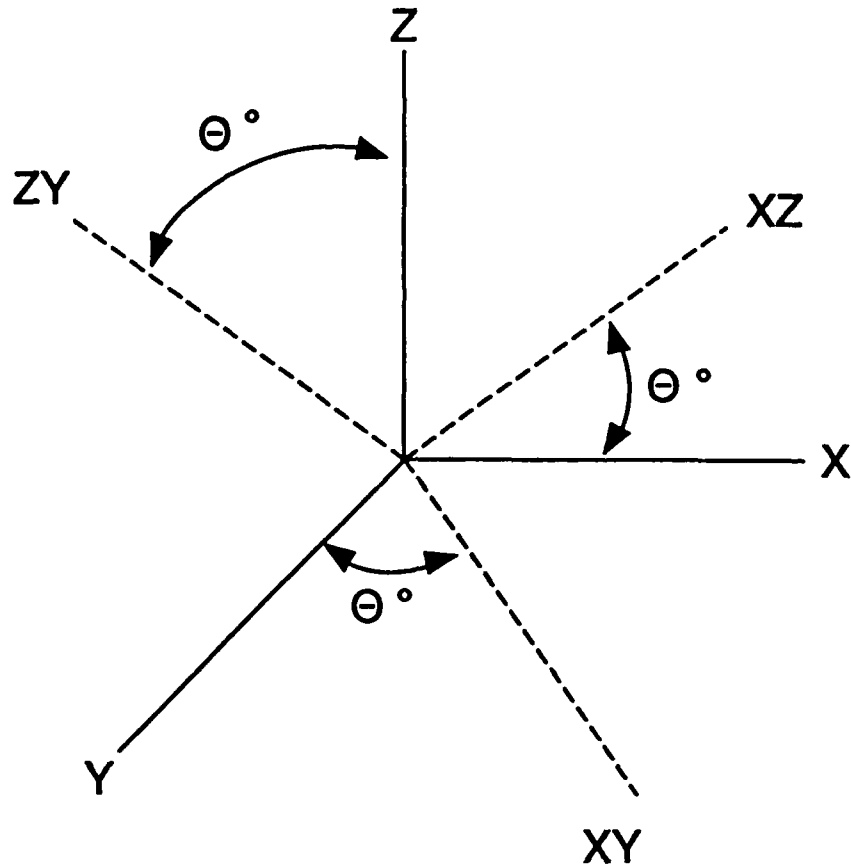


Figure 5. Arbitrary axis with respect to the coordinate axis.

d. Having calculated I_{xx} , I_{yy} , I_{zz} , I_{xy} , I_{xz} , and I_{yz} , the following transformation equations provide a means of solving for the products of inertia (Beer and Johnston, 396-397).

$$\begin{aligned}
 I'_{xy} &= \frac{I_{xx} + I_{yy}}{2} + \frac{I_{xx} - I_{yy}}{2} \cos(2\theta) - P_{xy} \sin(2\theta) \\
 I'_{xz} &= \frac{I_{xx} + I_{zz}}{2} + \frac{I_{xx} - I_{zz}}{2} \cos(2\theta) - P_{xz} \sin(2\theta) \\
 I'_{yz} &= \frac{I_{yy} + I_{zz}}{2} + \frac{I_{yy} - I_{zz}}{2} \cos(2\theta) - P_{yz} \sin(2\theta)
 \end{aligned}$$

The calculations are simplified by using 45° as the rotation angle ($\cos(2\theta)=0$ and $\sin(2\theta)=1$). After rearranging terms and reducing, P_{xy} , P_{xz} , and P_{yz} are solved by substituting the appropriate values into the following equations:

$$P_{xy} = \frac{I_{xx} + I_{yy}}{2} - I_{xy}$$

$$P_{xz} = \frac{I_{xx} + I_{zz}}{2} - I_{xz}$$

$$P_{yz} = \frac{I_{yy} + I_{zz}}{2} - I_{yz}$$

Principal moments and axes of inertia

For every three-dimensional object, there is a maximum, intermediate, and minimum MOI about the part's CM. These three moments are known as the principal MOIs. The axes about which these moments act are called principal axes of inertia. These axes define a unique coordinate system about which the products of inertia are equal to zero.

The principal MOIs are calculated by substituting the basic MOI and products of inertia into the following determinant equation:

$$\begin{bmatrix} I_{xx} - I & -P_{xy} & -P_{xz} \\ -P_{yx} & I_{yy} - I & -P_{yz} \\ -P_{zx} & -P_{zy} & I_{zz} - I \end{bmatrix} = 0.$$

The solution to this determinant is a cubic equation, whose roots yield three values for I. These values, I_1 , I_2 , and I_3 , yield the maximum, intermediate, and minimum values, respectively, of the test part's MOI.

To determine the orientation of the principal axes, the direction cosines (a, b, and c) for each axis must be determined. Direction cosines are the cosines of the angles between the principal axes (1, 2, and 3) and the basic reference axes (x, y, and z) (Figure 6). The direction cosines are calculated by substituting the basic MOIs, the products of inertia, and the individual principal MOIs (I_1 , I_2 , and I_3) into the equations which follow. This substitution yields three equations with three unknowns. The calculations must be performed three times, once for each of the principal moments of inertia.

$$(I_{xx} - I)a - P_{xy}b - P_{xz}c = 0$$

$$-P_{yx}a + (I_{yy} - I)b - P_{yz}c = 0$$

$$-P_{zx}a - P_{zy}b + (I_{zz} - I)c = 0.$$

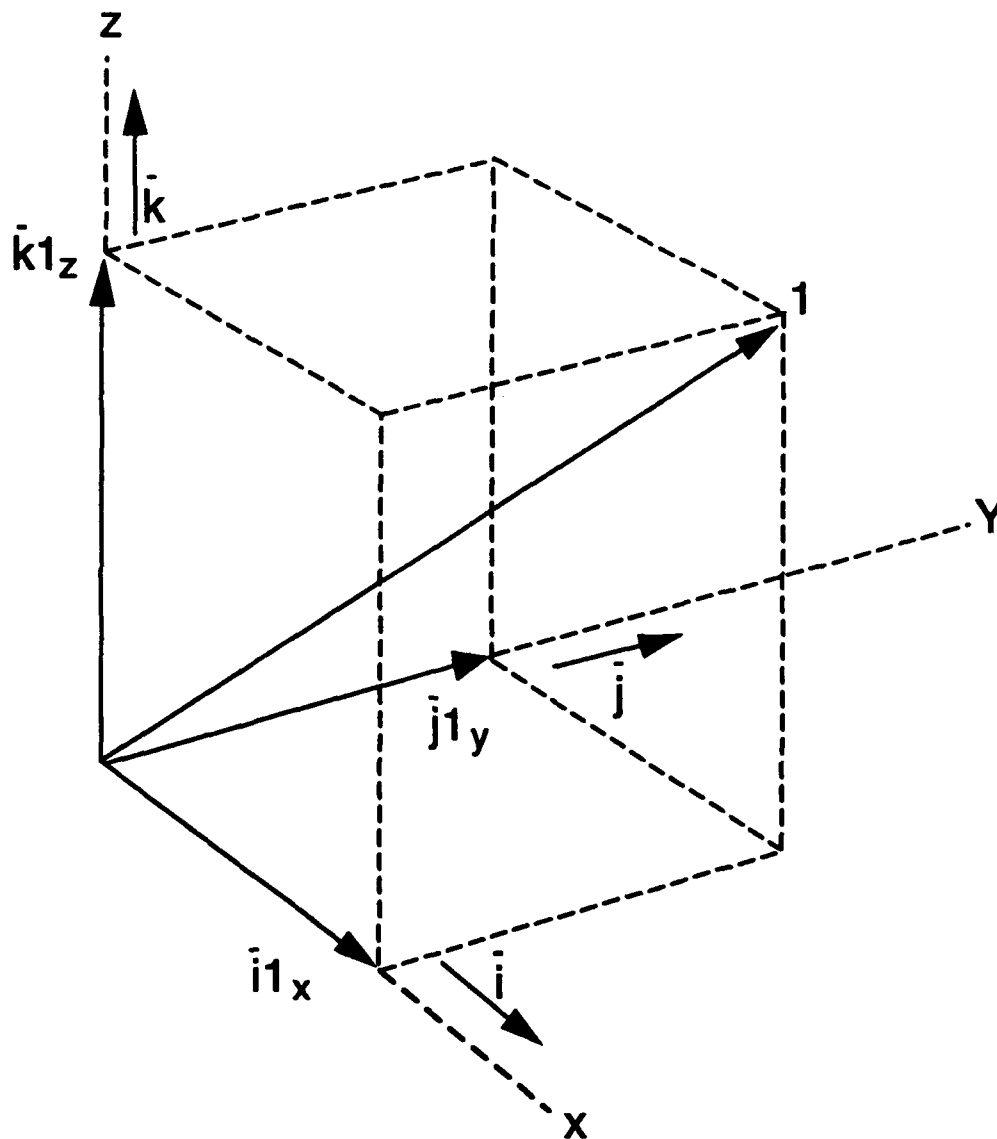


Figure 6. Direction cosines for principal axis 1.

The three calculations result in three sets of direction cosines:

- (1) $a_1, b_1, c_1,$
- (2) $a_2, b_2, c_2,$
- (3) $a_3, b_3, c_3.$

As a check, substitute each set of the direction cosines into the following equation:

$$a^2 + b^2 + c^2 = 1.$$

Each set of direction cosines defines a principal axis (1, 2, or 3) with respect to the basic coordinate axes. The relation between the principal axes and the coordinate axes is shown below in Figure 7 (Meriam and Kraige, 589-606).

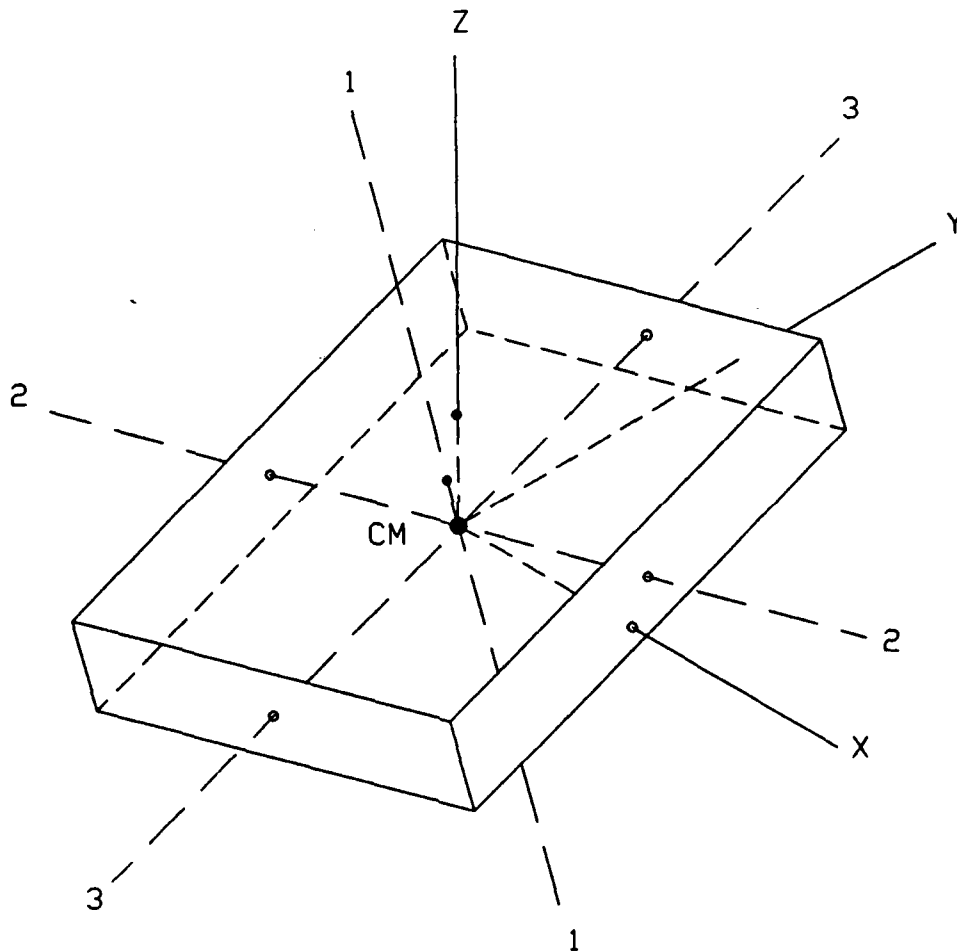


Figure 7. Relation between coordinate and principal axes.

Procedure of operation

The procedures to measure the CM and MOI of test parts are relatively simple, but must be carefully followed in order to ensure accuracy. Before operating the MPI, check for correct setup of the entire system as detailed in the MPI Model KSR330-60

instruction manual. The initial start up of the system is as follows:

1. Remove all objects from the MPI test platform.
2. Turn on the MPI, the gas bearing pressure, the air compressor, and the computer system. Allow at least 10 minutes for the console to warm up before operation of the MPI.
3. Execute the KSR program. If the program is not initiated, the rotary gas bearing system will not be supplied air by the compressor.
4. Press any key to display the main menu (Appendix C).
5. Select the update test information menu and enter the appropriate information. It is important to update the test part ID before each test performed.
6. Set the gas bearing pressure for the test part. To do so:
 - a. With the air off, mount the test fixture and the test part to the test platform.
 - b. Turn the air on and turn the pressure regulator counterclockwise until the bearing pressure gauge reads zero.
 - c. Slowly turn the regulator clockwise until the test platform and the test part float freely on the gas bearing.
 - d. Note the bearing pressure level.
 - e. Turn the regulator clockwise to get an additional 5 psi on the bearing pressure gauge. **WARNING:** Do not exceed 95 psi on the bearing pressure gauge.
 - f. Turn the air off at the toggle switch to lock the test platform in place and remove the test part from the fixture in order to prepare for the test sequence.

Calibration

Individual calibrations must be performed for the CM and the MOI. Both require the use of a calibration beam and weights (Figure 8). The calibration beam must be installed with holes A, B, and C in the quadrant between the +X and +Y axes. The computer will prompt the operator to place specific weights in particular holes on the calibration beam (Figure 9). The MPI system calibration essentially is a comparison of the measured CM

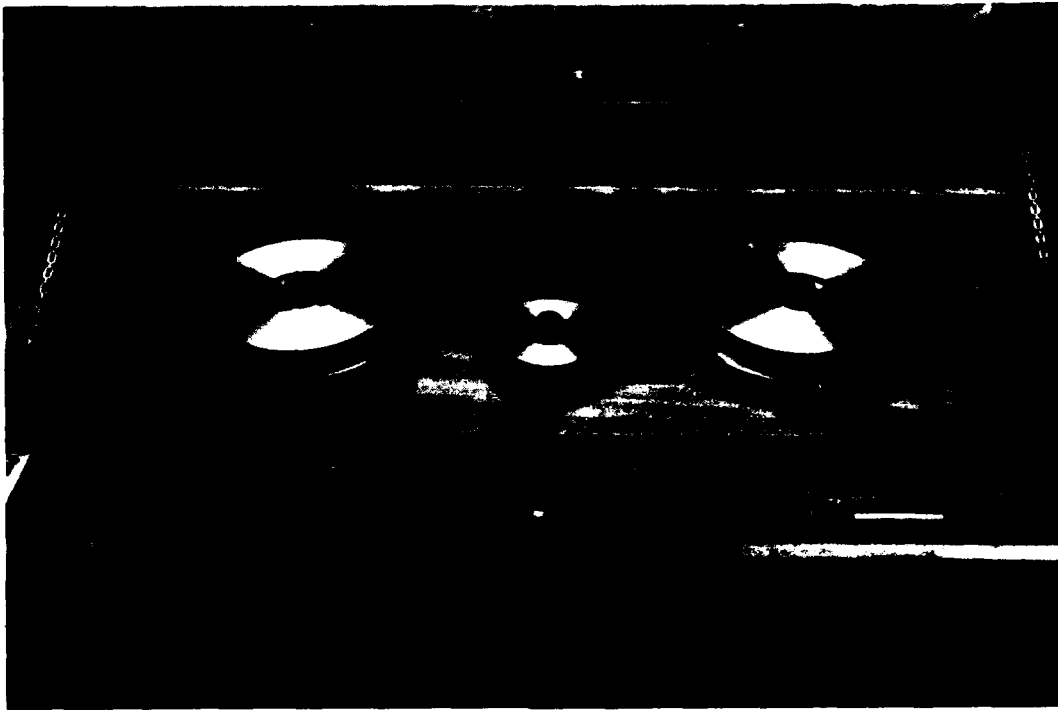


Figure 8. Calibration beam and weights.



Figure 9. MPI setup for CM calibration.

and MOI values of the simple weight against their theoretical CM and MOI values. After the sequence is complete, the computer will display that the calibration was accepted or rejected. An accepted calibration indicates the machine is operating correctly and will provide reliable data. If rejected, the system is faulty and a diagnosis must be performed to find the source of error. Possible sources include: floor vibration, air drafts, touching the MPI during testing, improper system setup, etc.

Measurements

In order to measure the CM and MOI of a test part, the weight of the part and the estimated height of the part's CM above the test platform must be entered into the computer. The weight of the test part is crucial in obtaining accurate results and should be checked before each calculation is made. The estimated CM height has only a minor affect on the actual calculations of the CM. The measurements consist of three steps:

1. Tare measurement.
2. Part measurement.
3. Calculations.

Step 1: Tare measurement

The tare measurement is performed to eliminate the test fixture from the test part's CM and MOI calculations. During this measurement, the fixture to be used to mount the test part onto the platform should be attached. Test fixtures should be fabricated to mount to the test platform (Appendix D). It is important to mount the test fixture so that the defined test part coordinate system aligns with the test platform system. If an offset is present, it should be noted and be accounted for in the CM and MOI test results.

As with the calibration, individual tare measurements must be performed for the CM and the MOI. With the fixture in place, press the moment display button off and then on to zero the display. The operator then should follow the screen prompts to initiate the tare measurement. The computer will take readings from each quadrant and print the results. All subsequent part measurements will compensate for the test fixture's properties. If the fixture is to be placed in a different configuration, another tare measurement must be performed for the new fixture orientation.

Step 2: Part measurement

The test part now should be installed in the fixture and part measurement selected from the main menu. Again, the CM and MOI tests are two separate functions. For each test, the computer will prompt the operator through the sequence and print the raw data.

Step 3: Calculations

To obtain the test results, the calculations menu should be selected. The program provides the option to calculate the CM and MOI separately or together. After selecting which calculations are to be made, be sure to verify that the actual weight of the test part and the estimated CM height are entered properly. These values are crucial to obtaining correct data. The computer will perform the calculations and print the data. The MOI is given about the test part's CM and the platform's coordinate system.

Conclusions

After following the procedure to measure the mass properties of a test part, several conclusions were made:

1. The positioning of the fixture is crucial to obtain accurate results. It is necessary for the test part to be referenced to the pivot axis of the test platform. If the defined axis of the test part is not aligned with the pivot axis of the test platform, an offset distance will have to be figured into the results from the CM calculations.

2. Vertical components of the test fixture must be perpendicular to the test platform. The MPI is incapable of factoring the lean of the vertical components into the calculations. Therefore, the results will be erroneous if the fixture is not perpendicular.

3. It is important to update the test information before each test. By changing the test part ID each time, the results will be labeled appropriately on the printouts. This is very important when performing a series of tests.

4. Each time calculations are done, the mass of the test part should be entered. This will minimize the possibility of having the calculations based on the wrong test part weight.

5. The MPI must be absolutely isolated from any vibration in the surrounding environment. Vibration isolators were added to the described MPI system.

6. A test fixture may require the part to be attached in such a way that would cause an offset between the center axis of the test part and that of the test platform. If this is the case, be sure to attach the part the same way each time so that the offset is always in the same direction.

7. The maximum moment allowed due to the offset of the CM of a test part is 60 lb/in. When designing test fixtures, this restriction must be kept in mind.

8. Accurate results are obtained by running both the CM and MOI tests together and then performing the calculations for both. The MPI calculates the MOI results about the test part's CM. If the CM was never determined, the MOI results will be based on the CM of the last part tested and will be wrong.

References

- Beer, Ferdinand P., and Johnston, Jr., E. Russell. 1977. Vector mechanics for engineers: Statics and dynamics. 3rd ed. New York: McGraw-Hill Book Company.
- Meriam, J. L., and Kraige, L. G. 1986. Engineering mechanics: Dynamics. Vol 2. 2nd ed. New York: John Wiley and Sons.
- Space Electronics, Inc. n.d. Mass Properties Instrument Model KSR330-60 instruction manual. Berlin, CT: Space Electronics, Inc.

Appendix A.

Accuracy of the MPI

The accuracies of the MPI are as follows:

CM error for 100-lb part with 0.5-in. CM offset

Bearing/machine tolerance.	0.001000 in.
Sensitivity error (0.003/part weight in lbs.	0.000030 in.
Linearity error (0.03% of full scale).	0.000180 in.
Max uncertainty this example.	<u>0.001210 in.</u>

Moment of inertia

Basic accuracy	$\pm 0.25\%$
Tare MOI (typical)	130 lb-in ²
MOI error for 800 lb-in ² part:	
Basic accuracy	2.0 lb-in ²

Appendix B.

Determination of torsional pendulum physical constant

The calibration constant primarily is a function of the spring rate of the torsion rod. A calibration weight is placed at a known distance from the center of the test platform and the period of oscillation (T_c) is measured. The weight then is placed at the center of the platform and a second period (T_o) is found. The difference between the two readings is completely due to the change in MOI. The MOI change, I , is defined as $I = WR^2$, where W is the calibration weight and R is the offset. The calibration constant then is determined:

$$C = I / (T_c^2 - T_o^2).$$

Appendix C.

KSR main menu

The main menu of the KSR program appears as follows:

>>>> MAIN MENU <<<<<

F1: Update test information
F2: Part measurement
F3: Tare measurement
F4: Calibration
F5: System utilities
F6: Calculations
<ESC>: Quit

Select function from list

Appendix D.

Mounting of test fixture

Any part to be tested must be interfaced with the test platform by means of a test fixture. The fixture must be fabricated to fit the hole pattern of the test platform. This pattern is detailed in Figure D-1 shown below.

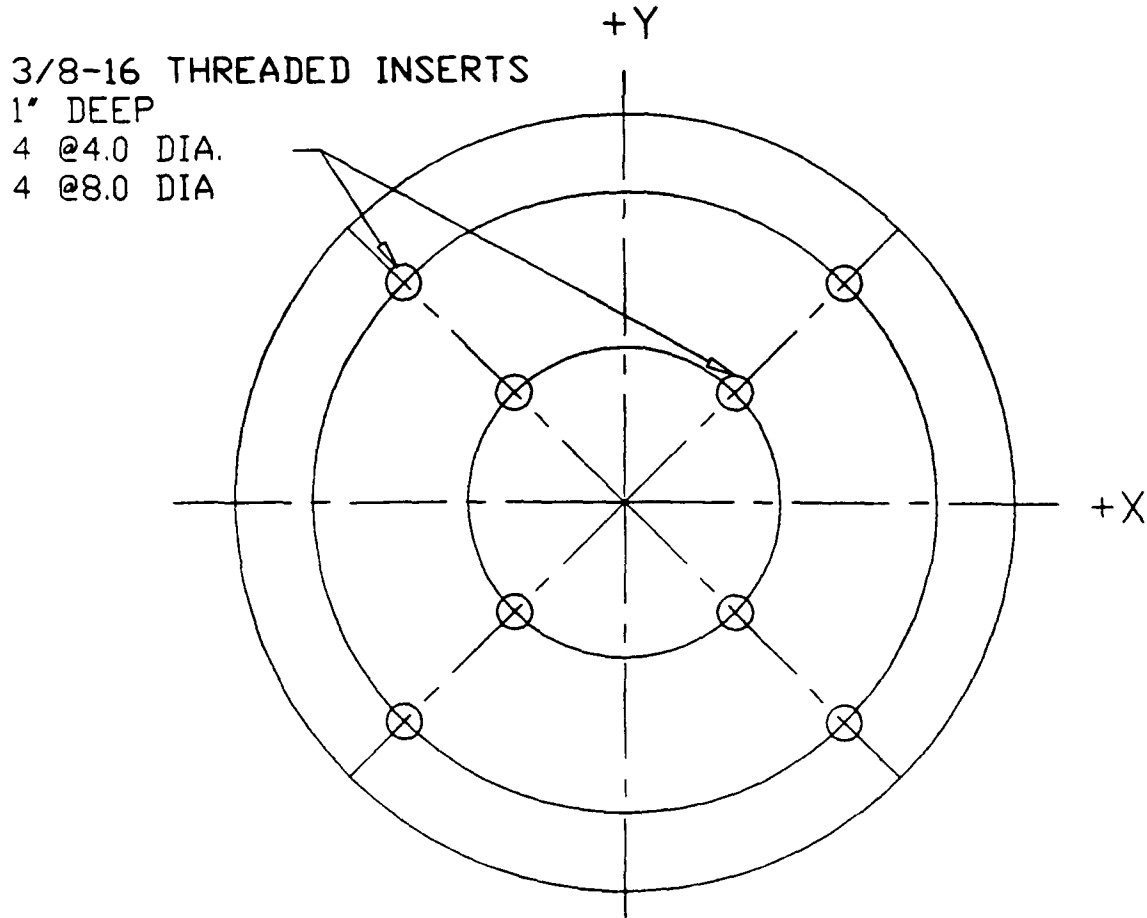


Figure D-1. Mounting pattern of test platform.

Appendix E.

Manufacturers' list

Space Electronics, Inc.
81 Fuller Way
Berlin, CT 06037

JUN-AIR (USA) Inc.
1303 Barclay Blvd.
Buffalo Grove, IL 60089

Initial distribution

Commander, U.S. Army Natick Research,
Development and Evaluation Center
ATTN: STRNC-MIL (Documents
Librarian)
Natick, MA 01760-5040

U.S. Army Communications-Electronics
Command
ATTN: AMSEL-RD-ESA-D
Fort Monmouth, NJ 07703

Commander/Director
U.S. Army Combat Surveillance
and Target Acquisition Lab
ATTN: DELCS-D
Fort Monmouth, NJ 07703-5304

Commander
10th Medical Laboratory
ATTN: Audiologist
APO New York 09180

Naval Air Development Center
Technical Information Division
Technical Support Detachment
Warminster, PA 18974

Commanding Officer, Naval Medical
Research and Development Command
National Naval Medical Center
Bethesda, MD 20814-5044

Deputy Director, Defense Research
and Engineering
ATTN: Military Assistant
for Medical and Life Sciences
Washington, DC 20301-3080

Commander, U.S. Army Research
Institute of Environmental Medicine
Natick, MA 01760

U.S. Army Avionics Research
and Development Activity
ATTN: SAVAA-P-TP
Fort Monmouth, NJ 07703-5401

Library
Naval Submarine Medical Research Lab
Box 900, Naval Sub Base
Groton, CT 06349-5900

Director, U.S. Army Human
Engineering Laboratory
ATTN: Technical Library
Aberdeen Proving Ground, MD 21005

Commander
Man-Machine Integration System
Code 602
Naval Air Development Center
Warminster, PA 18974

Commander
Naval Air Development Center
ATTN: Code 602-B (Mr. Brindle)
Warminster, PA 18974

Commanding Officer
Armstrong Laboratory
Wright-Patterson
Air Force Base, OH 45433-6573

Director
Army Audiology and Speech Center
Walter Reed Army Medical Center
Washington, DC 20307-5001

Commander, U.S. Army Institute
of Dental Research
ATTN: Jean A. Setterstrom, Ph. D.
Walter Reed Army Medical Center
Washington, DC 20307-5300

Naval Air Systems Command
Technical Air Library 950D
Room 278, Jefferson Plaza II
Department of the Navy
Washington, DC 20361

Commander, U.S. Army Test
and Evaluation Command
ATTN: AMSTE-AD-H
Aberdeen Proving Ground, MD 21005

Director
U.S. Army Ballistic
Research Laboratory
ATTN: DRXBR-OD-ST Tech Reports
Aberdeen Proving Ground, MD 21005

Commander
U.S. Army Medical Research
Institute of Chemical Defense
ATTN: SGRD-UV-AO
Aberdeen Proving Ground,
MD 21010-5425

Commander, U.S. Army Medical
Research and Development Command
ATTN: SGRD-RMS (Ms. Madigan)
Fort Detrick, Frederick, MD 21702-5012

Director
Walter Reed Army Institute of Research
Washington, DC 20307-5100

HQ DA (DASG-PSP-O)
5109 Leesburg Pike
Falls Church, VA 22041-3258

Harry Diamond Laboratories
ATTN: Technical Information Branch
2800 Powder Mill Road
Adelphi, MD 20783-1197

U.S. Army Materiel Systems
Analysis Agency
ATTN: AMXSY-PA (Reports Processing)
Aberdeen Proving Ground
MD 21005-5071

U.S. Army Ordnance Center
and School Library
Simpson Hall, Building 3071
Aberdeen Proving Ground, MD 21005

U.S. Army Environmental
Hygiene Agency
Building E2100
Aberdeen Proving Ground, MD 21010

Technical Library Chemical Research
and Development Center
Aberdeen Proving Ground, MD
21010--5423

Commander
U.S. Army Medical Research
Institute of Infectious Disease
SGRD-UIZ-C
Fort Detrick, Frederick, MD 21702

Director, Biological
Sciences Division
Office of Naval Research
600 North Quincy Street
Arlington, VA 22217

Commander
U.S. Army Materiel Command
ATTN: AMCDE-XS
5001 Eisenhower Avenue
Alexandria, VA 22333

Commandant
U.S. Army Aviation
Logistics School ATTN: ATSQ-TDN
Fort Eustis, VA 23604

Headquarters (ATMD)
U.S. Army Training
and Doctrine Command
ATTN: ATBO-M
Fort Monroe, VA 23651

Structures Laboratory Library
USARTL-AVSCOM
NASA Langley Research Center
Mail Stop 266
Hampton, VA 23665

Naval Aerospace Medical
Institute Library
Building 1953, Code 03L
Pensacola, FL 32508-5600

Command Surgeon
HQ USCENTCOM (CCSG)
U.S. Central Command
MacDill Air Force Base FL 33608

Air University Library
(AUL/LSE)
Maxwell Air Force Base, AL 36112

U.S. Air Force Institute
of Technology (AFIT/LDEE)
Building 640, Area B
Wright-Patterson
Air Force Base, OH 45433

Henry L. Taylor
Director, Institute of Aviation
University of Illinois-Willard Airport
Savoy, IL 61874

Chief, Nation Guard Bureau
ATTN: NGB-ARS (COL Urbauer)
Room 410, Park Center 4
4501 Ford Avenue
Alexandria, VA 22302-1451

Commander
U.S. Army Aviation Systems Command
ATTN: SGRD-UAX-AL (MAJ Gillette)
4300 Goodfellow Blvd., Building 105
St. Louis, MO 63120

U.S. Army Aviation Systems Command
Library and Information Center Branch
ATTN: AMSAV-DIL
4300 Goodfellow Boulevard
St. Louis, MO 63120

Federal Aviation Administration
Civil Aeromedical Institute
Library AAM-400A
P.O. Box 25082
Oklahoma City, OK 73125

Commander
U.S. Army Academy
of Health Sciences
ATTN: Library
Fort Sam Houston, TX 78234

Commander
U.S. Army Institute of Surgical Research
ATTN: SGRD-USM (Jan Duke)
Fort Sam Houston, TX 78234-6200

AAMRL/HEX
Wright-Patterson
Air Force Base, OH 45433

John A. Dellinger,
Southwest Research Institute
P. O. Box 28510
San Antonio, TX 78284

Product Manager
Aviation Life Support Equipment
ATTN: AMCPM-ALSE
4300 Goodfellow Boulevard
St. Louis, MO 63120-1798

Commander
U.S. Army Aviation
Systems Command
ATTN: AMSAV-ED
4300 Goodfellow Boulevard
St. Louis, MO 63120

Commanding Officer
Naval Biodynamics Laboratory
P.O. Box 24907
New Orleans, LA 70189-0407

Assistant Commandant
U.S. Army Field Artillery School
ATTN: Morris Swott Technical Library
Fort Sill, OK 73503-0312

Commander
U.S. Army Health Services Command
ATTN: HSOP-SO
Fort Sam Houston, TX 78234-6000

HQ USAF/SGPT
Bolling Air Force Base, DC 20332-6188

U.S. Army Dugway Proving Ground
Technical Library, Building 5330
Dugway, UT 84022

U.S. Army Yuma Proving Ground
Technical Library
Yuma, AZ 85364

AFFTC Technical Library
6510 TW/TSTL
Edwards Air Force Base,
CA 93523-5000

Commander
Code 3431
Naval Weapons Center
China Lake, CA 93555

Aeromechanics Laboratory
U.S. Army Research and Technical Labs
Ames Research Center, M/S 215-1
Moffett Field, CA 94035

Sixth U.S. Army
ATTN: SMA
Presidio of San Francisco, CA 94129

Commander
U.S. Army Aeromedical Center
Fort Rucker, AL 36362

U.S. Air Force School
of Aerospace Medicine
Strughold Aeromedical Library Technical
Reports Section (TSKD)
Brooks Air Force Base, TX 78235-5301

Dr. Diane Damos
Department of Human Factors
ISSM, USC
Los Angeles, CA 90089-0021

U.S. Army White Sands
Missile Range
ATTN: STEWS-IM-ST
White Sands Missile Range, NM 88002

U.S. Army Aviation Engineering
Flight Activity
ATTN: SAVTE-M (Tech Lib) Stop 217
Edwards Air Force Base, CA 93523-5000

Ms. Sandra G. Hart
Ames Research Center
MS 262-3
Moffett Field, CA 94035

Commander, Letterman Army Institute
of Research
ATTN: Medical Research Library
Presidio of San Francisco, CA 94129

Commander
U.S. Army Medical Materiel
Development Activity
Fort Detrick, Frederick, MD 21702-5009

Commander
U.S. Army Aviation Center
Directorate of Combat Developments
Building 507
Fort Rucker, AL 36362

U. S. Army Research Institute
Aviation R&D Activity
ATTN: PERI-IR
Fort Rucker, AL 36362

Commander
U.S. Army Safety Center
Fort Rucker, AL 36362

U.S. Army Aircraft Development
Test Activity
ATTN: STEBG-MP-P
Cairns Army Air Field
Fort Rucker, AL 36362

Commander U.S. Army Medical Research
and Development Command
ATTN: SGRD-PLC (COL Schnakenberg)
Fort Detrick, Frederick, MD 21702

MAJ John Wilson
TRADOC Aviation LO
Embassy of the United States
APO New York 09777

Netherlands Army Liaison Office
Building 602
Fort Rucker, AL 36362

British Army Liaison Office
Building 602
Fort Rucker, AL 36362

Italian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Directorate of Training Development
Building 502
Fort Rucker, AL 36362

Chief
USAHEL/USAAVNC Field Office
P. O. Box 716
Fort Rucker, AL 36362-5349

Commander U.S. Army Aviation Center
and Fort Rucker
ATTN: ATZQ-CG
Fort Rucker, AL 36362

Chief
Test & Evaluation Coordinating Board
Cairns Army Air Field
Fort Rucker, AL 36362

MAJ Terry Newman
Canadian Army Liaison Office
Building 602
Fort Rucker, AL 36362

German Army Liaison Office
Building 602
Fort Rucker, AL 36362

LTC Patrice Cottebrune
French Army Liaison Office
USAAVNC (Building 602)
Fort Rucker, AL 36362-5021

Australian Army Liaison Office
Building 602
Fort Rucker, AL 36362

Dr. Garrison Rapmund
6 Burning Tree Court
Bethesda, MD 20817

Commandant, Royal Air Force
Institute of Aviation Medicine
Farnborough Hampshire GU14 6SZ UK

Commander
U.S. Army Biomedical Research
and Development Laboratory
ATTN: SGRD-UBZ-I
Fort Detrick, Frederick, MD 21702

Defense Technical Information OCC
Selection
Cameron Station
Alexandra, VA 22313

Commander, U.S. Army Foreign Science
and Technology Center
AIFRTA (Davis)
220 7th Street, NE
Charlottesville, VA 22901-5396

Director,
Applied Technology Laboratory
USARTL-AVSCOM
ATTN: Library, Building 401
Fort Eustis, VA 23604

U.S. Air Force Armament
Development and Test Center
Eglin Air Force Base, FL 32542

Commander, U.S. Army Missile
Command
Redstone Scientific Information Center
ATTN: AMSMI-RD-CS-R
/ILL Documents
Redstone Arsenal, AL 35898

Dr. H. Dix Christensen
Bio-Medical Science Building, Room 753
Post Office Box 26901
Oklahoma City, OK 73190

Director
Army Personnel Research Establishment
Farnborough, Hants GU14 6SZ UK

U.S. Army Research and Technology
Laboratories (AVSCOM)
Propulsion Laboratory MS 302-2
NASA Lewis Research Center
Cleveland, OH 44135

Dr. Christine Schlichting
Behavioral Sciences Department
Box 900, NAVUBASE NLON
Groton, CT 06349-5900

Dr. Eugene S. Channing
7985 Schooner Court
Frederick, MD 21701-3273

LTC Gaylord Lindsey (5)
USAMRDC Liaison at Academy
of Health Sciences
ATTN: HSHA-ZAC-F
Fort Sam Houston, TX 78234

Aviation Medicine Clinic
TMC #22, SAAF
Fort Bragg, NC 28305

Dr. A. Kornfield, President
Biosearch Company
3016 Revere Road
Drexel Hill, PA 29026

NVEOD
AMSEL-RD-ASID
(Attn: Trang Bui)
Fort Belvoir, VA 22060

CA Av Med
HQ DAAC
Middle Wallop
Stockbridge Hants S020 8DY UK

**Commander and Director
USAE Waterways Experiment Station
ATTN: CEWES-IM-MI-R
Alfrieda S. Clark, CD Dept.
3909 Halls Ferry Road
Vicksburg, MS 39180-6199**

**Col. Otto Schramm Filho
c/o Brazilian Army Commission
Office-CEBW
4632 Wisconsin Avenue NW
Washington, DC 20016**

**Mr. Richard Thornley
ILS Manager, APACHE
Box 84
Westland Helicopters Limited
Yeovil, Somerset BA202YB UK**