

# DNA Cleavage Using Artificial Hydrolases Based on Cobalt (III) Cyclen Complexes


Eddie Chang\*, Rajeev Soni\*, Ellen R. Goldman\*, James B. Delehanty\*, Jason Bongard#, Frederick Strich# and, D. Andrew Knight#

\*Center for Bio/Molecular Science and Engineering, Naval Research Laboratory, Washington, D.C. 20375

#Department of Chemistry, Loyola University, 6363 St. Charles Avenue, New Orleans, LA 70115



**BAMHI**

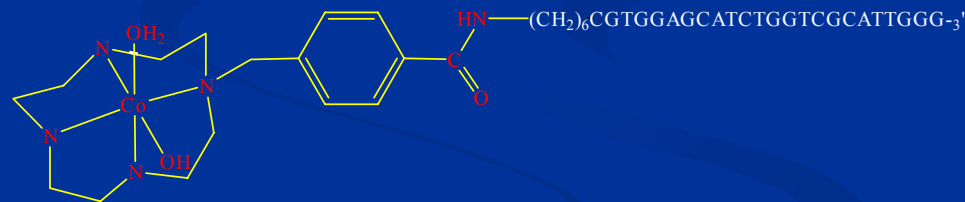

**Artificial Hydrolase**



# Report Documentation Page

*Form Approved*  
*OMB No. 0704-0188*

Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.

1. REPORT DATE <b>01 OCT 2005</b>	2. REPORT TYPE <b>N/A</b>	3. DATES COVERED <b>-</b>	
4. TITLE AND SUBTITLE <b>DNA Cleavage Using Artificial Hydrolases Based on Cobalt (III) Cyclen Complexes</b>		5a. CONTRACT NUMBER	
		5b. GRANT NUMBER	
		5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)		5d. PROJECT NUMBER	
		5e. TASK NUMBER	
		5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) <b>Center for Bio/Molecular Science and Engineering, Naval Research Laboratory, Washington, D.C. 20375 Washington, D.C. 20375</b>		8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)		10. SPONSOR/MONITOR'S ACRONYM(S)	
		11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT <b>Approved for public release, distribution unlimited</b>			
13. SUPPLEMENTARY NOTES <b>See also ADM001851, Proceedings of the 2003 Joint Service Scientific Conference on Chemical &amp; Biological Defense Research, 17-20 November 2003. , The original document contains color images.</b>			
14. ABSTRACT			
15. SUBJECT TERMS			
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT
a. REPORT <b>unclassified</b>	b. ABSTRACT <b>unclassified</b>	c. THIS PAGE <b>unclassified</b>	<b>UU</b>
			18. NUMBER OF PAGES <b>20</b>
			19a. NAME OF RESPONSIBLE PERSON

# ABSTRACT

We are developing small, robust catalytic systems to efficiently degrade DNA/RNA with a high degree of specificity. Here, we describe the use of tetraazamacrocyclic-chelated Co(III) complexes for the effective hydrolysis of double stranded plasmid DNA. These chelator molecules have been functionalized with a coupling group and we examine whether the modification has affected the ability of the complex to degrade DNA. Hydrolysis of DNA does appear to be slowed by the intermediate functionalization, but the eventual form of the coupled chelator-complex shows very high activity against a model pBluescript supercoiled DNA. We discuss these results, future tests, and plans for using these systems for sensitive equipment decontamination.




# Introduction

Naturally occurring nucleases, while potent as native defenses against foreign DNA, are limited by their short sequence-recognition capability and temperature instability. Development of "artificial nucleases" can overcome these limitations.

We seek to develop strongly binding artificial nucleases, with enhanced sequence recognition/ specificity and stability and to determine their suitability for anti-bacterial and anti-viral applications. These systems will be developed specifically to disrupt or destroy viruses and/or bacteria either through direct contact of the nucleases with the target genetic material or by delivering the nucleases through the protective envelopes of microorganisms using suitable carriers. The novel nucleases will be tested as decontamination agents.

# Components of Our System


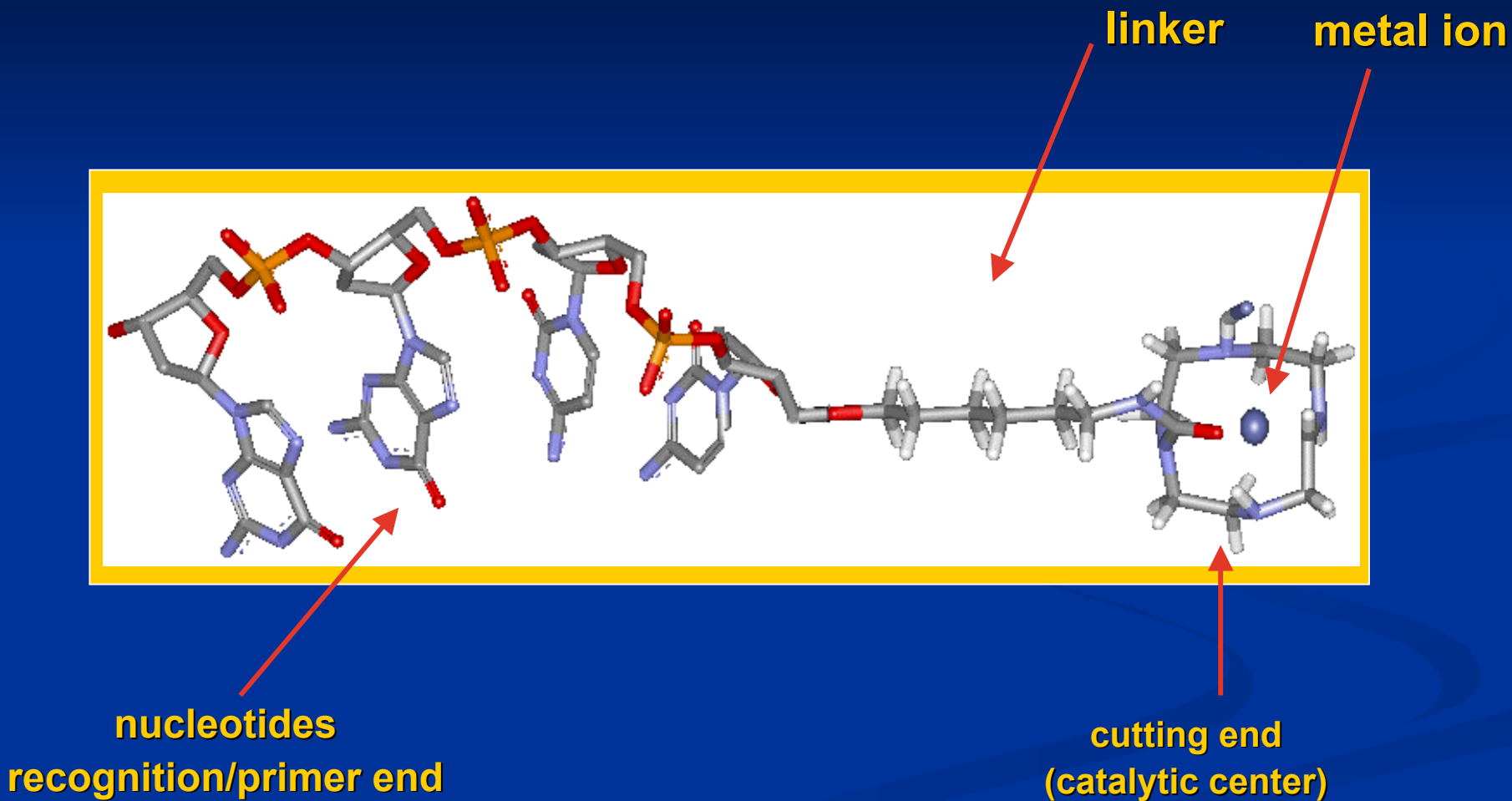


# Co(III)-Chelator Characteristics

- Among the fastest in hydrolysis of phosphodiester systems
- Kinetically inert, i.e., no turnover at RT
- "Promotes" hydrolysis, rather than "catalytic" at RT
- Inertness to help suppress viral growth



# How Fast?

	Cyclen	TRPN
BNPP	$4.6 \times 10^{-3} \text{ s}^{-1}$	$2.5 \times 10^{-2} \text{ s}^{-1}$
DNA	$2 \times 10^{-4} \text{ s}^{-1}$	$5 \times 10^{-5} \text{ s}^{-1}$

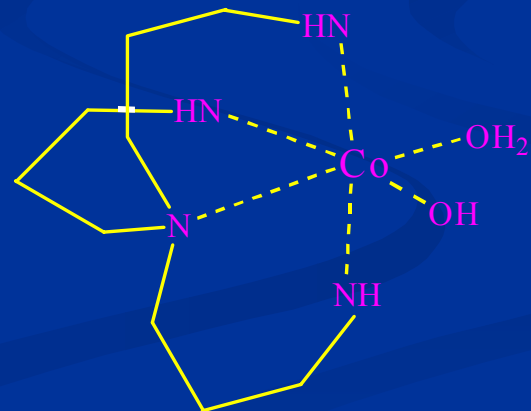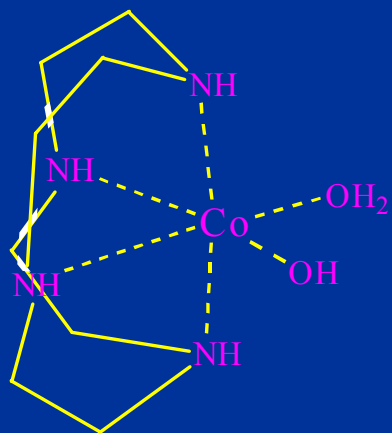


# Kinetic Inertness

Substitutional Inertness:

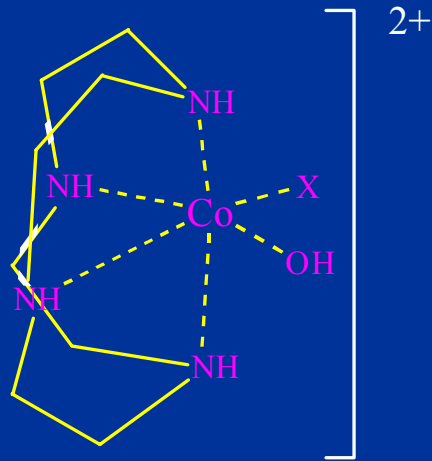
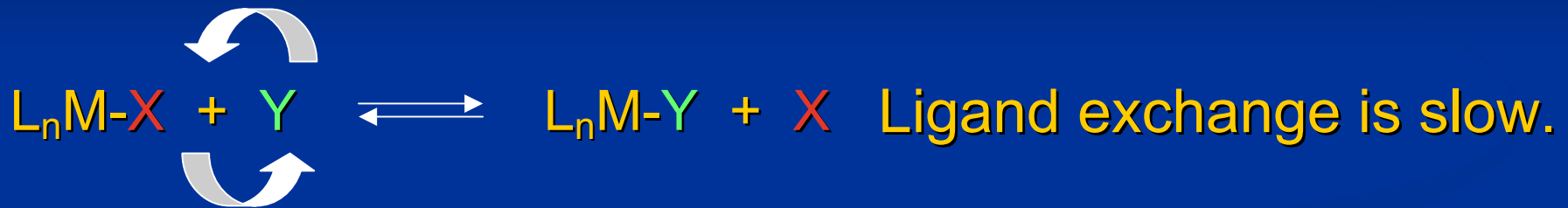


Example of a kinetically inert  
chelator-Co(III) complex


# Anti-Sense Paradigm

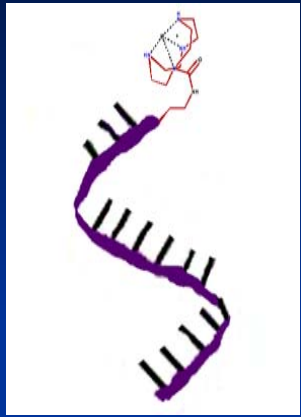

+

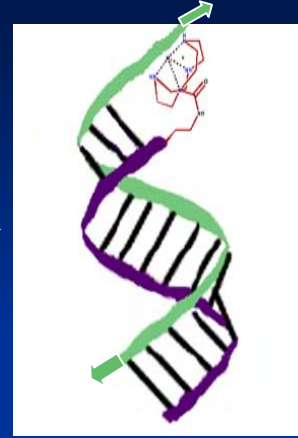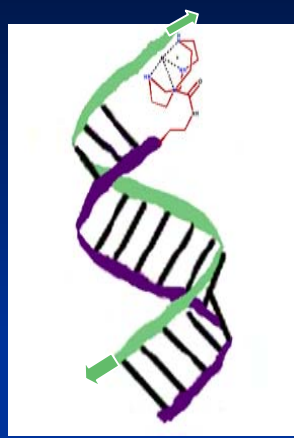

+

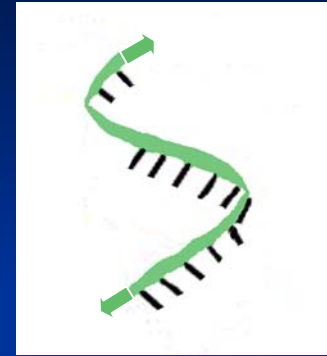

ODN probe + catalyst

RNA

probe + target RNA  
+ catalyst

probe + target RNA  
fragment + catalyst

RNA fragment




Hydrolysis




probe + target RNA + catalyst

kinetic inertness  
used for added  
steric blockade  
advantage

# CoCyclen and Its Functionalized Forms

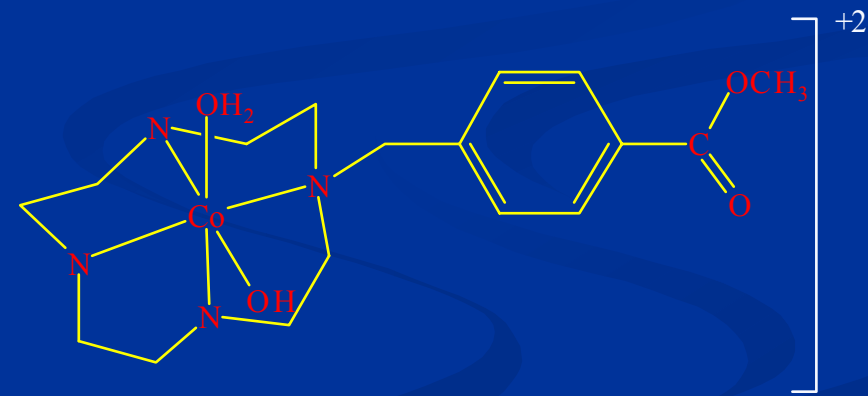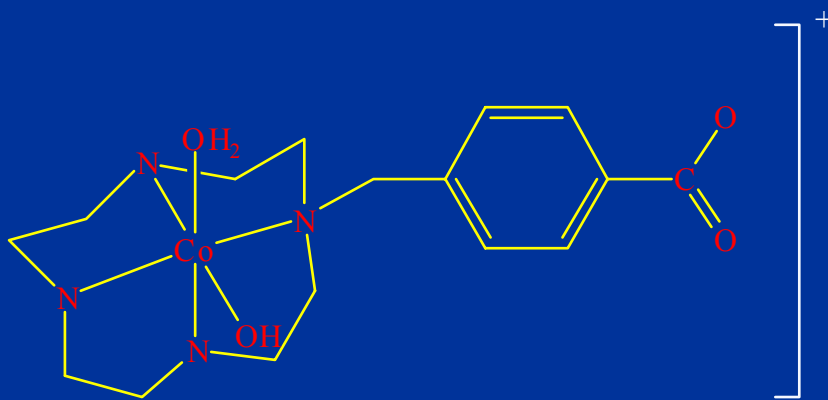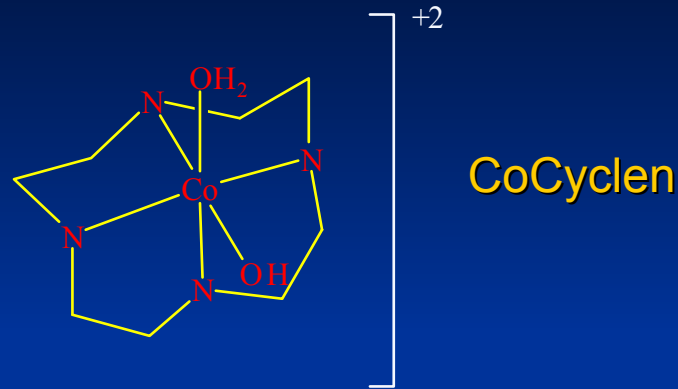

# CoCyclen Hydrolysis of dsDNA

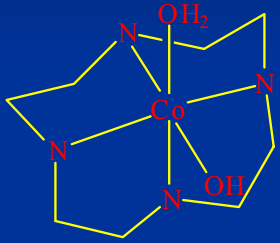

5mM Cocyc + DNA 37°C

2.5mM Cocyc + DNA 37°C

1.25mM Cocyc + DNA 37°C

DNA 37°C

5mM Cocyc + DNA 37°C

2.5mM Cocyc + DNA 37°C

1.25mM Cocyc + DNA 37°C

DNA 37°C


H  
E  
A  
T  
80°C, 5 min  
25°C, 5min

uncut dsDNA ↑

↑ nicked DNA

Unmodified chelator hydrolysis of dsDNA at 37°C for 2 h.  
Increasing nicked DNA as a function of Co(III) concentration.

# CoCyclenmba Hydrolysis of dsDNA

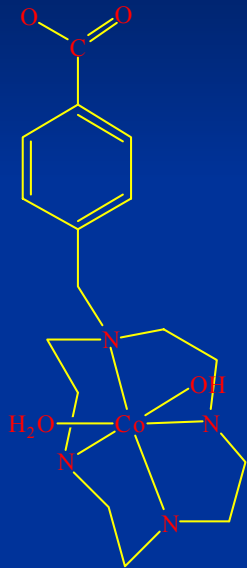

pH 10 control

pH 10

pH 8 control

pH 8

pH 7 control

pH 7

Nickase control

Nickase

pH 7 fastest for functionalized hydrolysis center. Possible slight increase at pH 10.

# Cocyclenmba Results

Lane	Supercoiled	Linear	Nicked	Comments
8	93%	3%	4%	pH 10 cntl
7	84%	5%	12%	pH 10
6	92%	5%	3%	pH 8 cntl
5	87%	5%	9%	ph8
4	90%	4%	5%	pH 7 cntl
3	72%	3%	25%	pH 7
2	93%	7%	0%	nickase cntl
1	10%	4%	87%	nickase



Cocyclenmba is most active at pH 7 with about 25% nicked population.

# CoCyclenmba Concentration-Dependence



Lane	[cocycmba]	[DNA]
1.	1.7 mM (ctrl)	0.139 mM bp
2.	1.7 mM	0.139 mM bp
3.	3.4 mM (ctrl)	0.139 mM bp
4.	3.4 mM	0.139 mM bp
5.	5 mM	0.139 mM bp
6.	5 mM (ctrl)	0.139 mM bp
7.	5 mM	0.0696 mM bp
8.	5 mM (ctrl)	0.0696 mM bp
9.	5 mM	0.0348 mM bp
10.	5 mM (ctrl)	0.0348mM bp

# Cocyclenmmb Hydrolysis of dsDNA

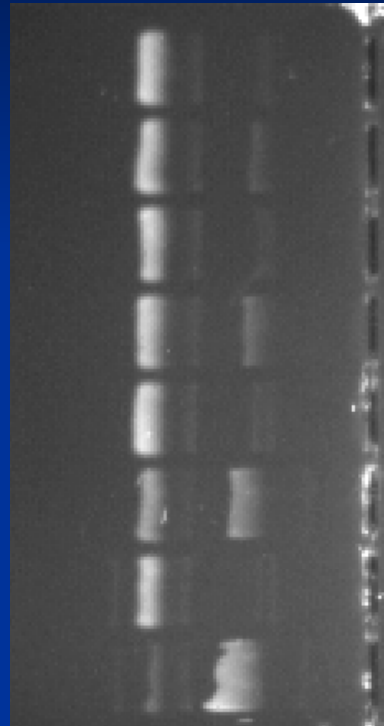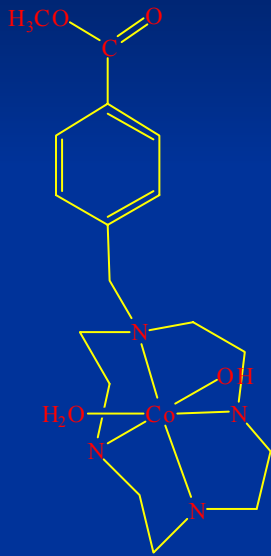

pH 10 control

pH 10

pH 8 control

pH 8

pH 7 control

pH 7

Nickase control

Nickase

pH 7 fastest for "coupled" form of the hydrolysis center--about 50% of nickase activity.

# Cocyclenmme Results

Lane	Supercoiled	Linear	Nicked	Comments
8	93%	3%	4%	pH 10 cntl
7	89%	4%	7%	pH 10
6	90%	5%	5%	pH 8 cntl
5	81%	4%	15%	ph8
4	91%	4%	5%	pH 7 cntl
3	54%	2%	43%	pH 7
2	94%	2%	4%	nickase cntl
1	9%	3%	88%	nickase


Cocyclenmme is most active at pH 7 with 43% nicked population-- almost half of the nickase rate.

# Next Steps

## Sequence-Specific Cut of a ssDNA

5' GCGCTGACCC GCACCTCGTA GACCAGCGTA  
ACCCAGTGGT CGTTTAGCGC GACAATCGCC  
CGGGTAATTC AAGACAGAGC CGCGCAGACG  
CAGACCGACC 3'

100-mer target  
ssDNA (substrate)

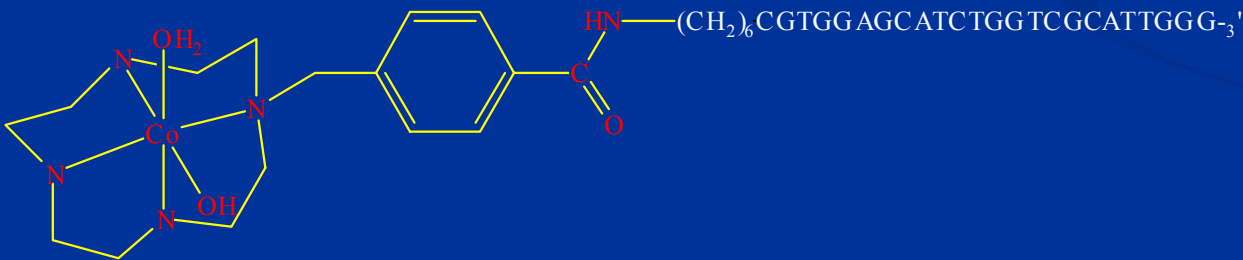


24-mer CoCyclenmba  
complex

2 products around 34- and 66-  
bases long




# Cell-Free Demonstration

## Inhibition of a transcription/translation system

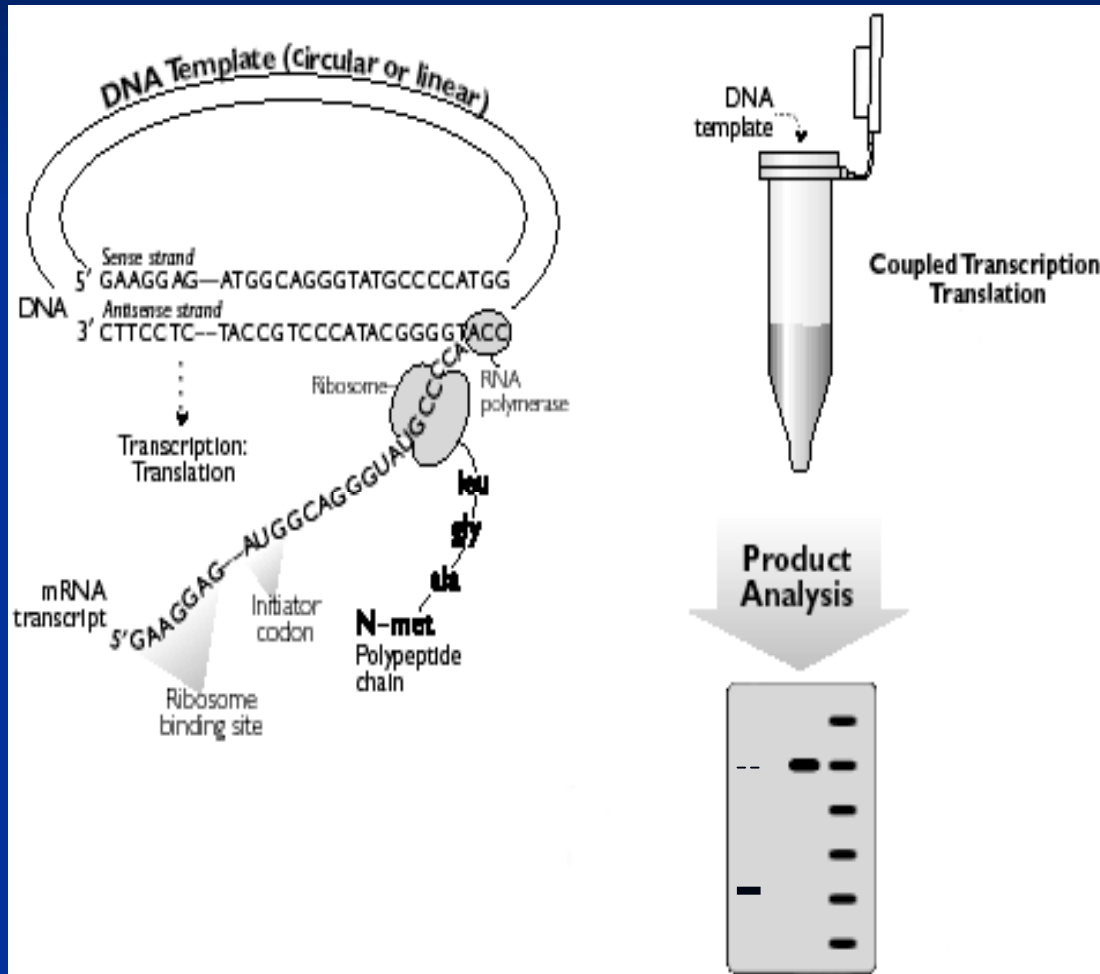


# Decontamination

Package into delivery platform, e.g., fusion-enabled vesicles:

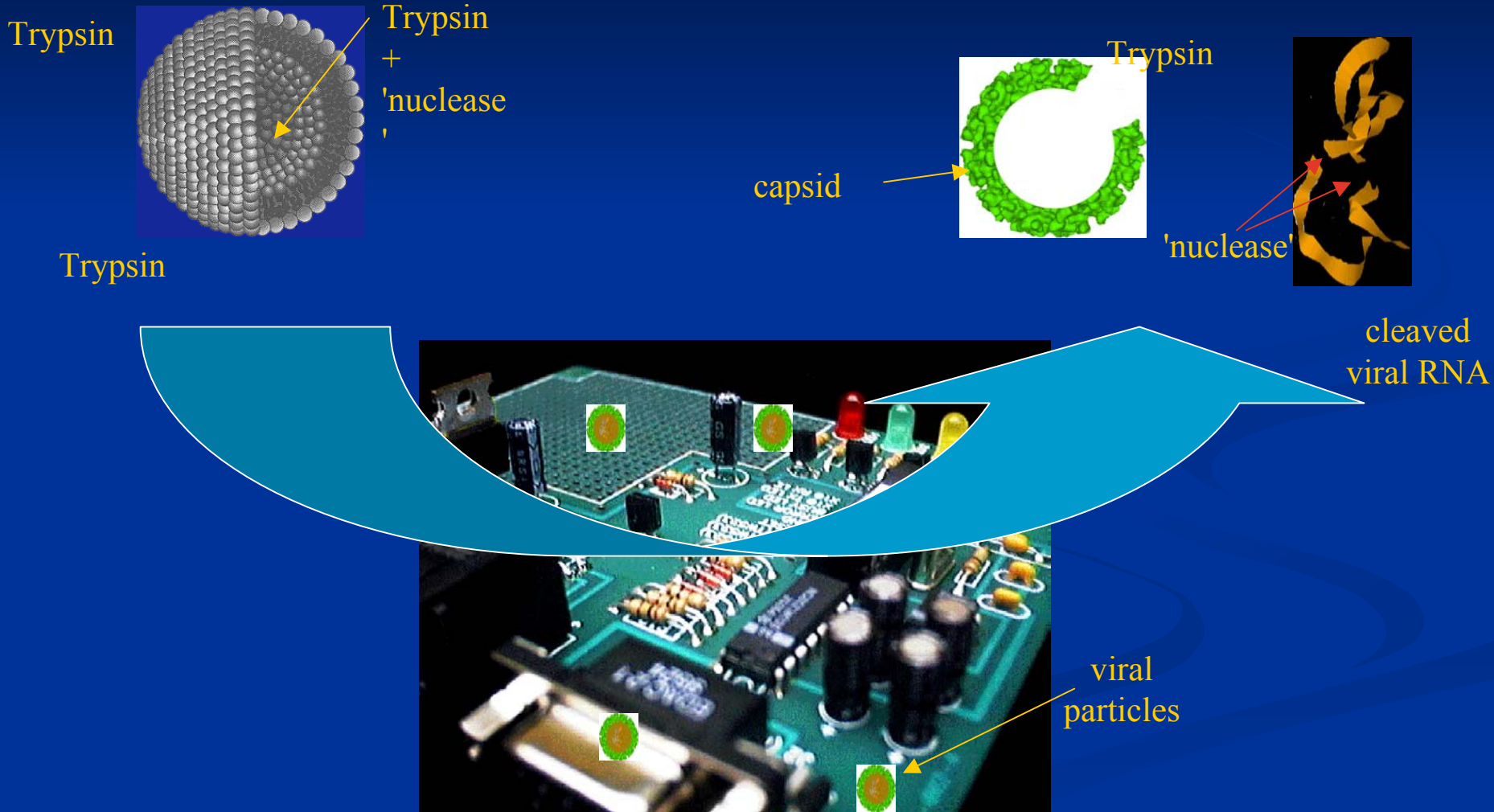

# Coupled to Beads

We have also coupled the functionalized complex to solid-supports, such as agarose beads. These can either serve as filtration elements or as a way to produce ssDNA-chelator complexes.

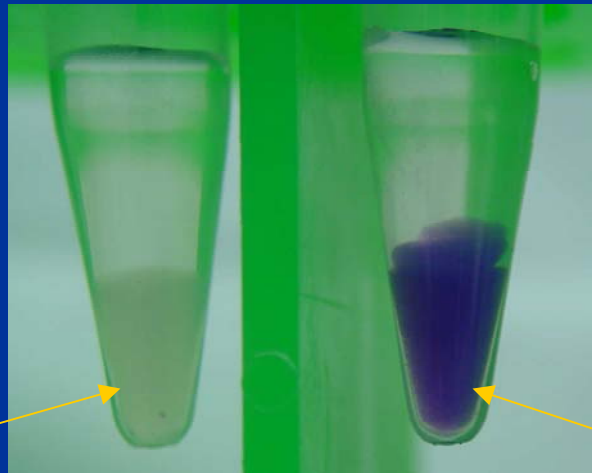

Control agarose  
beads

Cocycmba-  
coupled agarose  
beads

# Conclusion

We have synthesized a Co(III)-chelator complex with a coupling functionality and tested the system against dsDNA, using the amount of nicked dsDNA produced as a measure of activity.

The functionalized complex was active against pBluescript supercoiled DNA with the methyl-ester most active at close to 50% activity of a control nickase enzyme.

Our next steps are to couple the system to a ss-oligo "primer" for sequence-specific cutting and then to a demonstration of cell-free inhibition of protein production via nucleic acid interdiction.

This work was partially supported by ONR and associated DoD organizations.