

Field Portable Concentration and Extraction of Pathogen Analytes from Large Volume Blood Samples



Matt Ewert

mewert@bt.usf.edu

Report Documentation Page				Form Approved OMB No. 0704-0188	
Public reporting burden for the collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to a penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number.					
1. REPORT DATE 18 NOV 2003		2. REPORT TYPE N/A		3. DATES COVERED -	
4. TITLE AND SUBTITLE Field Portable Concentration and Extraction of Pathogen Analytes from Large Volume Blood Samples				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S)				5d. PROJECT NUMBER	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) University of South Florida				8. PERFORMING ORGANIZATION REPORT NUMBER	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release, distribution unlimited					
13. SUPPLEMENTARY NOTES See also ADM001851, Proceedings of the 2003 Joint Service Scientific Conference on Chemical & Biological Defense Research, 17-20 November 2003. , The original document contains color images.					
14. ABSTRACT					
15. SUBJECT TERMS					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT UU	18. NUMBER OF PAGES 34	19a. NAME OF RESPONSIBLE PERSON
a. REPORT unclassified	b. ABSTRACT unclassified	c. THIS PAGE unclassified			

Claim

Levels of bacteremia associated with common bacterial pathogens in adults and children may be less than 1 CFU per ml (3,6). The use of a low volume of blood for a blood culture is known to reduce the detection of bacteremia in adults (5). In an effort to understand how bacterial concentration could contribute clinically significant capabilities to PCR detection systems, blood culture techniques provide relevant information. Tenney et al have shown that 7 ml blood samples resulted in 29 % more positive cultures than found with 2 ml samples (7). A 15 ml culture input sample provided a 25 % greater positive blood culture yield than 5 ml inputs in another study by Hall et al (1). Still others have found improvements for positive blood culture yield when moving from 10 to 20 and 30 ml blood sample inputs (2,4).

1. Hall, M. M., D. M. Ilstrup, J. A. II. Washington. Effect of volume of blood cultured on detection of bacteremia. J. Clin. Microbiol. 3:643-645.
2. Ilstrup, D. M., J. A. II. Washington. The importance of volume of blood cultured in the detection of bacteremia and fungemia. Diagn. Microbiol. Infect. Dis. 1:107-110.
3. Kellogg, J.A., J. P. Manzella, and D. A. Bankert. 2000. Frequency of low level bacteremia in children from birth to fifteen years of age. J. Clin. Microbiol. 38:2181-2185.
4. Li, J., J.J. Plorde, and L. G. Carlson. 1994. Effects of volume and periodicity in blood cultures. J. Clin. Microbiol. 32:2829-2831.
5. Mermel, L. A., and D. G. Maki. Detection of bacteremia in adults: consequences of culturing an inadequate volume of blood. 1993. Ann. Intern. Med. 119:270-272.
6. Reimer, L.G., M. L. Wilson, M. P. Weinstein. 1997. Update on detection of bacteremia and fungemia. Clin. Microbiol. Rev. 10:444-465.
7. Tenney, J. H., L. R. Reller, S. Mirrett, W-L. L. Wang, M. P. Weinstein. 1982. Controlled evaluation of the volume of blood cultured in detection of bacteremia and fungemia. J. Clin. Microbiol. 15:558-561

Collaborations and Coordination

CDC Bioterrorism RRAT Lab

USAMRIID

MIT LL


AFIP

UT

Roche Applied Science, Penzburg, Germany

Path Centre, Perth, Australia

Clinical Testing Scenarios

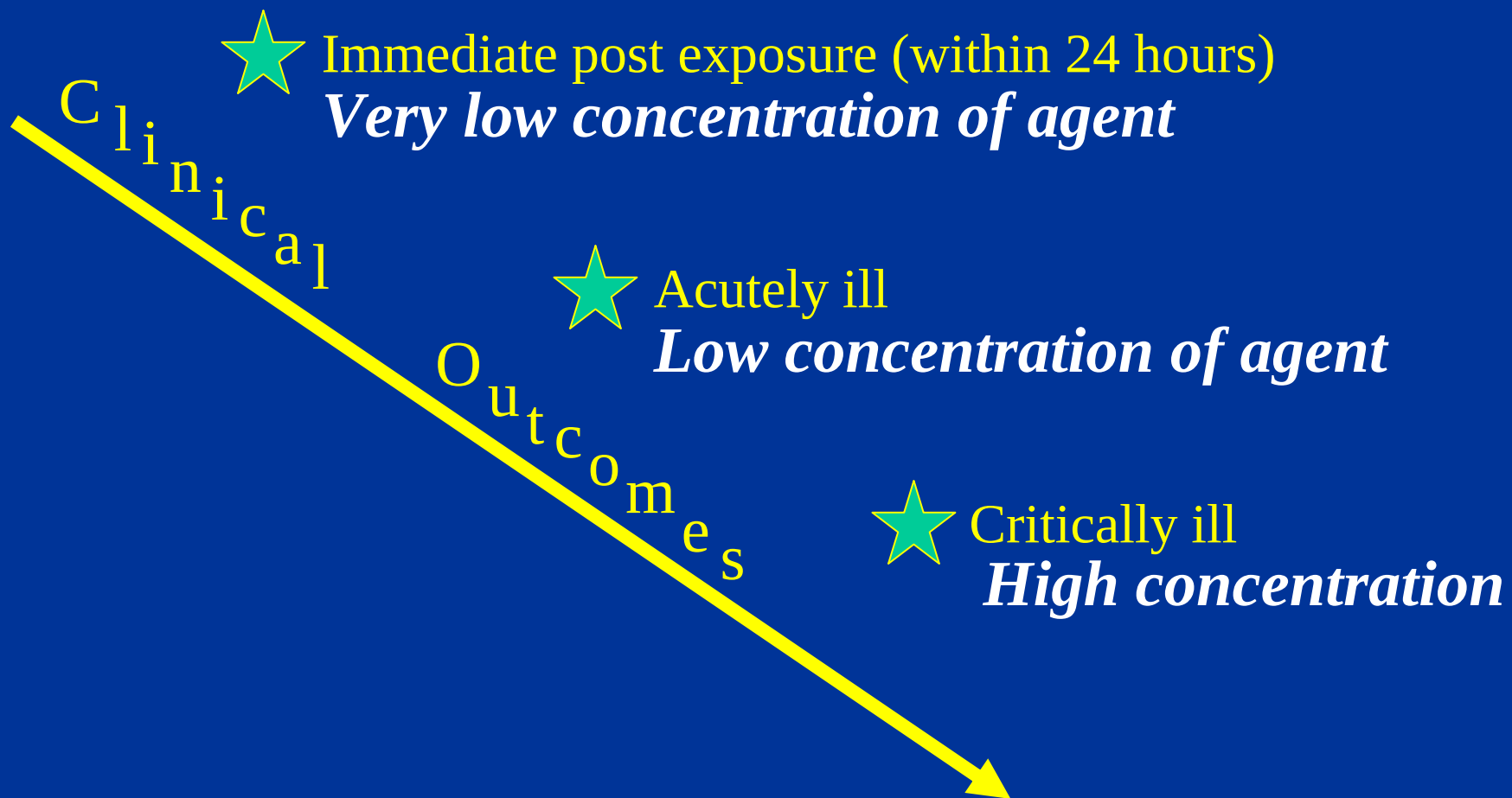

Benefits of a Rapid and Sensitive Protocol for PCR Detection of Bacteria in Blood

B. anthracis – Possible Early ID of Mediastinitis

Y. pestis – Early ID of Septicemic Bubonic Patients

F. tularensis – Slight Aerosol Danger from Culture

Brucella spp. – Blood Culture Can Take 4 Weeks

Burkholderia – Rapid Onset of Sepsis Within 24-48 h

U.N. Report: Efforts to Halt al-Qaida Failing Terror Group Moving Toward Chemical, Biological Weapons By EDITH M. LEDERER, AP

UNITED NATIONS (Nov. 15) - The al-Qaida terror network is determined to use chemical and biological weapons and is restrained only by the technical difficulties of doing so, a U.N. expert panel said in a confidential report.

Date: 14 Nov 2003

From: ProMED-mail <promed@promedmail.org>

Source: News 24 website / Xinhua [edited]

Food laced with rat poison killed one and sent 24 to the hospital during a wedding party in east China's Jiangsu province, state media said on Friday.

The mass poisoning happened on Thursday in Sihong county, as several guests started feeling sick shortly after the wedding banquet had started, the Xinhua news agency reported.


Typical Clinical Specimen Collection Time Line

Time Period	<i>Bacillus anthracis</i>	<i>Yersinia pestis</i>
0 – 24 h	nasal swab	nasal swab
 72 h	serum	serum sputum


Performance Specifications Objectives

- **Field Portable**
- **Closed System Architecture**
- **Minimal End User Input / Fully Automated**
- **RNA and DNA from Bacteria, Virus, and Host Response Using 1 – 8 ml of Whole Blood**
- **Protection of Host RNA as well as Analyte RNA, DNA and Protein**
- **Minimize PCR Inhibitory Host DNA and Blood Matrix Biomass**
- **Rapid Processing Time (analyte capture and elution within 15 minutes / total nucleic acid extraction within 45 minutes)**
- **Plug and Play Compatibility with All Nucleic Acid, Biosensor, Mass Spec, or Culture Based Detection Systems, and Automated Host Response Array Technologies**

Challenge for Rare Analyte PCR Detection in Clinical Samples

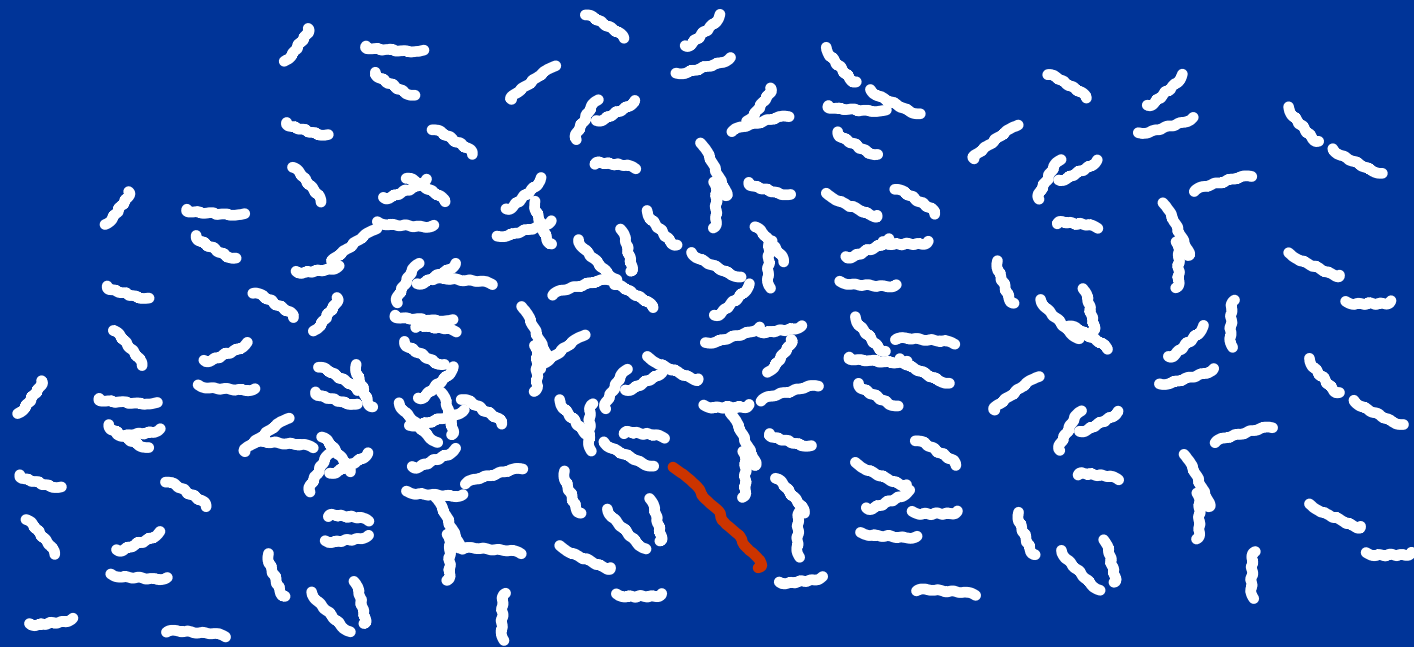


B. anthracis *Gram stain*

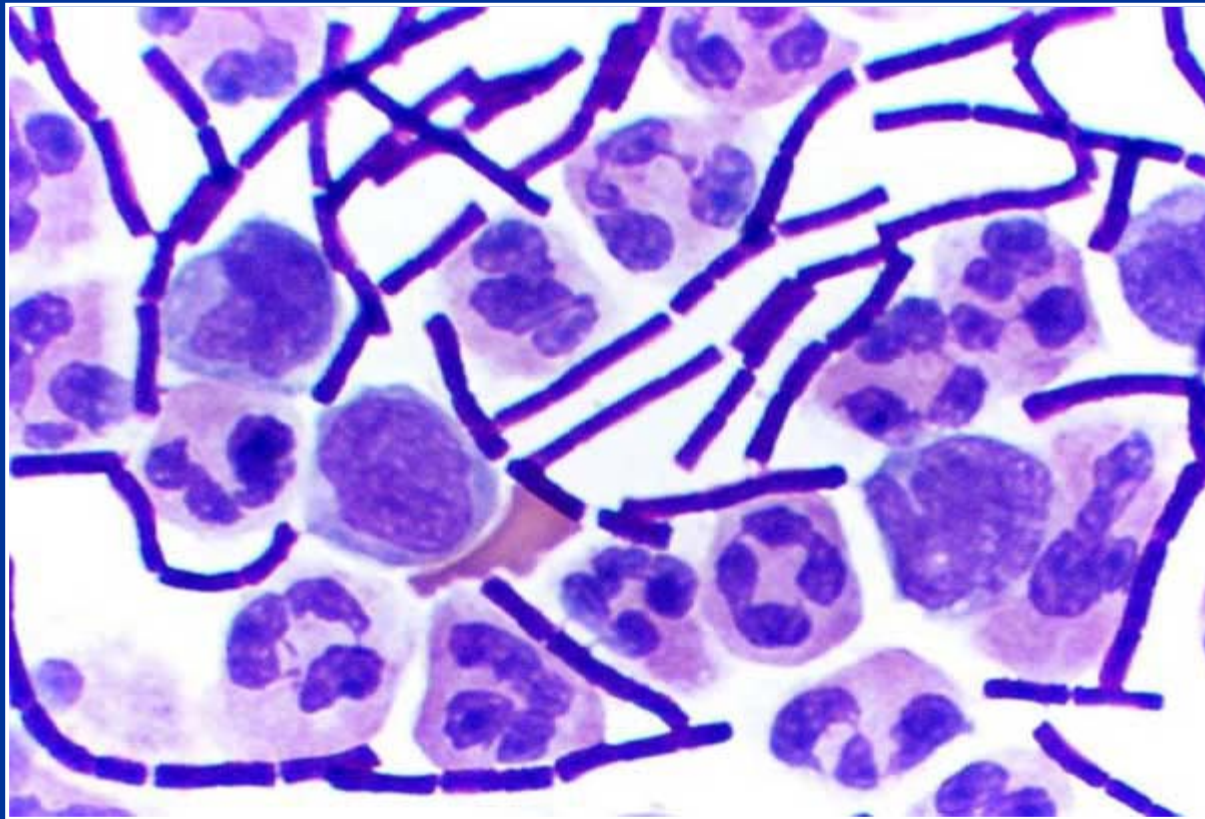

Light Cycler Crossing Points for *Burkholderia pseudomallei*


Sample Number	CFU / 6 ml Whole Blood	Crossing Point
1	0.1	0.00
2	0.1	0.00
3	5.0	34.10
4	5.0	37.77
5	10.0	36.19
6	10.0	37.27
7	20.0	38.69
8	20.0	34.96
9	40.0	37.00
10	40.0	34.34

**Laboratory of Dr. Tim Inglis and PCR
Reagents from USAMRIID**

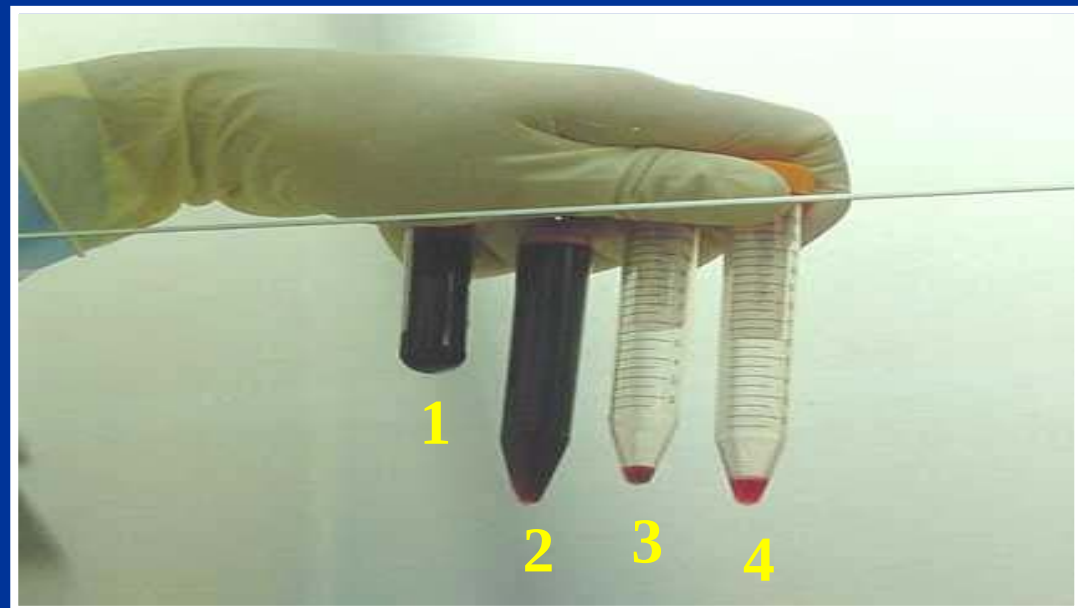


**Prep Method: 1) spiked blood, 2) diluted blood,
3) pellet, 4) enzyme treated pellet.**


DNase Induced Protein Solubility

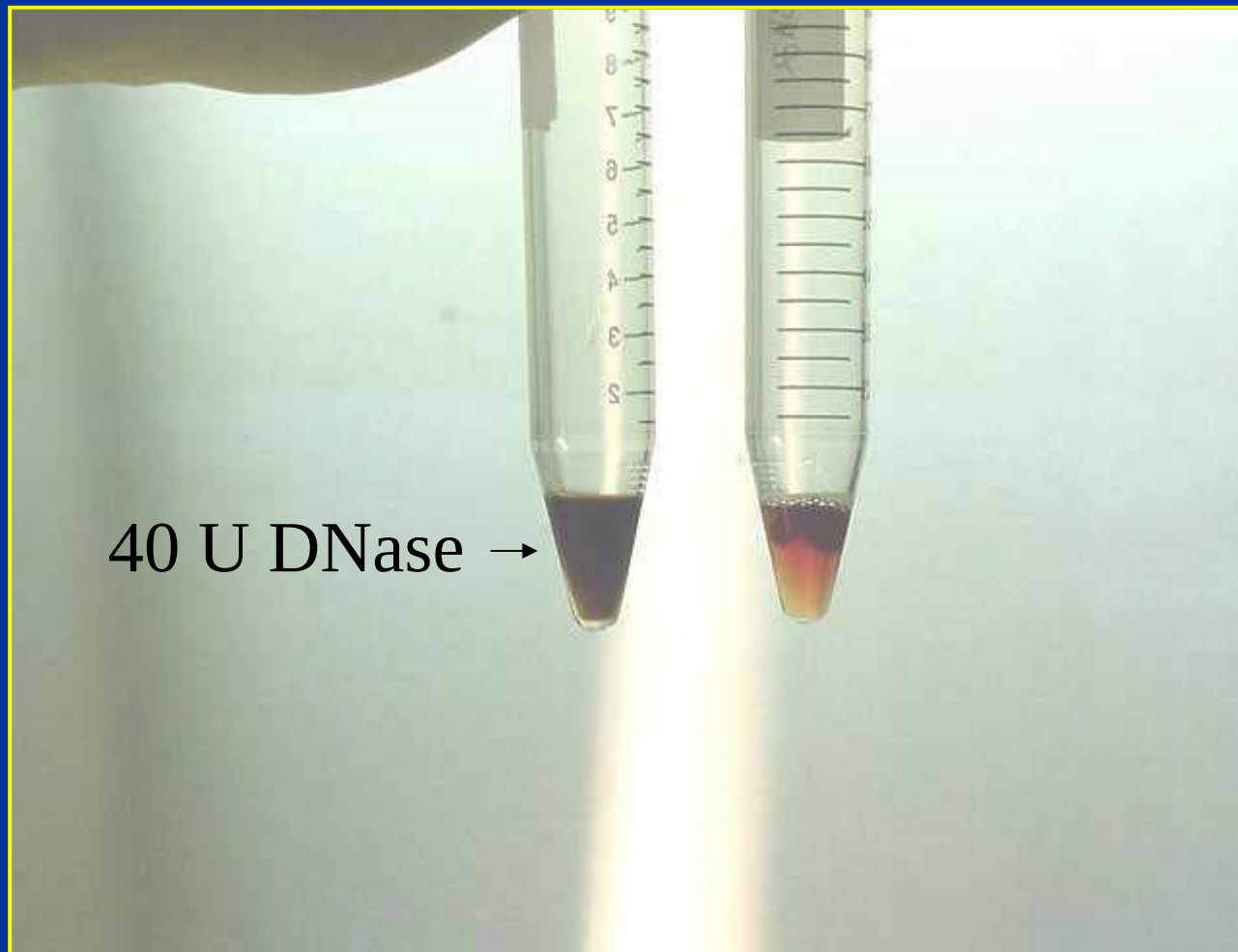

DNase Induced Protein Solubility

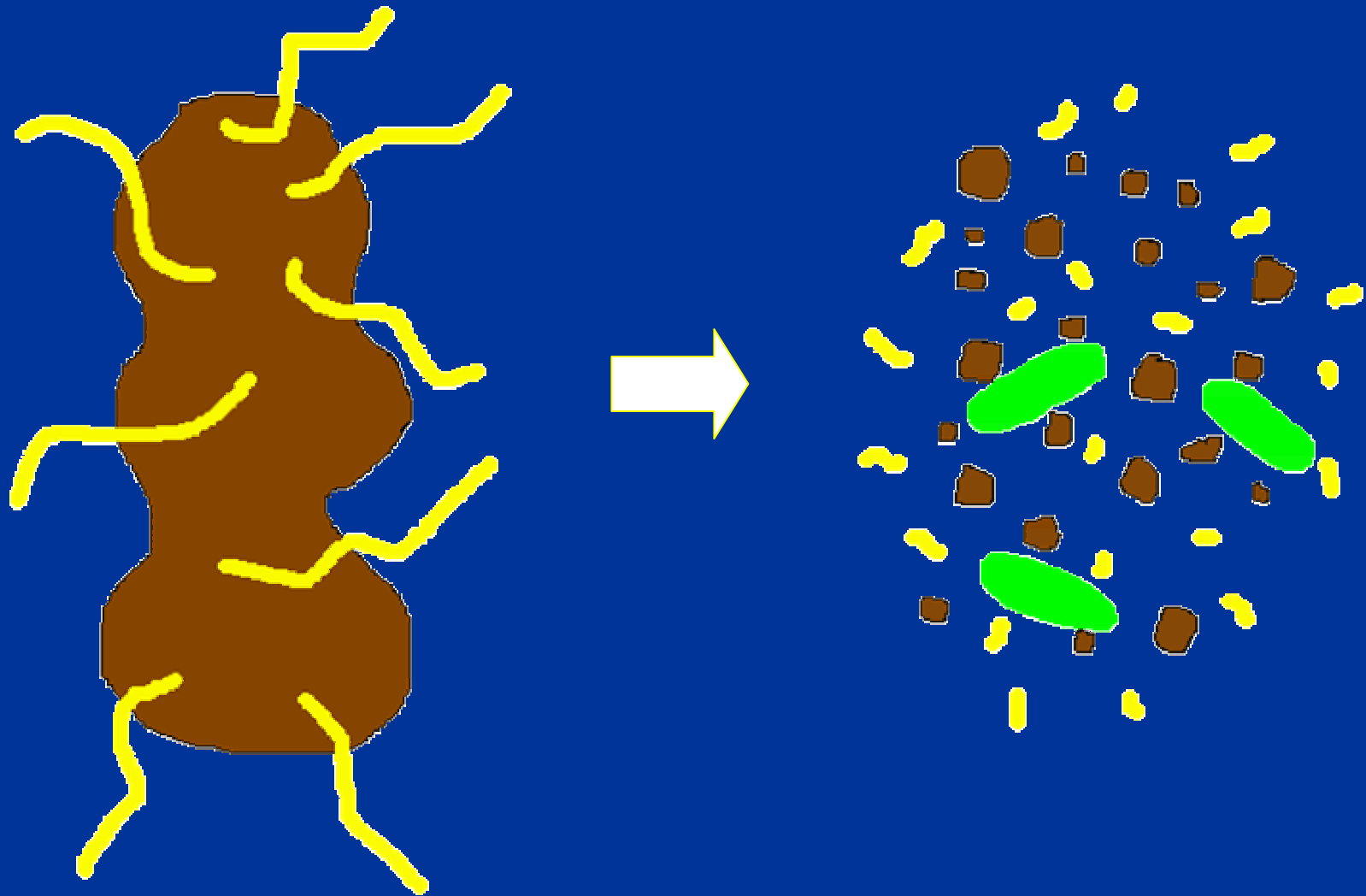

Sample Lysis and Nucleic Acid Extraction

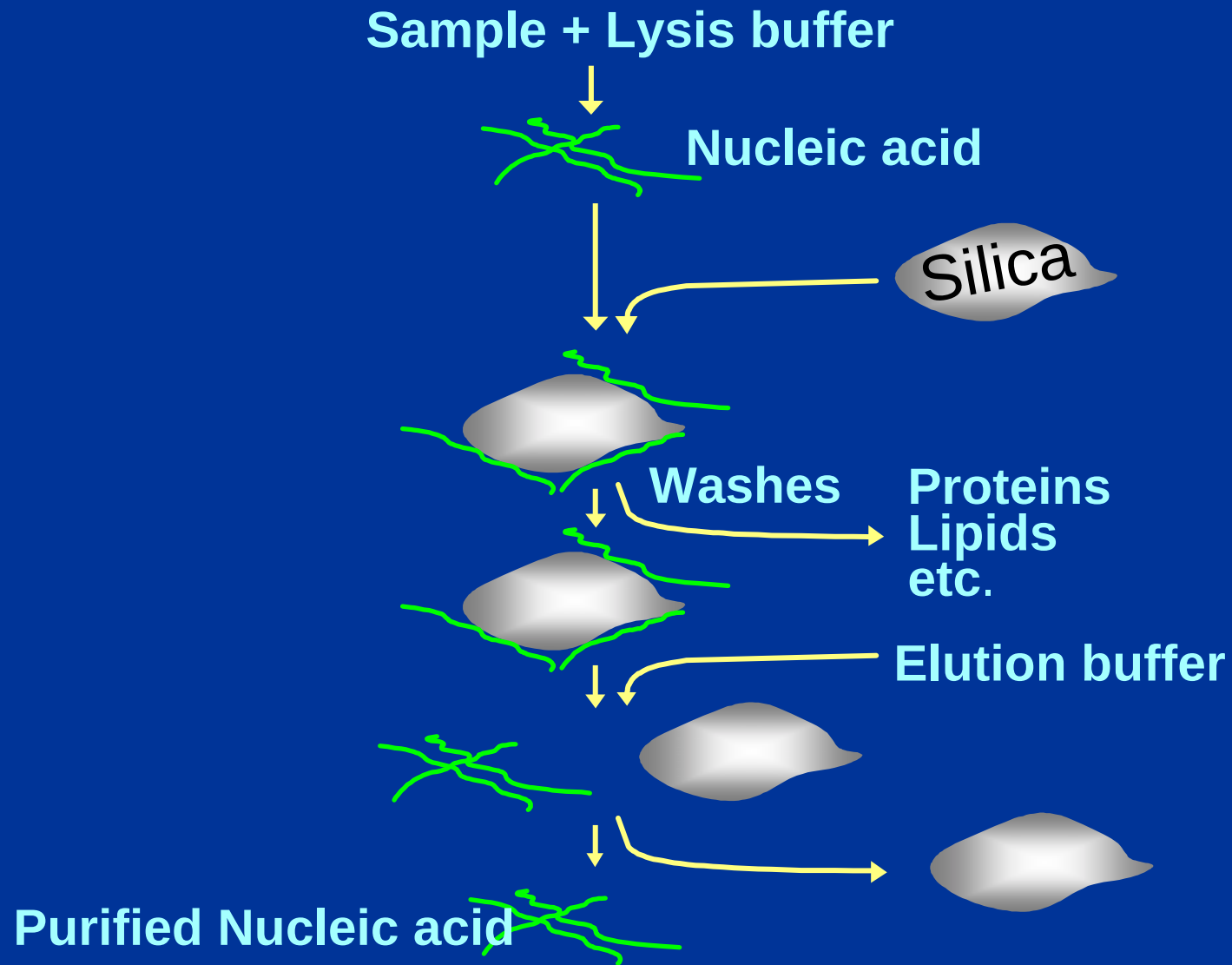

Whole Blood Seeded with 45 - 75 CFU *Brucella suis* / 6ml


Light Cycler Crossing Points

Spiked Conc. (cfu/6ml)	Sample I.D.*	Oligo Set BRU 1	Oligo Set BRU 2	Oligo Set BRU 3
45	1	38.69		
	2	39.09	42.41	
	3		42.26	
	4			38.12
	5	39.66	37.8	39.56
	6	38.66	39.71	35.67
60	7	37.97	37.36	
	8	41.77	38.36	
	9	39.06	38.06	38.91
	10	42.31	38.85	
	11	39.63	38.91	
	12	39.31	36.32	35.67
75	13	41.28	38.21	
	14	44.72	36.84	
	15	38.57	37.41	38.61
	16	39.03	37.82	38.06

Crossing Points for *B. suis* Samples Processed with and with out Human DNA Removal (HDR)


CFU / 6 ml Blood	Sample Number	No HDR Treatment	With HDR Treatment
20	1	0.00	0.00
20	2	0.00	0.00
50	3	42.11	36.68
50	4	0.00	0.00
100	5	40.19	36.88
100	6	0.00	39.24
200	7	44.17	38.58
200	8	0.00	36.91
400	9	40.45	34.48
400	10	43.29	35.67

Only 50 % of the untreated spiked samples produced positive signals. Ninety percent of the HDR treated spiked samples produced positive signals. The HDR treated samples produced earlier crossing points and 50 to 100 % higher fluorescence signals.

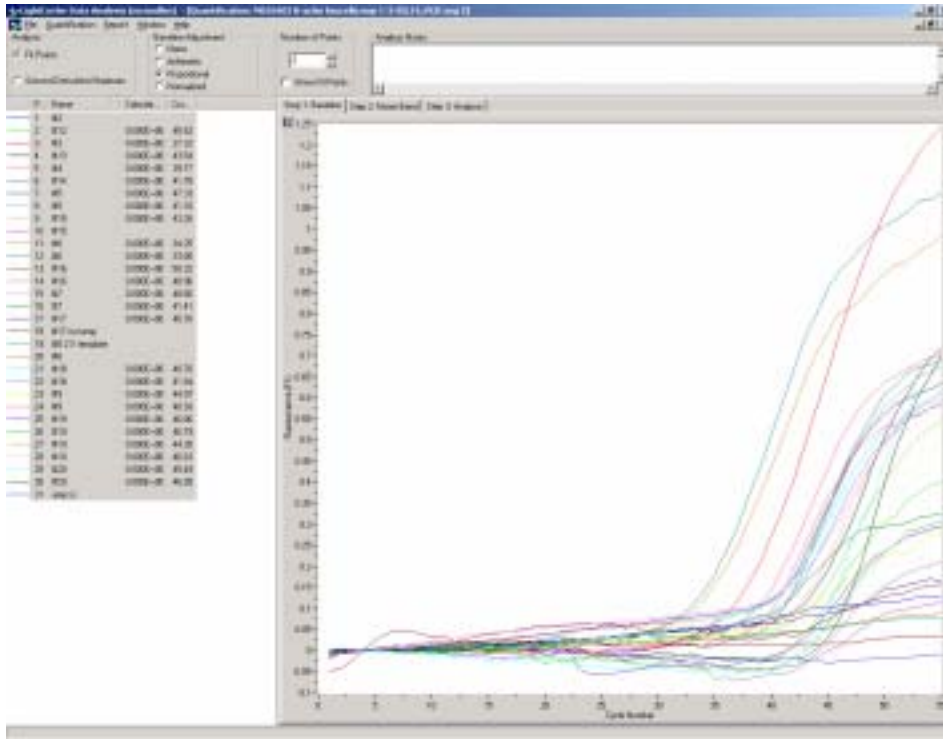

Right. Human DNA Amplification with CDC's Beta actin oligo set using CBD BB005.0 Protocol.

The crossing points ranged from 20.41 to 21.92 with a standard deviation of .052 and a C.V. of 2.42. Note the narrow range of total fluorescence values.

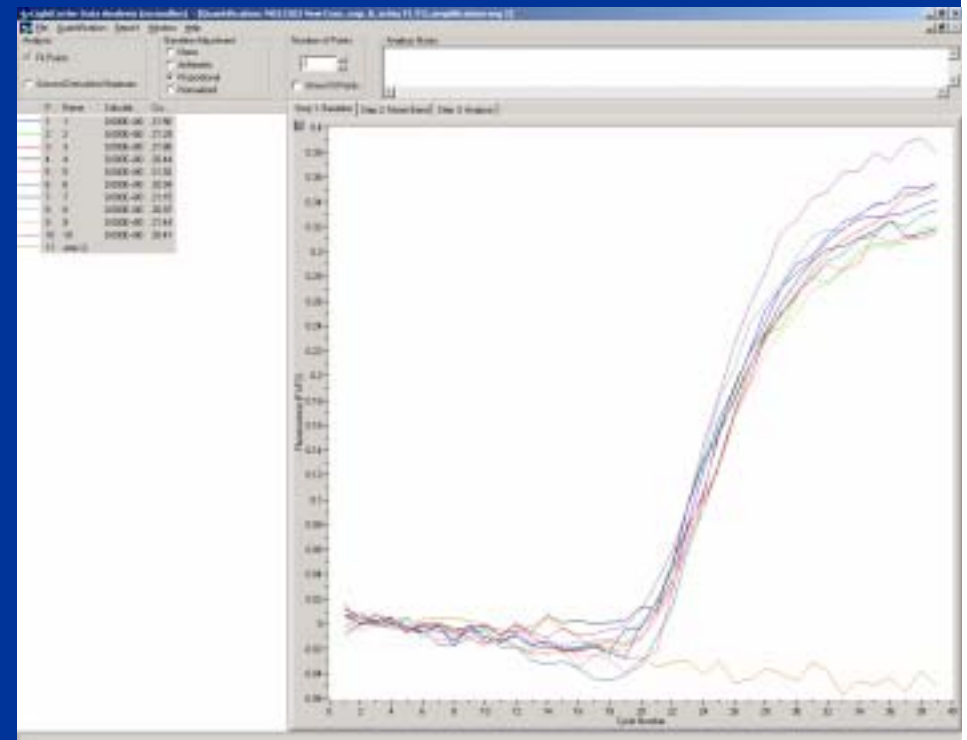

Whole Blood Seeded with 10-80 CFU *Yersinia pestis* / 6ml

Light Cycler Crossing Points

Spiked Conc. cfu/6ml	Sample I.D.	YP2 Event A			YP2 Event B			YP9 Event A			YP9 Event B		
		AVG.	S.D.	C.V.	AVG.	S.D.	C.V.	AVG.	S.D.	C.V.	AVG.	S.D.	C.V.
1	1 2												
10	3 4							42.90 41.99	0.82 0.87	1.91 2.07	41.90 42.26	0.35 0.49	0.84 1.16
20	5 6	37.98 37.04	1.78 0.30	4.69 0.81	38.09 39.01	0.62 1.93	1.63 4.95	41.19 42.18			42.59 42.89	0.22 1.33	0.52 3.10
40	7 8	36.91 38.51	0.35 0.08	0.10 0.21	36.94 37.33	0.13 0.37	0.35 0.99	41.03 42.39	0.78 0.42	1.90 0.99	41.80 40.94	0.59 0.31	1.41 0.76
80	9 10	38.01 36.61	0.05 0.40	0.13 1.10	36.63 36.87	0.27 0.14	0.74 0.38	40.25 40.18	0.11 0.40	0.27 1.00	39.91 39.39	0.21 0.48	0.53 1.22
Spiked Conc. cfu/6ml	Sample I.D.	YP12 Event A			YP12 Event B			YP16 Event A			YP16 Event B		
		AVG.	S.D.	C.V.	AVG.	S.D.	C.V.	AVG.	S.D.	C.V.	AVG.	S.D.	C.V.
1	1 2												
10	3 4	35.82 36.00			36.48 35.82	0.23 0.06	0.63 0.17	37.30 35.23	0.63 0.24	1.69 0.68	36.10 37.45	0.40 0.16	1.11 0.46
20	5 6	35.97 36.15	0.21 0.31	0.58 0.86	36.43 36.47	0.45 0.14	1.24 0.38	36.35 35.17	1.30 0.36	3.58 1.02	36.17 34.10	0.64 0.20	1.77 0.59
40	7 8	35.46 36.42	0.27 0.15	0.76 0.41	35.88 35.68	0.09 0.06	0.25 0.17	33.11 33.74	0.26 0.04	0.79 0.12	33.97 35.22	0.23 0.51	0.68 1.45
80	9 10	35.72 35.15	0.12 0.31	0.34 0.88	35.32 34.78	0.32 0.10	0.91 0.29	33.19 33.58	0.02 0.41	0.06 1.22	33.37 32.60	0.33 0.25	0.99 0.77

Dilution Tubes for Vegetative *B. anthracis*:
Used to Achieve S.D. < 5 % on Plate Counts



Whole Blood Seeded with ≤ 10 CFU *B. anthracis* / 6ml

Sample #	CDC pX02 Oligo Set Avg. Noise Band Crossing Point	CDC pX02 Oligo Set Standard Deviation	CDC Genomic Oligo Set Avg. Noise Band Crossing Point	CDC Genomic Oligo Set Standard Deviation
1				
2				
3	37.86	0.93	40.39	
4	39.07		40.00	
5				
6				
7			40.77	
8			39.86	
9	38.42	0.54	38.53	0.44
10	36.90	0.27	37.75	0.11
11				
12			39.10	
13			38.61	
14				
15				
16	37.79	0.80	38.39	
17	38.27	0.37	38.68	1.20
18	37.52		39.60	0.63
19	36.19	0.36	37.47	
20	37.56	0.73	38.80	0.10
21	36.71	0.44	38.74	1.02
22	37.09		38.47	
23	38.07	0.80	39.06	
24	37.36		38.37	0.55
25			40.24	
26	38.16	0.73	38.34	0.31
27	37.15	0.30	38.28	0.54
28	39.44	0.29	38.31	0.29
29	37.42	0.33	37.69	0.29

PES Filter Mechanism




Concentration of *B. anthracis* via Filtration and Centrifugation Followed by PCR Analysis

Centrifugation				Filtration			
Amount <i>B. anthracis</i> Dev. Seeded (cfu/6ml)	Noise Band Crossing Points	Mean	Std. Dev.	Noise Band Crossing Points	Mean	Std.	
< 0.01							
< 0.01							
< 10.0				40.33	39.89		40.11
< 10.0					37.79		37.79
< 10.0				40.36		37.69	39.03
< 10.0	41.93	40.31	41.12				
< 10.0		40.47	40.47	37.90	37.70	37.79	37.80
< 10.0	38.11	40.36	39.24	36.45	36.09	36.81	36.45
							0.10
							0.36
< 50.0	37.15	36.14	37.98	37.23	37.81		
< 50.0	36.45	38.15	38.49	37.70	1.09	35.24	34.18
						34.68	34.70
							0.82
							0.53

Filtration Method May Be Automated

Six ml Whole Blood With and Without Detergent and Enzyme Treatment Followed by Centrifugation @ 5,500 RCF

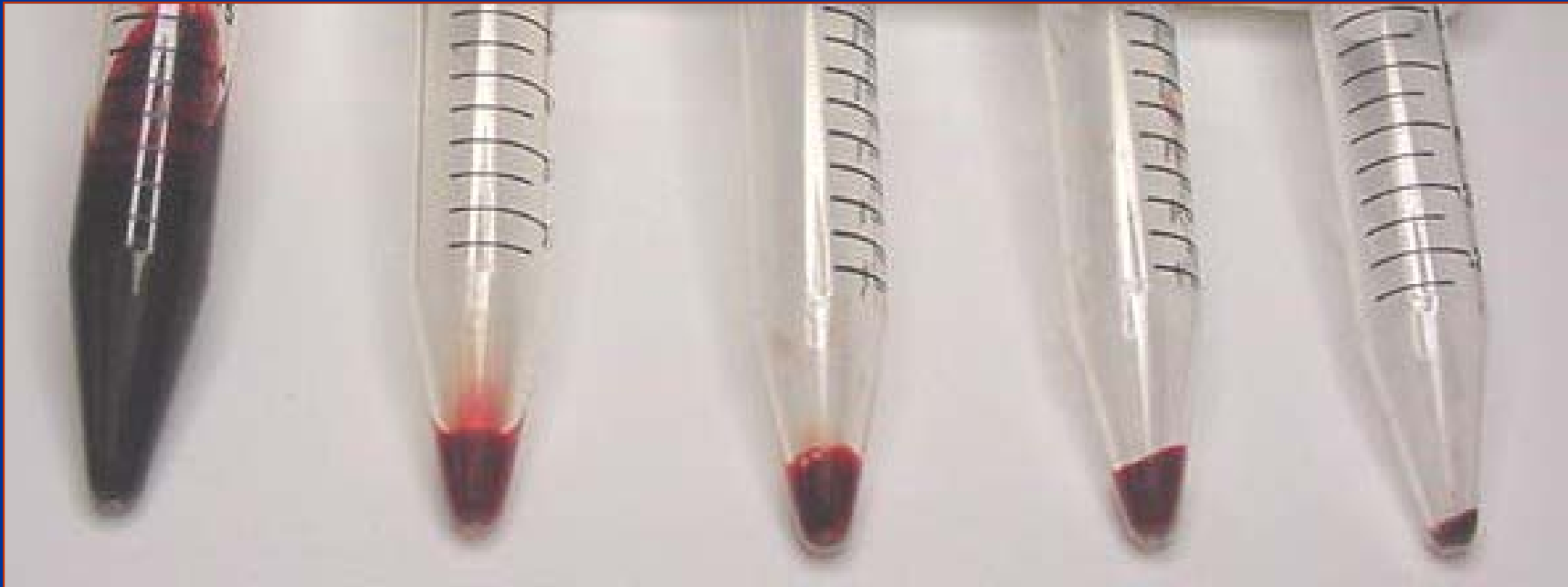

Untreated

1 %

1 % Triton

Exp.

Exp.

Blood

Saponin


X – 100

Detergent


Detergent

+ Enzymes


Successful CRP dose response for reagent processed whole blood

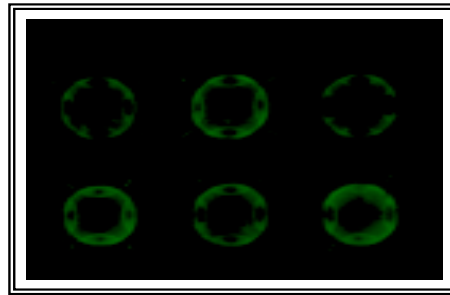


0ng/ml

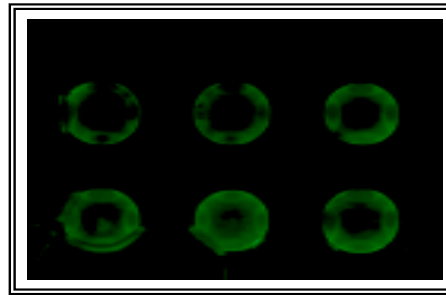


1ng/ml

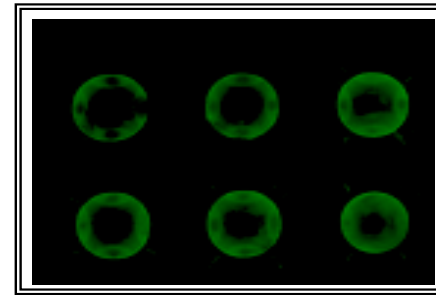


10ng/ml

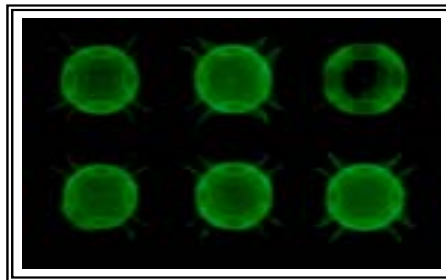

100ng/ml

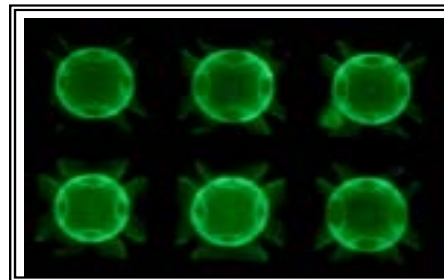

1000ng/ml


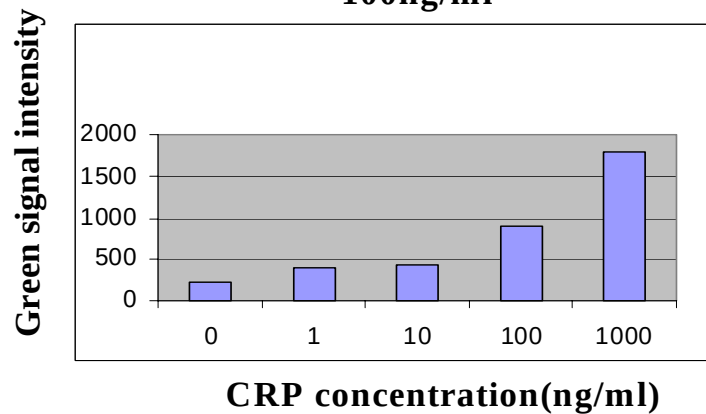


*Data acquired at the University of Texas at Austin by the McDevitt laboratory on 6-30-03;
Sample provided by Matt Ewert,
USF*

CANARY Sensor Concept

MIT Lincoln Laboratories

- Responds in seconds
- High sensitivity
- B-cell specificity is engineered with recombinant antibodies
- Aequorin amplifies signal and provides simple interface

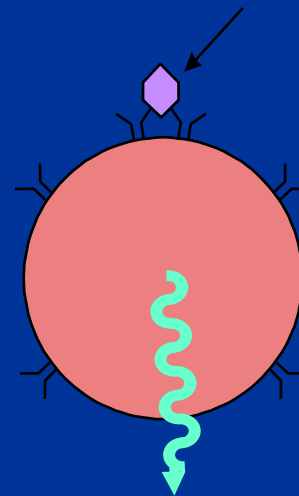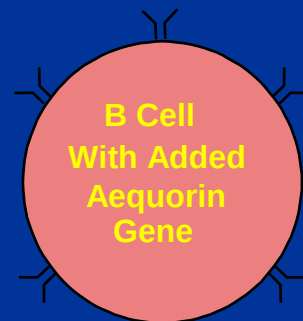


(1) Pathogens cross link membrane-bound recombinant antibodies

(2) Biochemical signal amplification releases Ca^{2+}


(3) Ca^{2+} makes aequorin emit photons

(4) Detect photons from cells


Analyte Capture

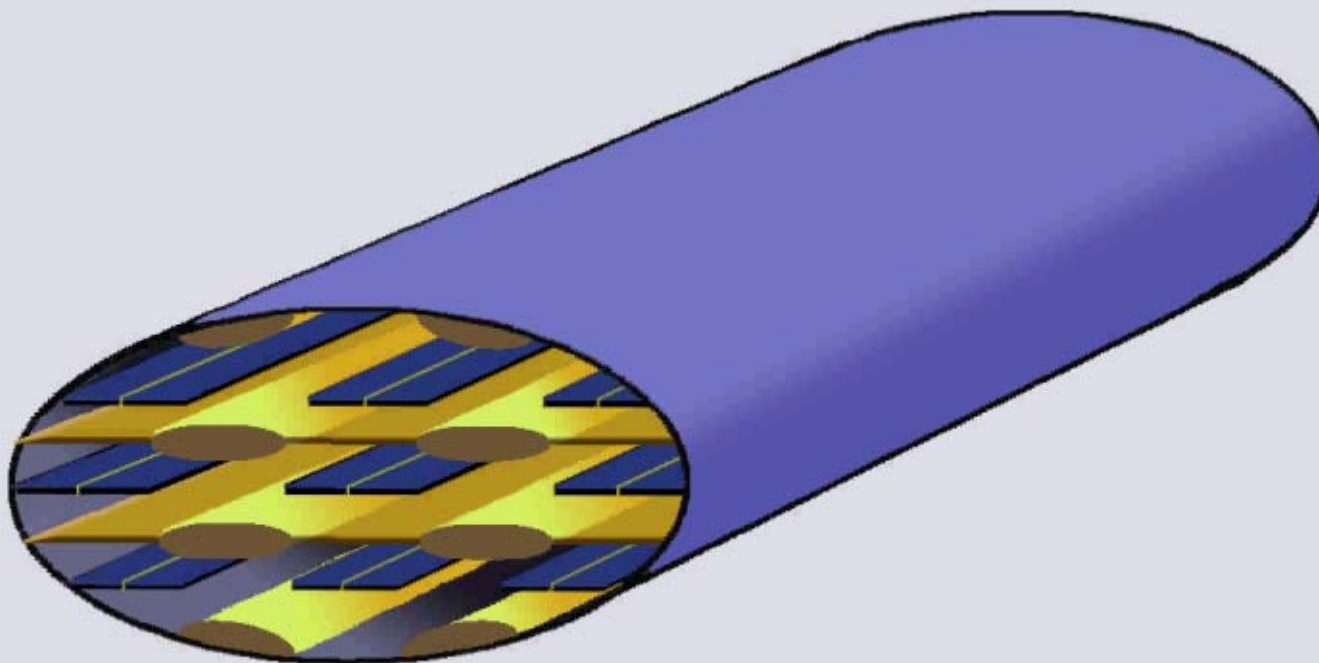

Analyte Elution

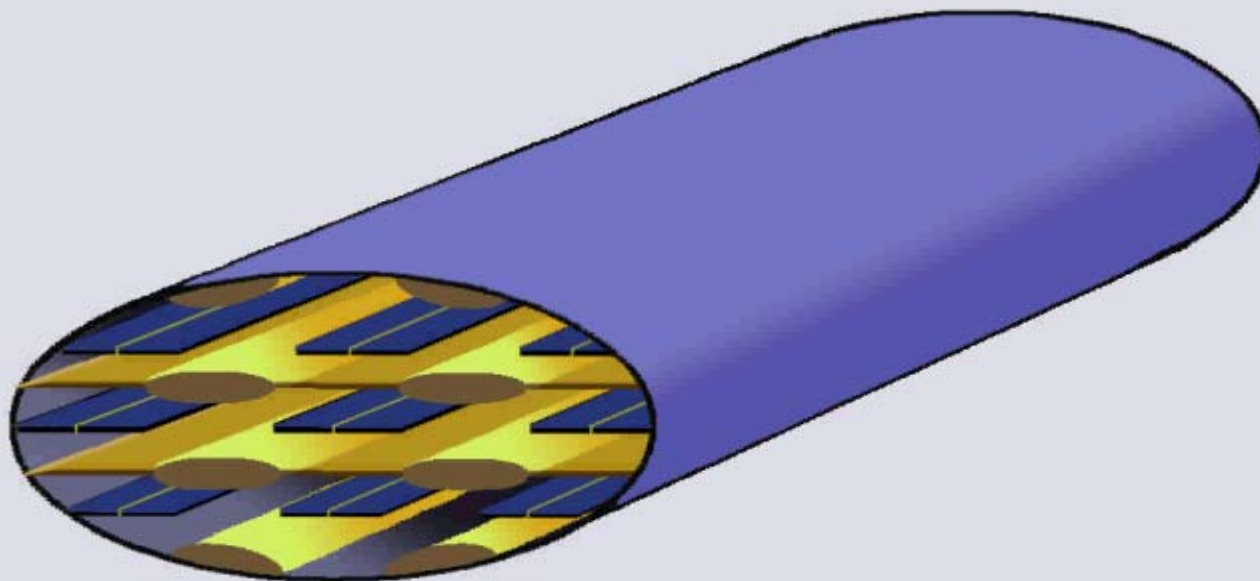


Crude Prototype Platform for Developing Capture Chemistries

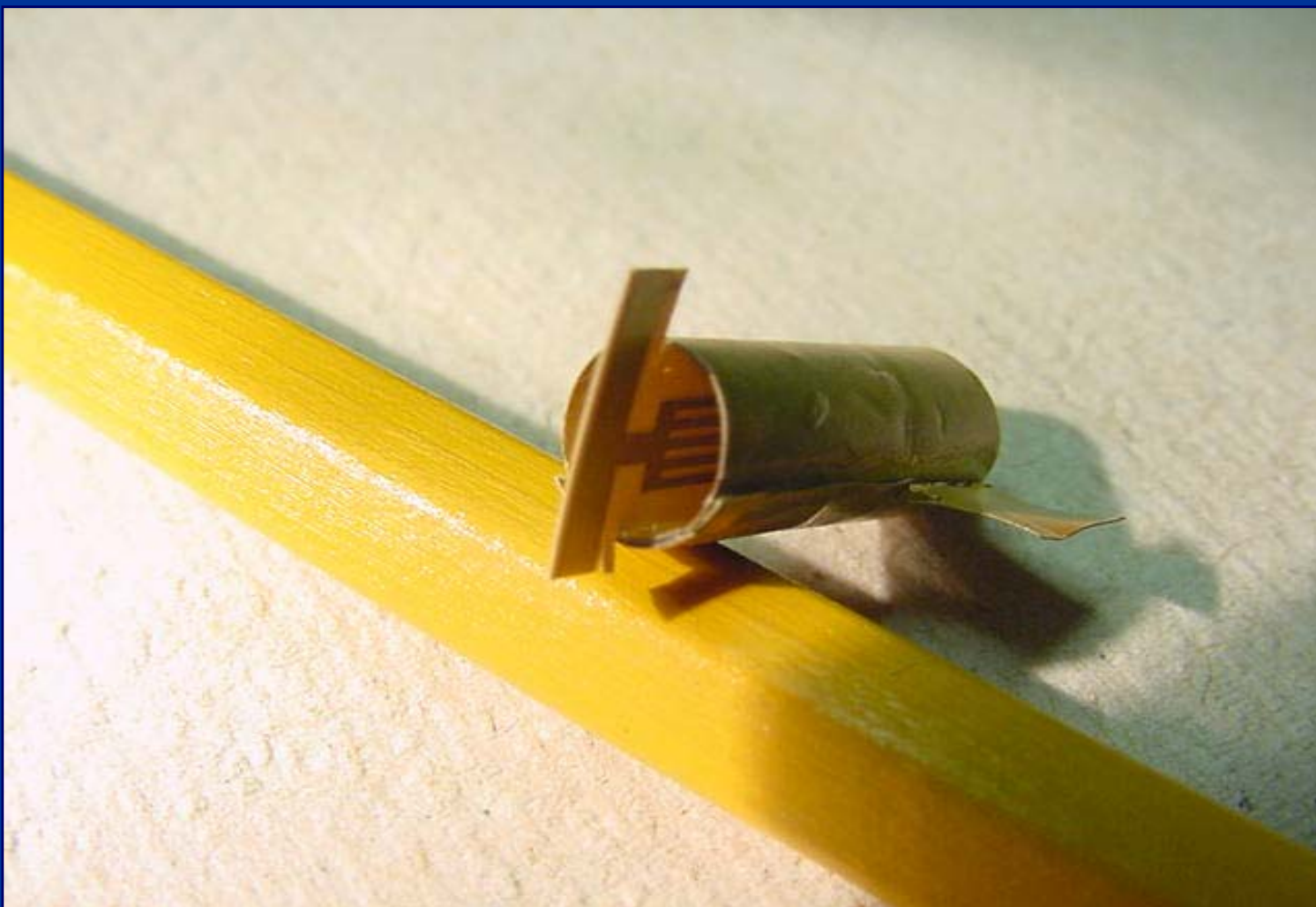




Blow Up View

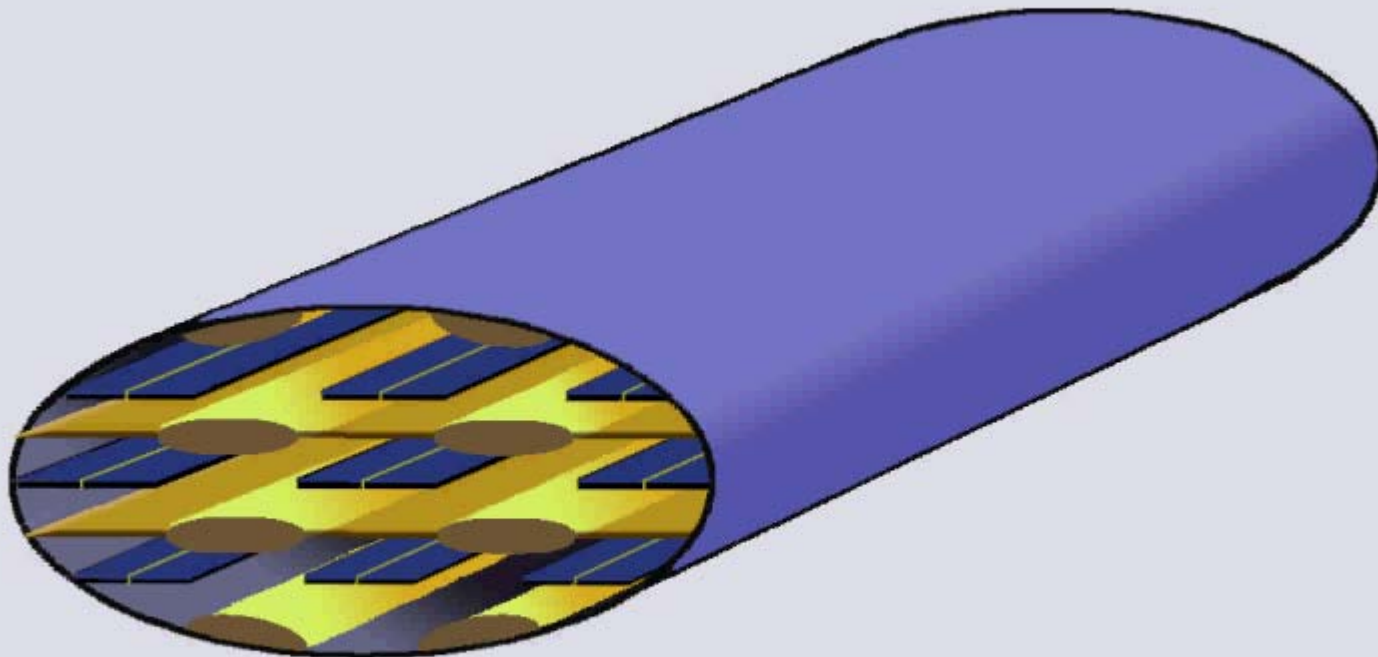

Analyte Elution: Chemical

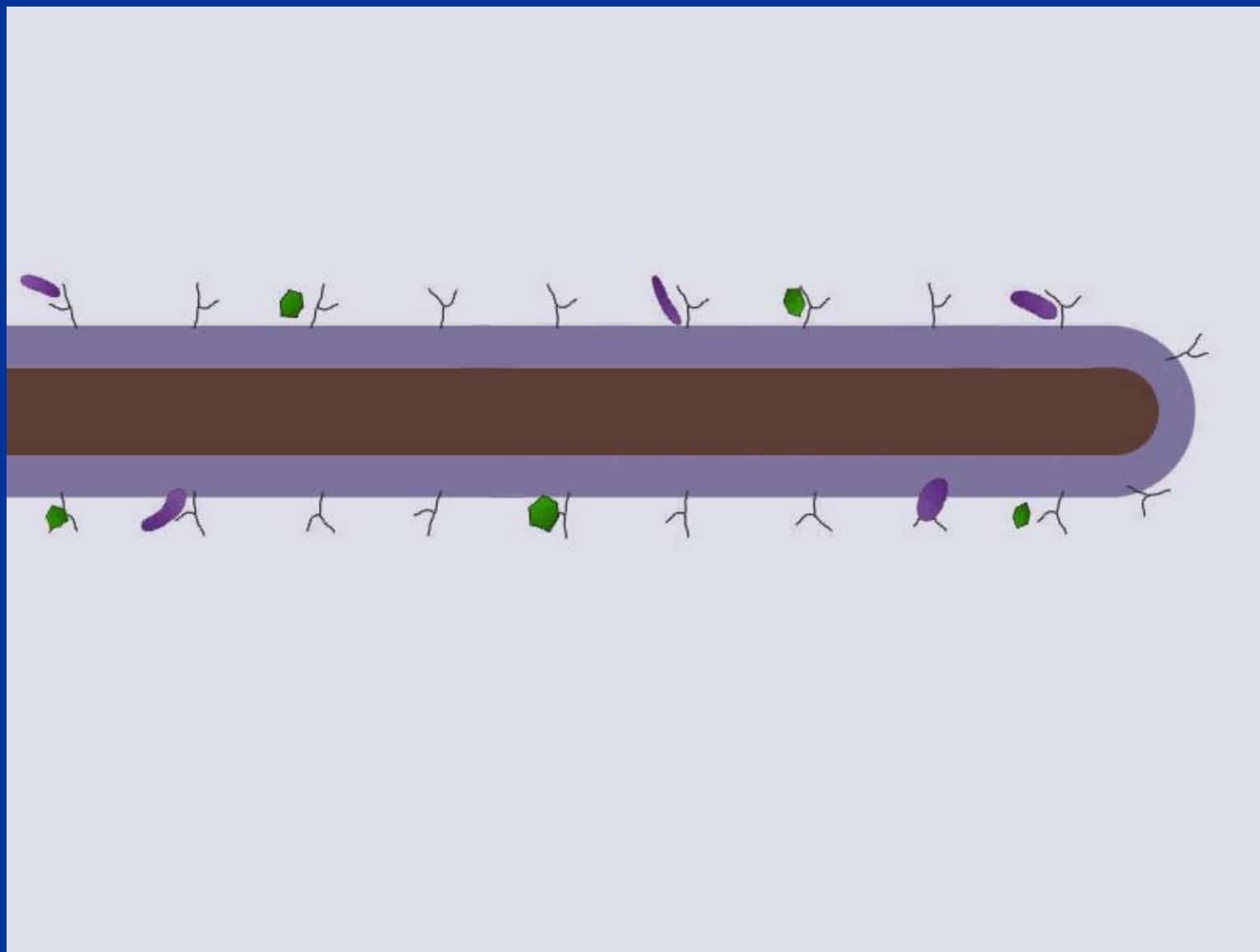

Florida Department of Health Paula Mooty and Selena Grant



Acknowledgements

**Philip Amuso, Andrew Cannons,
Jacqueline Cattani, Selena Grant, Tim
Inglis, Phil Lee, John McDevitt, Patrick
McMullen, Rich Meyer, Paula Mooty,
Dave Norwood, and Martha Petrovick**