



EDGEWOOD

CHEMICAL BIOLOGICAL CENTER

U.S. ARMY RESEARCH, DEVELOPMENT AND ENGINEERING COMMAND

ECBC-TR-749

MALVERN NANO ZS PARTICLE SIZE COMPARISON WITH THE INTEGRATED VIRUS DETECTION SYSTEM (IVDS)

Charles H. Wick

RESEARCH AND TECHNOLOGY DIRECTORATE

Patrick E. McCubbin



OptiMetrics, Inc.
Research & Engineering

OPTIMETRICS, INC.
Abingdon, MD 21009

March 2010

Approved for public release;
distribution is unlimited.



ABERDEEN PROVING GROUND, MD 21010-5424

Disclaimer

The findings in this report are not to be construed as an official Department of the Army position unless so designated by other authorizing documents.

REPORT DOCUMENTATION PAGE

Form Approved
OMB No. 0704-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Department of Defense, Washington Headquarters Services, Directorate for Information Operations and Reports (0704-0188), 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302. Respondents should be aware that notwithstanding any other provision of law, no person shall be subject to any penalty for failing to comply with a collection of information if it does not display a currently valid OMB control number. **PLEASE DO NOT RETURN YOUR FORM TO THE ABOVE ADDRESS.**

1. REPORT DATE (1-January-2005) XX-03-2010		2. REPORT TYPE Final		3. DATES COVERED (From - To) Aug 2009 - Sep 2009	
4. TITLE AND SUBTITLE Malvern Nano ZS Particle Size Comparison with the Integrated Virus Detection System (IVDS)				5a. CONTRACT NUMBER	
				5b. GRANT NUMBER	
				5c. PROGRAM ELEMENT NUMBER	
6. AUTHOR(S) Wick, Charles H. (ECBC); and McCubbin, Patrick E. (OMI)				5d. PROJECT NUMBER None	
				5e. TASK NUMBER	
				5f. WORK UNIT NUMBER	
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) DIR, ECBC, ATTN: RDCB-DRD-D, APG, MD 21010-5424 OMI, 100 Walter Ward Boulevard, Abingdon, MD 21009				8. PERFORMING ORGANIZATION REPORT NUMBER ECBC-TR-749	
9. SPONSORING/MONITORING AGENCY NAME(S) AND ADDRESS(ES)				10. SPONSOR/MONITOR'S ACRONYM(S)	
				11. SPONSOR/MONITOR'S REPORT NUMBER(S)	
12. DISTRIBUTION/AVAILABILITY STATEMENT Approved for public release; distribution is unlimited.					
13. SUPPLEMENTARY NOTES					
14. ABSTRACT The Malvern Nano ZS particle size and zeta potential equipment was tested for particle size accuracy and compared to the Integrated Virus Detection System. Polystyrene latex standard reference materials were used in this comparison. The materials tested were all in the sub-micron particle size region.					
15. SUBJECT TERMS Particle size Malvern Nano ZS Polystyrene latex Sub-micron Integrated Virus Detection System					
16. SECURITY CLASSIFICATION OF:			17. LIMITATION OF ABSTRACT	18. NUMBER OF PAGES	19a. NAME OF RESPONSIBLE PERSON Sandra J. Johnson
a. REPORT U	b. ABSTRACT U	c. THIS PAGE U			19b. TELEPHONE NUMBER (include area code) (410) 436-2914

20100402023

Blank

PREFACE

The work described in this report was started in August 2009 and completed in September 2009.

The use of either trade or manufacturers' names in this report does not constitute an official endorsement of any commercial products. This report may not be cited for purposes of advertisement.

This report has been approved for public release. Registered users should request additional copies from the Defense Technical Information Center; unregistered users should direct such requests to the National Technical Information Service.

Blank

CONTENTS

1.	INTRODUCTION	7
2.	CONCLUSION.....	8
	LITERATURE CITED.....	17

FIGURES

1.	20 nm PSL (Malvern)	9
2.	20 nm PSL (IVDS).....	9
3.	30 nm PSL (Malvern)	10
4.	30 nm PSL (IVDS).....	10
5.	70 nm PSL (Malvern)	11
6.	70 nm PSL (IVDS).....	11
7.	90 nm PSL (Malvern)	12
8.	90 nm PSL (IVDS).....	12
9.	30 nm Plus 90 nm PSL (Malvern)	13
10.	30 nm Plus 90 nm PSL (IVDS).....	13
11.	30 nm Plus 130 nm PSL by Intensity (Malvern)	14
12.	30 nm Plus 130 nm PSL by Volume (Malvern)	14
13.	30 nm Plus 130 nm PSL, Full Scale (IVDS)	15
14.	30 nm Plus 130 nm PSL, Expanded Scale (IVDS).....	15

MALVERN NANO ZS PARTICLE SIZE COMPARISON WITH THE INTEGRATED VIRUS DETECTION SYSTEM (IVDS)

1. INTRODUCTION

The Malvern Nano ZS measures zeta potential or electrophoretic mobility for fine particle systems in solution, including bacteria and viruses. The characterization of the surface of microorganisms is determined by analyzing the surface charge and layer effects of bacteria and viruses. The system can assess the stability of biological organisms in suspension and emulsions and explore metabolic changes in organisms due to either matrix or environmental variables. Various parameters such as pH, ion interactions, solution composition, and other in vitro parameters can affect the surface charge of the species in question. The non-disruptive analysis leaves organisms intact for further sample analyses.

The system also uses a non-invasive back scatter (NIBS) technology for particle analysis and size distribution. The NIBS technology tests biological materials in solution, using a small sample size, without disrupting the sample.¹

Particle Size Measurements

A series of samples were analyzed for particle size with the Malvern instrument. The following polystyrene latex (PSL) particles, from Bangs Laboratories, were measured: 20 nm, 30 nm, 70 nm, 90 nm, 30 nm plus 90 nm, and 30 nm plus 130 nm. The samples were also measured with the Integrated Virus Detection System (IVDS)² as a comparison.

The Malvern instrument analyzed the 20 nm, 30 nm, 70 nm, and 90 nm particles, and the scans are shown in Figures 1, 3, 5, and 7, respectively. The value for Z-Avg in the table in each figure is the diameter result in nanometers. The IVDS measured the 20 nm, 30 nm, 70 nm, and 90 nm particles, and the scans are shown in Figures 2, 4, 6, and 8, respectively. The resultant size from the IVDS is listed on the graph. The IVDS and the Malvern instrument are in close agreement for the major peak of the PSL standards. The IVDS instrument does show a smaller extra peak in the 70 and the 90 nm size standard that is not shown in the Malvern analysis.

To examine the ability to differentiate several peaks in both instruments, two mixed standards were prepared. A 30 nm plus 90 nm sample was mixed in a 1:3 ratio by volume (where 30 nm=1). The Malvern analysis of the 30 nm plus 90 nm sample is shown in Figure 9. The Z-Avg was determined to be 79.23 nm with only one peak shown. The IVDS detected both standards in Figure 10. Conversations with Malvern technical support indicate that there needs to be an approximate 3X difference in size for their instrument to differentiate peaks.

A 30 nm plus 130 nm sample was mixed in a 10:1 ratio by volume (where 30 nm=10) to see if the Malvern instrument could differentiate the peaks. As shown in Figure 11, the separate peaks of the intensity of the standards are shown. Figure 12 shows the size distribution by volume percent. The IVDS analysis of the standards shows a full scale scan

in Figure 13, and an expanded scale scan in Figure 14. The IVDS instrument was able to show both PSL size standards.

2. CONCLUSION

The Malvern instrument was able to analyze monodispersed particles accurately. The Malvern instrument needs an approximate 3X difference in size to differentiate particles. The size determination in the instrument is by measuring Brownian motion in solution by Dynamic Light Scattering. The 3X difference is a theoretical limit for sizing for this technology.

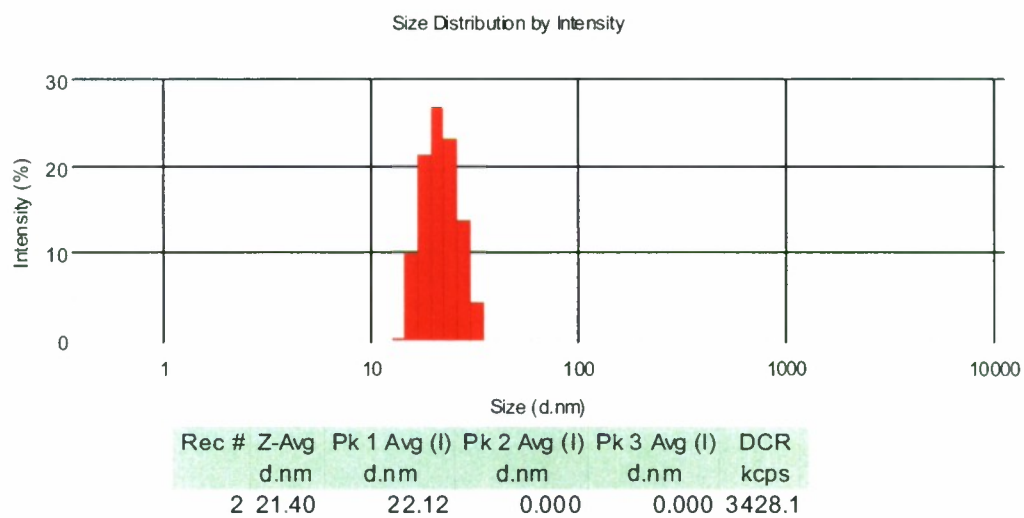


Figure 1. 20 nm PSL (Malvern)

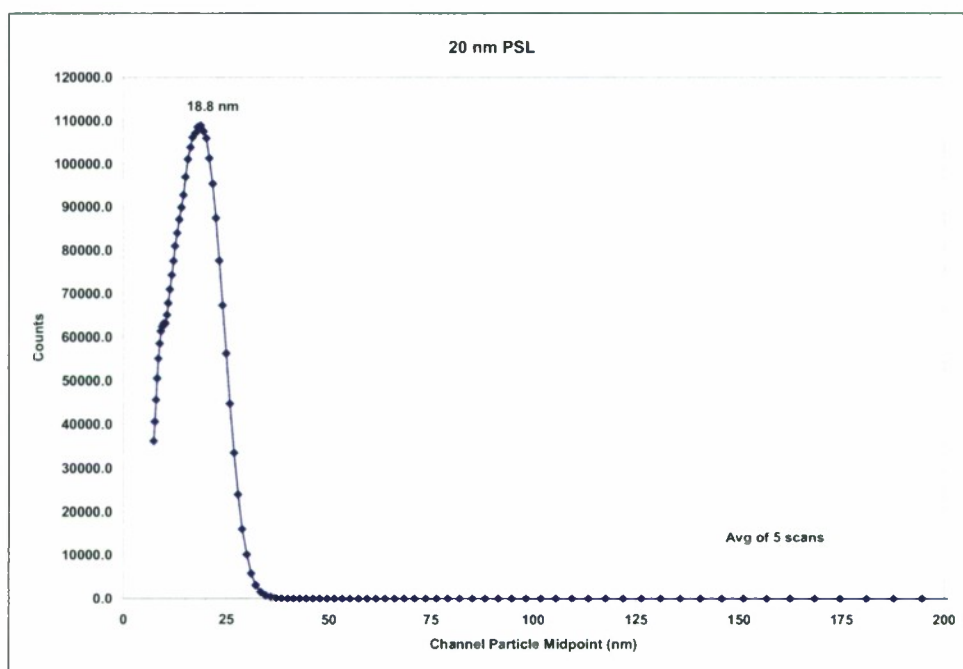


Figure 2. 20 nm PSL (IVDS)

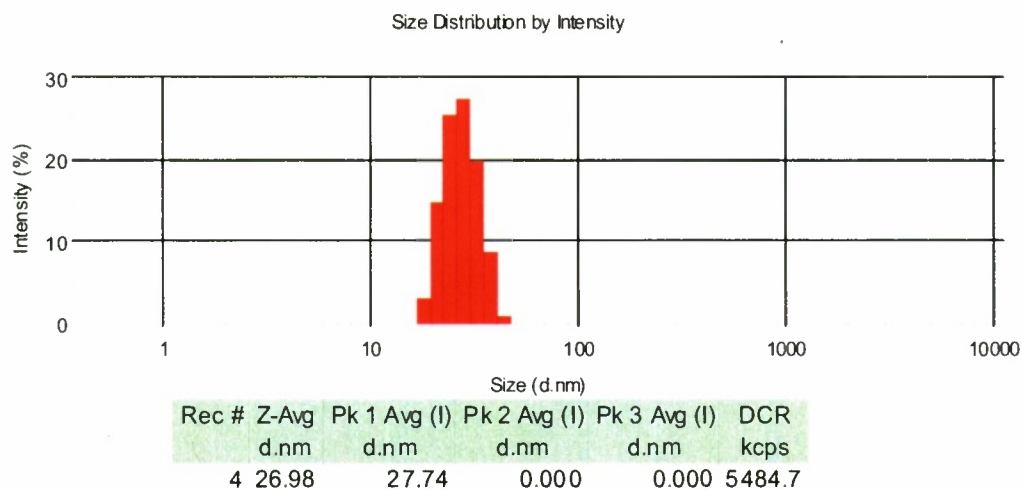


Figure 3. 30 nm PSL (Malvern)

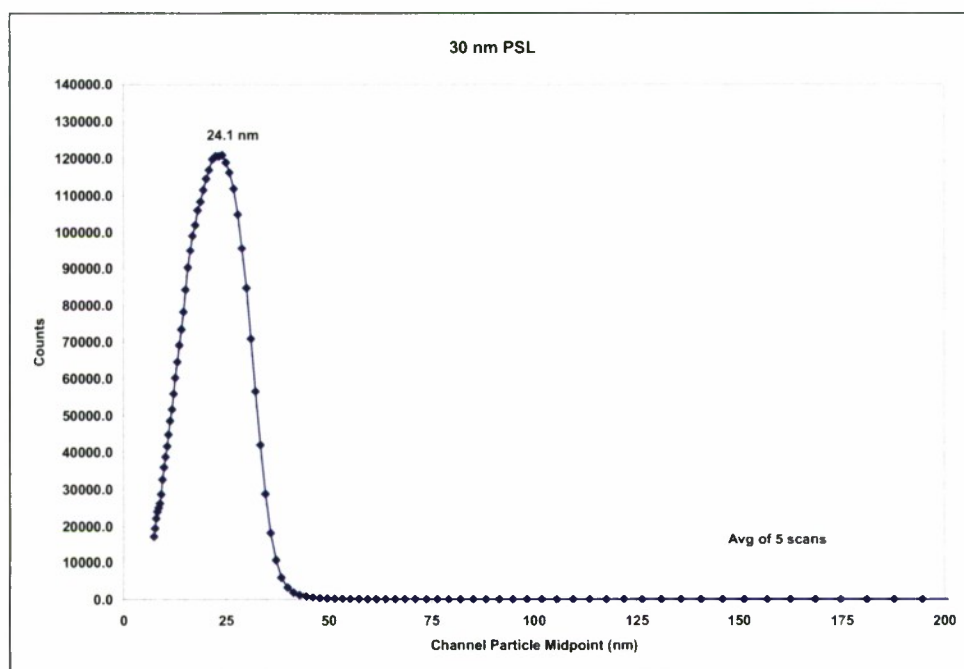


Figure 4. 30 nm PSL (IVDS)

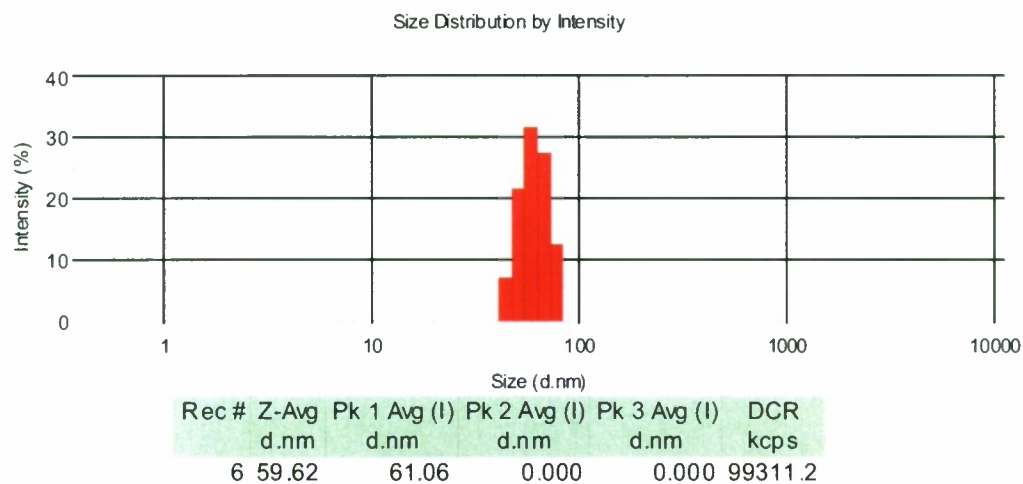


Figure 5. 70 nm PSL (Malvern)

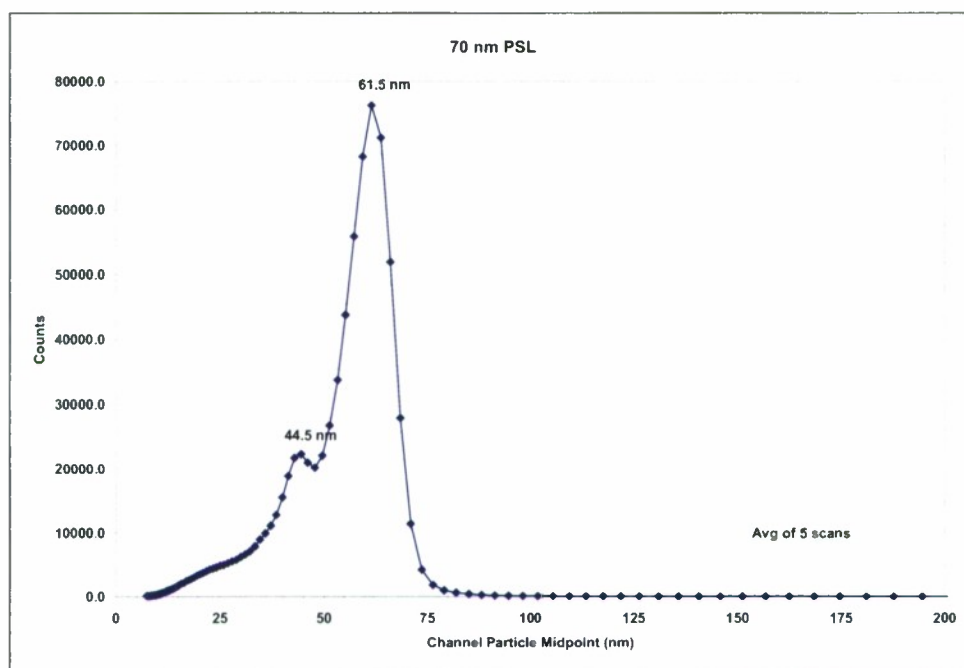


Figure 6. 70 nm PSL (IVDS)

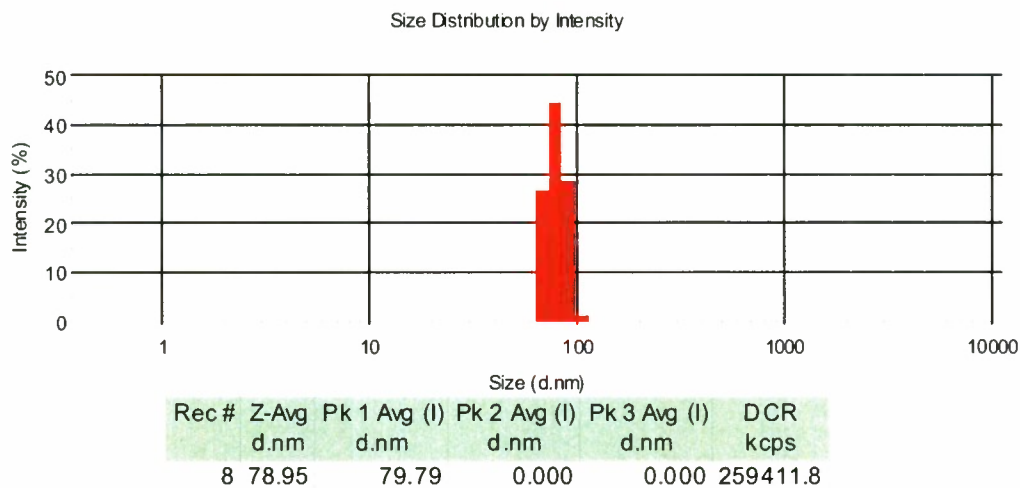


Figure 7. 90 nm PSL (Malvern)

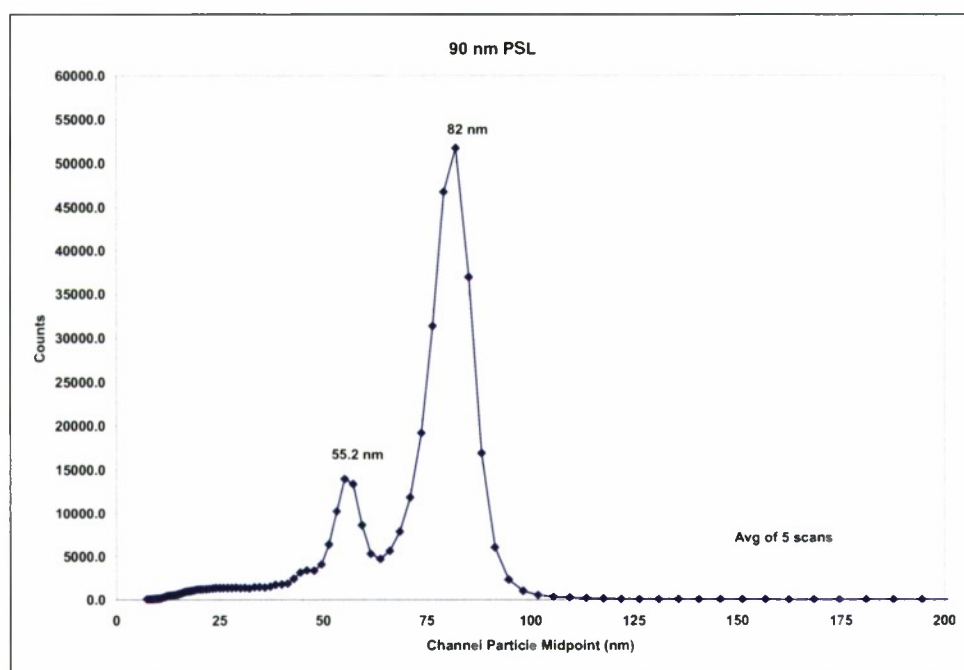


Figure 8. 90 nm PSL (IVDS)

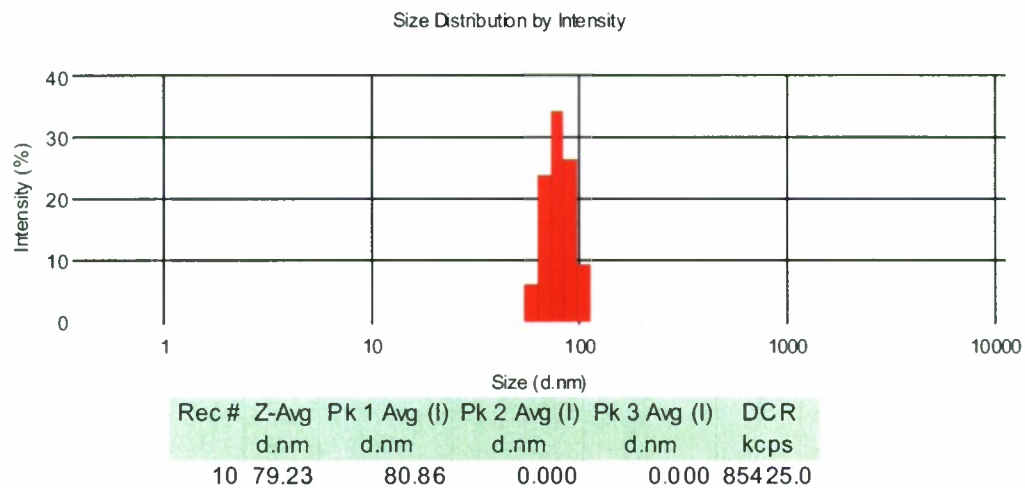


Figure 9. 30 nm Plus 90 nm PSL (Malvern)

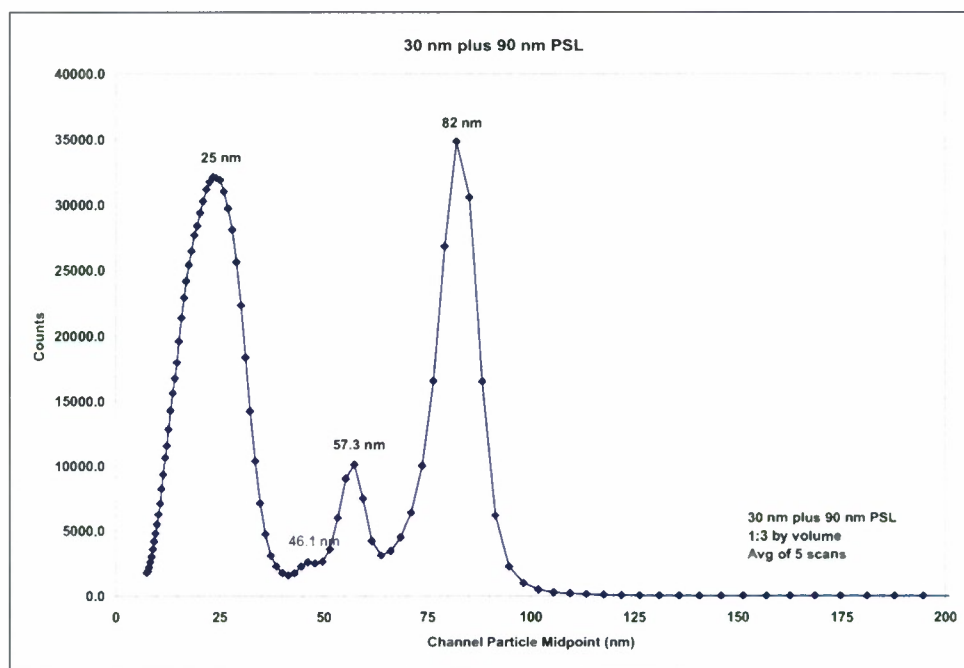


Figure 10. 30 nm Plus 90 nm PSL (IVDS)

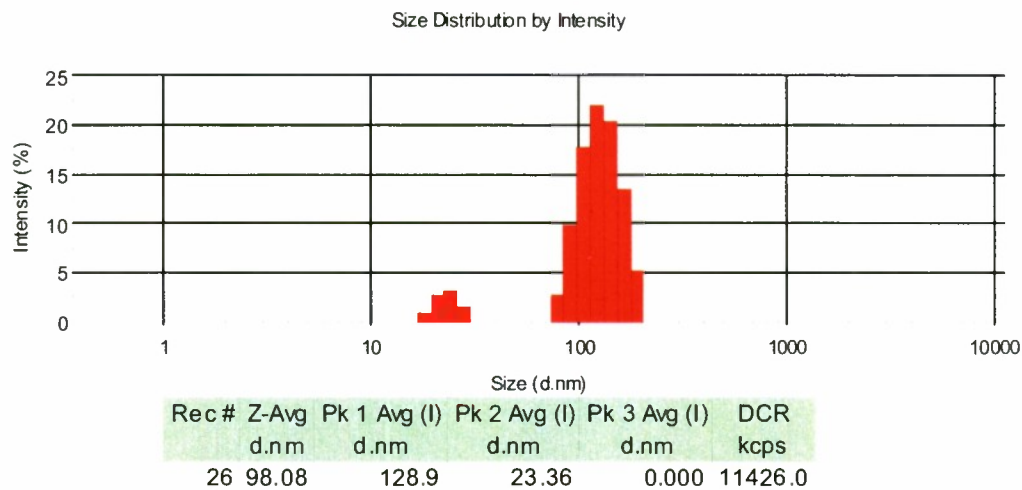


Figure 11. 30 nm Plus 130 nm PSL by Intensity (Malvern)

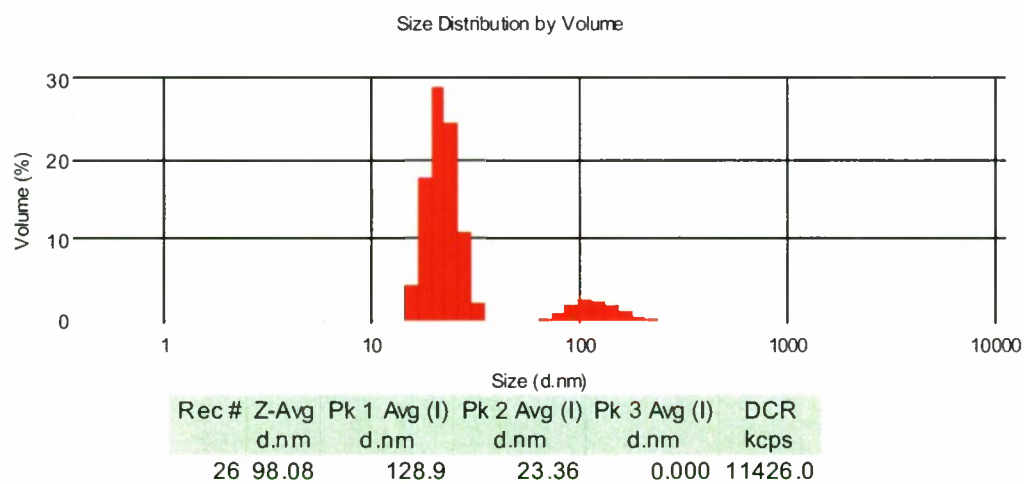


Figure 12. 30 nm Plus 130 nm PSL by Volume (Malvern)

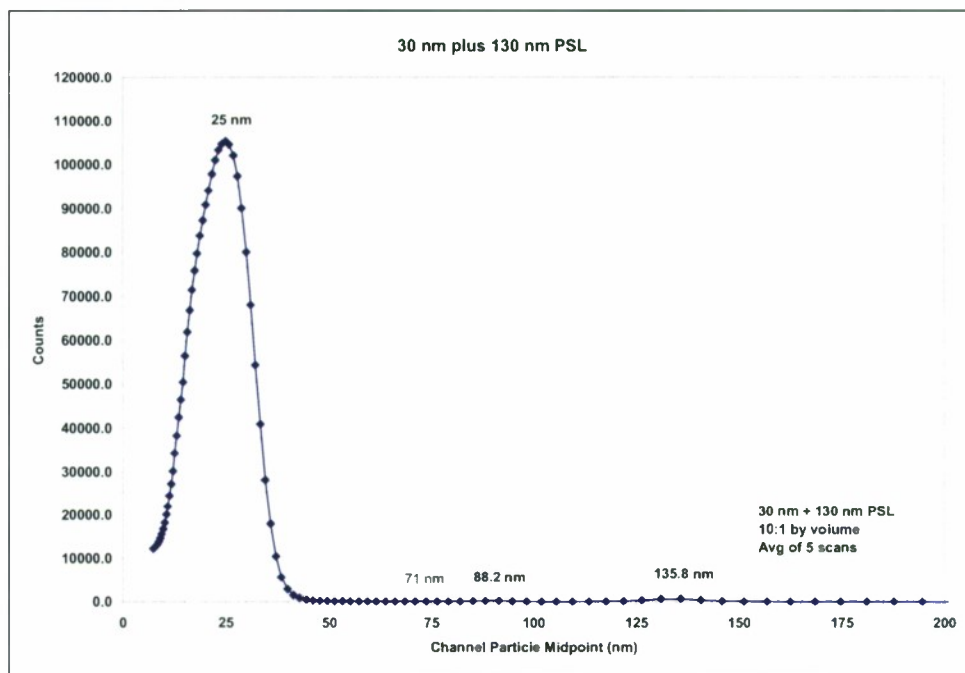


Figure 13. 30 nm Plus 130 nm PSL, Full Scale (IVDS)

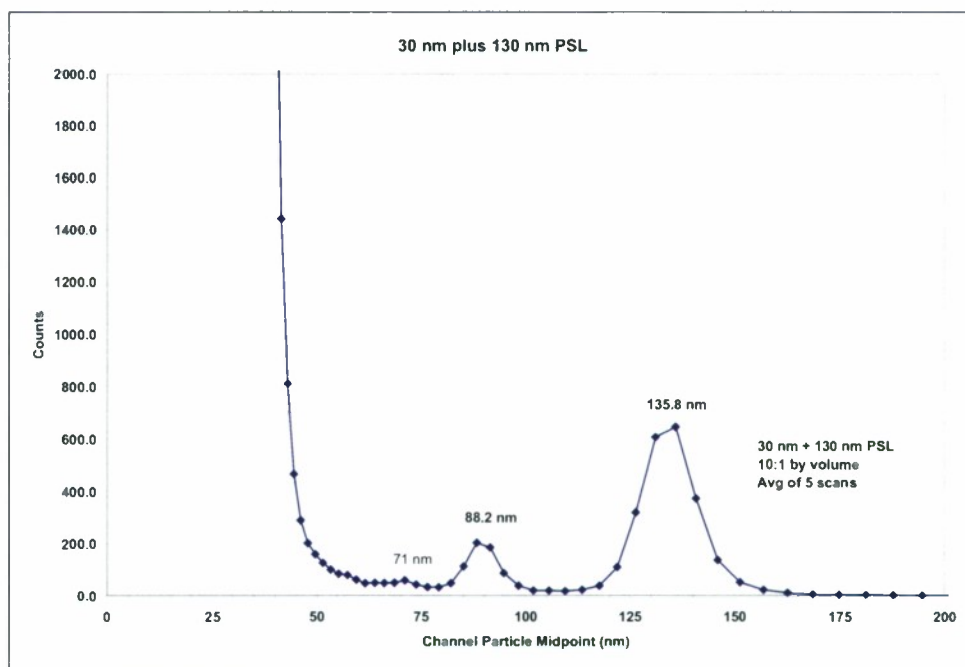


Figure 14. 30 nm Plus 130 nm PSL, Expanded Scale (IVDS)

Blank

LITERATURE CITED

1. Kaszuba, M.; McKnight, D.; Connah, M.T.; McNeil-Watson, F.K.; Nobbmann, U. Measuring Sub Nanometre Sizes Using Dynamic Light Scattering. *J Nanopart Res.* **2008**, *10*, 823–829.
2. Wick, C.H.; McCubbin, P.E.; Birenzvice, A. *Detection and Identification of Viruses Using the Integrated Virus Detection System (IVDS)*; ECBC-TR-463; U.S. Army Edgewood Chemical Biological Center: Aberdeen Proving Ground, MD, 2005; UNCLASSIFIED Report (AD-A454 377).